

METRIA M92 Benchtop pH/ORP Meter

Introduction

Thank you for selecting the METRIA M92 benchtop pH/ORP meter. This manual provides a step-by-step guide to help you operate the meter, please carefully read the following instructions before use.

Unpacking

Before unpacking, ensure that the current work environment meets following conditions.

- Relative humidity is less than 80%.
- Ambient temperature is greater than 0°C and less than 60°C.
- No potential electromagnetic interference.


The following list describes the standard components of the meter. After the unpacking, please check all components are complete. If any are damaged or missing, please contact nearest distributor.




METRIA M92 pH/ORP Meter




Electrode Arm

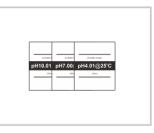

pH Buffer Pouches

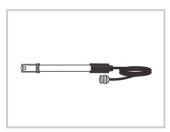

E201 pH Electrode

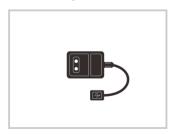

DC5V Power Adapter



TP-10K Temperature Probe



USB Cable

Display

The METRIA M92 benchtop pH/ORP meter is equipped with an easy-read LCD display that used to show the measured values and mode icons. The following table describes the function of each icon.




INDEX:

Measure	Measurement mode icon: Indicates the meter is in the measurement mode.	Stable	Stable icon: Indicates the measuring value has stabilized.
Calibration	Calibration mode icon: Indicates the meter is in the calibration mode.	Hold	Hold icon: Indicates the measuring value has been locked.
Setup	Setup mode icon: Indicates the meter is in the setting mode.	A	Calibration Due Alarm: Prompts the user to calibrate the meter.
Memory	Memory icon: Indicates the data is stored into memory.	ATC	Automatic Temperature Compensation: Indicates the temperature compensation is enabled.
Slope	Electrode slope icon: Indicates the average slope of the pH electrode.		

Keypad

KEY	FUNCTION
Measle	 Switches the meter ON/OFF. Locks the measured value, press the key again to resume measuring. Exits the calibration or setting and returns to measurement.
Mode I°C	Toggles between available measurement modes.Sets the temperature (Press and hold the key for 3 seconds).
Cal	 Starts calibration. Enters the setup menu (Press and hold the key for 3 seconds).
	 Stores current reading to memory. Increase value or scroll up through the menu item.
MRIV	 Views the calibration report or data logs. Decrease value or scroll down through the menu item.
Enter	 Confirms the calibration, settings or displayed options. Turn on/off the backlight (Press and hold the key for 3 seconds).

Connectors

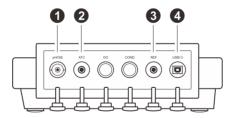


INDEX:

NO.	CONNECTOR	DESCRIPTION	
1	pH/ISE	Used for connecting the pH or ORP electrode	
2	ATC	Used for connecting the temperature probe	
3	REF	Used for connecting the reference electrode	
4	USB ()	Used for connecting the USB cable and DC5V power adapter	

Installing the Electrode Holder


Take out the electrode arm from the packaging. The base plate of the electrode holder has a circular hole, the electrode arm has a connecting rod. Insert the connecting rod into the circular hole and swivel the electrode arm 90°. Electrode holder is now ready to swing into desired position.

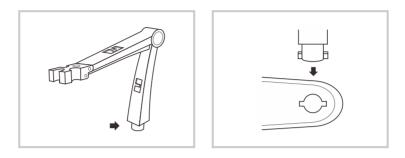

Adjustment of electrode arm


After installation, if the electrode arm automatically rises or falls, you need to adjust the screws until arm locate at any position.

- 1. Remove the plastic cover from the electrode arm.
- 2. Use the screwdriver to tighten the screw moderately.
- 3. Insert the plastic cover to previous position. Installation is completed.

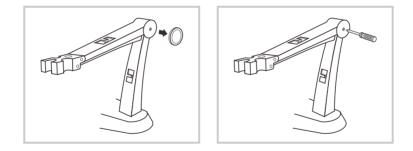


Connecting the Electrode

1. Take out the pH electrode from the packaging. Place the electrode into left or right side of the electrode arm.

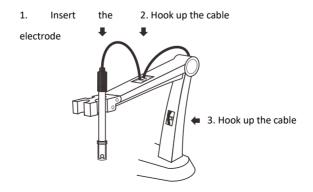

2. Insert the BNC connector into the connector socket labeled pH/ISE. Rotate and push the connector clockwise until it locks. After the connection is completed, DO NOT pull on the cable. Always make sure that the connector is clean and dry.

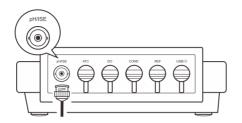


Connecting the Temperature Probe

1. Place the temperature probe into the circular hole of the electrode arm.




2. Insert the phone plug to the connector socket labeled ATC.

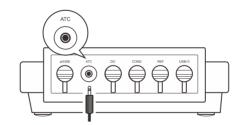

Connecting the Power Adapter

- 1. Connect the USB cable to power adapter.
- 2. Insert the other side of cable into the power socket. The meter is now ready for use.

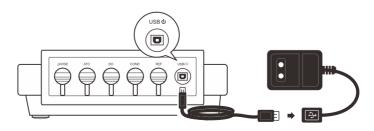

Prior to Use

Remove the protective cap from the bottom of the electrode.

pH Electrode:

If the glass sensitive membrane has dried out, soak the electrode in 3M KCL solution (pH adjusted to 4.0) for at least 30 minutes.



 ORP Electrode (purchase separately): If the sensing element has dried out, soak the electrode in 4M KCL solution for at least 20 minutes.

Switching the Meter On and Off

- Press the Meas key to switch on the meter, the display shows the measured value.
- Press and hold the **Meas** key for 5 seconds, the meter will switch off.

Setup Menu

The METRIA M92 benchtop pH/ORP meter contains an integrated setup menu that is used to customize the displayed option to meet measurement requirements. In the different modes, the meter will show the corresponding options. The following table describes the functions of the menu items.

MENU	DESCRIPTION	OPTIONS	DESCRIPTION	DEFAULT
		USR	USA (pH1.68/4.01/7.00/10.01/12.45)	
	pH Buffer: Set the pH buffer group for calibration and	N ISE	NIST (pH1.68/4.01/6.96/9.18/12.45)	USA
	auto-recognition.		DIN	
		9 IU	(pH1.09/4.65/6.79/9.23/12.75)	
		USEr	User-Defined	
		1	1 point	
		2	2 points	
Calibration Points:	Calibration Points: Set the number of calibration points.	3	3 points	3 points
		Ч	4 points	
		5	5 points	
rE50	Resolution: Set the resolution of the pH measurement.	0.00 (0.001pH	
		0.0 1	0.01pH	0.001pH
		0. 1	0.1pH	

pH MODE:

חט וד	Measurement Unit:	°۲	Degrees Celsius	ŝ
UT 1C	Set the default temperature unit.	°F	Degrees Fahrenheit	C

METRIA M92 Benchtop pH/ORP Meter

ORP MODE:

MENU	DESCRIPTION	OPTIONS	DESCRIPTION	DEFAULT
r E 5 D	Resolution:	0. 1	0.1mV	0.1)/
rc30	Set the resolution of the ORP measurement.	1	1mV	0.1mV

GENERAL OPTIONS:

MENU	DESCRIPTION	OPTIONS	DESCRIPTION	DEFAULT
51.5	Stability Criteria: When the LO option is enabled, the Stable icon will guickly appear on the display.	LO	Low	
528	When the HI option is enabled, the icon will take longer to appear, but guarantees high accuracy of the measurement.	н	High	Low
נ וחט	Auto-Hold: When the option is enabled, the meter will	¥E5	Enable	Disable
XOLA	automatically sense a stable reading and lock the measurements.	по	Disable	Disable
	Auto-Power Off: When the option is enabled, the meter will automatically turn off if no key is pressed within a	10	10 minutes	
066		20	20 minutes	Disable
UFF		30	30 minutes	
	specified time period.	по	Disable	
	Calibration Due: When the option is enabled, if the meter does not	131	1 to 31 days	
CALL	calibrated within a specified time period, the meter will automatically show the \triangle icon.	OFF	Disable	Disable
937BP	Date and Time: Set the current date and time.		Year-month-day, hour-minutes	
5 .	Clear Stored Data:	YE 5	Enable	Disable
Elr	Delete all stored readings in the memory.	по	Disable	Disable
r St	Factory Reset: If enabled, all of the calibration data and	462	Enable	Disable


selected parameters will back to factory default settings,	по	Disable	
the meter must be recalibrated.			

METRIA M92 Benchtop pH/ORP Meter

Setting the default option

- 1.1 If necessary, press the Mode key until the display shows desired measurement mode (e.g., pH).
- 1.2 Press and hold the 🖹 key for 3 seconds to enter the setup menu and the ▲ or 🔻 key to select the menu item (e.g., BUF/P-01).
- 1.3 Press the **Enter** key, the display shows an option.
- 1.4 Press the \blacktriangle or \blacktriangledown key to select the desired option.
- 1.5 Press the Enter key to confirm, the meter returns to the measurement mode. Setting is completed.

() If you want to exit the setting, press the **Meas** key.

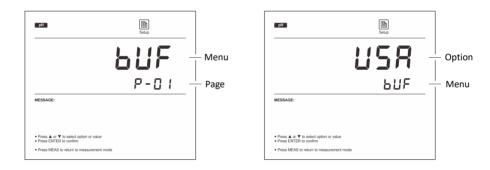

Setting the date and time

2.1 Press and hold the \exists key for 3 seconds to enter the setup menu and the \blacktriangle or \forall key until the display shows "DATE".

2.2 Press the Enter key, the meter shows current year (e.g., 2018).

2.3 Press the \blacktriangle or \checkmark key to set year and the **Enter** key to confirm, the display shows current date and time (Format: month-day, hour-minutes).

2.4 Press the \blacktriangle or \blacktriangledown key to set the date and time, press the **Enter** key to confirm until the meter returns to the measurement mode. Setting is

completed.



8

Temperature Compensation

For better accuracy, we recommend the use of either a sensor with a built-in or a separate temperature probe for the calibration or measurements.

Automatic Temperature Compensation

Connect the temperature probe to the meter (Refer to page 5 "Connecting the Temperature Probe"). The ATC icon immediately appears on the display, the meter is now switched to the automatic temperature compensation mode.



Manual Temperature Compensation

If the meter does not detect a temperature probe, the °C icon will show on the display indicating that the meter is switched to the manual temperature compensation mode. To set the temperature value of sample, follow the steps below.



- 1. Press and hold the °C key for 3 seconds to enter the temperature setting mode.
- 2. Press the \blacktriangle or \triangledown key to modify the temperature value.
- 3. Press the Enter key to confirm.
- Press the ▲ or ▼ key once, the setting value will increase or decrease by 0.1. Press and hold the ▲ or ▼ key, the setting value will increase or decrease by 1.

pH Calibration

The METRIA M92 benchtop pH/ORP meter allows 1 to 5 points calibration in the pH mode. We recommend that you perform at least 2 points calibration for high accuracy measurement. The meter will automatically recognize and calibrate to following standard buffer values.

USA Standard Buffers	pH1.68, 4.01, 7.00, 10.01, 12.45
NIST Standard Buffers	pH1.68, 4.01, 6.86, 9.18, 12.45
DIN Standard Buffers	pH1.09, 4.65, 6.79, 9.23, 12.75

If the USER option is selected, the meter will allow only 2 points calibration. Single point calibration should only be carried out with pH7.00, 6.86 or 6.79, otherwise calibration will not be accepted.

Make sure to calibrate the meter when attaching a new electrode or during first use. DO NOT reuse the calibration solution after calibration, contaminants in solution will affect the calibration and eventually the accuracy of the measurement.

Single point calibration



- 1.1 Press the **Mode** key until the meter shows **PH** icon and you have selected 1 point calibration in the setup menu.
- 1.2 Press the Cal key, the display shows pH7.00/CAL (or 6.86/CAL, or 6.79/CAL).

рН	- Calbration
	ERL
MESSAGE:	
Press MEAS to return to measurement mode	

1.3 Rinse the pH electrode with distilled water, place the electrode (and temperature probe) into the pH7.00 (or 6.86, or 6.79) buffer solution. The end of the electrode must be completely

immersed into the calibration solution. Stir the electrode gently to create a homogeneous

solution. Press the **Enter** key, the Calibration icon begins flashing.



1.4 Wait for the reading to stabilize, the meter automatically shows END and returns to the

measurement mode. Calibration is completed.

METRIA M92 Benchtop pH/ORP Meter

Multi-point calibration

pH	Calibration
-	
	CALS
MESSAGE:	
Immerse the sensor in the calibration solution	
Press MEAS to return to measurement mode	

2.1 Ensure that you have selected 2 to 5 points calibration in the setup menu.

2.2 Repeat the steps 1.2 to 1.3 above. When the first calibration point is completed, the display will show CAL2. The meter prompts you to continue with second point calibration.



2.3 Rinse the pH electrode with distilled water, place the electrode (and temperature probe) into the next buffer solution. The meter will automatically recognize the calibration solution (e.g., pH4.01) and begins the calibration, the Calibration icon continuously flashing.



2.4 Wait for the reading to stabilize, the display will show CAL3. The meter prompts you to

continue with third point calibration.

2.5 Repeat the step 2.3 above until the display shows END. The meter automatically returns to the measurement mode. Calibration is completed.

pH calibration with custom buffers

рН	Setup
L	JSEr
	ЬUF
MESSAGE:	
 Press A or V to select option or value Press ENTER to confirm 	
Press MEAS to return to measurement	mode

рН	Calibration
	5.00
	ERL I
MESSAGE:	
Immerse the sensor in the calibration solution Press ▲ or ▼ to select option or value Press ENTER to confirm	
Press MEAS to return to measurement mode	

рН	Calibration
	4 <u>.0</u> 0
	ERL2
MESSAGE:	
 Immerse the sensor in the calibration solution Press ▲ or ▼ to select option or value Press ENTER to confirm 	
Press MEAS to return to measurement mode	

3.1 Ensure that you have selected the USER option in the setup menu. The calibration solutions should be at least 1 pH unit apart from each other.

3.2 Rinse the pH electrode with distilled water, place the electrode (and temperature probe) into the custom buffer solution. Stir the electrode gently and wait until the measurement is stable.

3.3 Press the **Cal** key, the meter enters the calibration mode.

3.4 If necessary, press the \blacktriangle or \checkmark key to set the calibration value, press the **Enter** key to begin the calibration (e.g., 6.00pH).

- 3.5 Wait for the reading to stabilize, the display shows CAL2. The meter prompts you to
 - continue with second point calibration.
- 3.6 Rinse the pH electrode with distilled water, place the electrode (and temperature probe) into the next buffer solution and wait until the measurement is stable.

3.7 If necessary, press the \blacktriangle or \checkmark key to set the calibration value, press the Enter key to begin the calibration (e.g., 4.00pH).

estance:

3.8 Wait for the reading to stabilize, the meter automatically shows END and returns to the

measurement mode. Calibration is completed.

$\mathbf{\hat{U}}$

- During the calibration process, if the meter shows \mathcal{E}_{rr} , please check the pH electrode and ensure the pH buffers are fresh and uncontaminated.
- If the electrode slope is not within the normal range (< 70% or >110%), the Slope=II icon will disappear on the display.
- If you want to exit the calibration, press the Meas key.

Viewing the pH calibration report

рН	Memory
	818
	P-02
MESSAGE:	
Press ▲ or ▼ to select option or value Press ENTER to confirm	
Press MEAS to return to measurement mode	

рН	Memory
	12-31
MESSAGE:	
 Press ▲ or ▼ to select option or 	r value
Press MEAS to return to measure	rement mode

pH		H Memory
	0F	5
MESSAGE:		
 Press ▲ or ▼ to select option or value 		
Press MEAS to return to measurement mode		

pH	Memory
	4-7
	99.8 ×
MESSAGE:	
 Press ▲ or ▼ to select option or value 	
Press MEAS to return to measurement mode	

- 4.1 Press the **MR** key in the pH measurement mode, the meter shows LOC/P-01.
- 4.2 Press the \blacktriangle or \checkmark key until the meter shows ELE/P-02.

4.3 Press the Enter key, the meter shows the last calibration date (Format: month-day).

4.4 Press the $\mathbf{\nabla}$ key, the meter shows the zero-point offset (e.g., 2mV).

4.5 Press the \checkmark key, the meter shows the pH buffer group and slope (e.g., pH4~7, slope: 99.8%).

- 4.6 To exit the calibration report, press the **Meas** key.
- () If the meter does not calibrated, the display will only show "----".

ORP Calibration

ORP

MESSAGE:

· Press MEAS to return to measurement mode

The METRIA M92 benchtop pH/ORP meter allows 1 point calibration in the relative mV mode, but calibration is not necessary unless exact readout agreement with a work standard and at a specific ORP value is needed.

MESSAGE:
Press ▲ or ▼ to select cotion or value Press ENTER to confirm Press MEAS to return to measurement mode

6

105.0

1.1 Press the **Mode** key until the meter shows **ORP** icon.

1.2 Rinse the ORP electrode with distilled water, place the electrode into the calibration solution. Stir the electrode gently and wait until the measurement is stable.

1.3 Press the **Cal** key, the meter enters the calibration mode.

- 1.4 Press the \blacktriangle or \triangledown key to set the displayed value (e.g., 105 R.mV).
 - 1.5 Press the **Enter** key to confirm, the Calibration icon begins flashing.

ORP	Calibration
	End
MESSAGE:	

1.6 Wait for the reading to stabilize, the meter automatically shows END and returns to the

measurement mode. Calibration is completed.

Viewing the ORP calibration report

ORP	Calibration
	E: E
	P-02
MESSAGE:	
Press ▲ or ▼ to select option or value Press ENTER to confirm	
Press MEAS to return to measurement mode	

ORP	Calibration
	1231
MESSAGE:	
Press ▲ or ▼ to select optio	n or value
Press MEAS to return to mea	ssurement mode

ORP	Calibration	
	BF5	
MESSAGE:		
Press ▲ or ▼ to select option or value Press MEAS to return to measurement mod	e	

- 2.1 Press the **MR** key in the ORP measurement mode, the display shows LOC/P-01.
- 2.2 Press the \blacktriangle or \checkmark key until the display shows ELE/P-02.

2.3 Press the **Enter** key, the display shows the last calibration date (Format: month-day).

- 2.4 Press the $\mathbf{\nabla}$ key, the display shows the offset potential (e.g., 5mV).
- 2.5 To exit the calibration report, press the **Meas** key.

Temperature Calibration

During the measurement process, if the temperature reading displayed differs from that of an accurate thermometer, the meter needs to be calibrated.

- 1. Connect the temperature probe to the meter and place into a solution with a known accurate temperature.
- 2. Press and hold the °C key for 3 seconds to enter the temperature setting mode.
- 3. Press the \blacktriangle or \triangledown key to set the temperature value.
- 4. Press the **Enter** key to confirm. Calibrating is completed.



During the setting process, press the \blacktriangle or \blacktriangledown key once, the setting value will increase or decrease by 0.1. Press and hold the \blacktriangle or \blacktriangledown key, the setting value will increase or decrease by 1.

pH Measurement

1. Press the **Mode** key until the display shows **PH** icon.


2. Rinse the pH electrode with distilled water. Place the electrode (and temperature probe) into the sample solution, stir the electrode gently.

3. Record the measured value when the reading is stable.

ORP Measurement

The METRIA M92 benchtop pH/ORP meter provides two millivolt measurement modes.

- Press the Mode key until the display shows measurement unit "mV", the meter is now enters the absolute mV measurement mode.
- Press the Mode key until the display shows ORP icon, the meter is now enters the relative mV measurement mode.

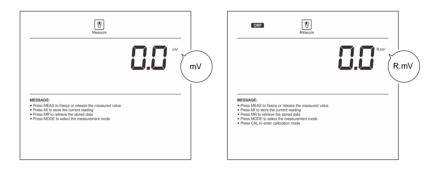

• Select one of the above modes. Place the ORP electrode into the sample. Record the measured value when the reading is stable.

Auto-Hold

The meter contains an Auto-Hold function. If enabled, the meter will automatically sense a stable reading and lock the measurements, the HOLD icon

appears on the display. If disabled, press the â key, the meter will immediately lock the displayed value. Press the **Meas** key to resume measuring.



Storing and Recalling Data

На	Memory Memory
7.0	
Stope and	35.0 Å
MESSAGE: • Press MEAS to freeze or release the measured value • Press MI to store the current reading • Press MI to retrieve the stored data • Press MODE to select the measurement mode • Press GAI to redre calization mode	

The METRIA M92 benchtop pH/ORP meter is capable of storing and recalling up to 500 data sets.

Storing readings into memory

During the measurement process, press the **MI** key to store the reading into the memory, the Memory icon appears on the display.

Viewing stored readings

- Press the **MR** key in the measurement mode, the meter shows LOC/P-01 (Data Log). 1.
- 2. Press the Enter key, the meter shows the serial number of the stored data.

Press the $\mathbf{\nabla}$ key, the meter shows the date and time of the stored data (Format: 3. month-day, hour-minutes).

- 4. Press the $\mathbf{\nabla}$ key, the meter shows the stored data. 5. Press the $\mathbf{\nabla}$ key again, the meter shows next data set.
 - 6. Press the Meas key, the meter returns to the measurement mode.

Clearing the memory

Please refer to page 6 SETUP MENU.

pH	Memory
	° 0.25
MESSAGE:	



12-3

Β

10:59

pН

MESSAGE:

 Press ▲ or ▼ to select option or value Press MEAS to return to measurement m 18

Electrode Care and Maintenance

pH electrode

Since pH electrode is susceptible to dirt and contamination, clean as necessary depending on the extent and condition of use.

- After measuring: rinse the electrode in distilled water, store the electrode into the 3M KCL solution.
- Salt deposits: soak the electrode in warm tap water to dissolve deposits, then thoroughly rinse with distilled water.

• Oil or Grease film: wash the glass sensitive membrane of electrode gently in some detergents and water. If necessary, using the alcohol to clean the sensitive membrane, then rinse with distilled water. Place the electrode in the 3M KCL solution for at least 30 minutes.

• Clogged reference junction: heat a diluted KCl solution to 60°C to 80°C. Place the electrode into the heated solution for about 10 minutes. Allow the electrode to cool in some unheated KCl solution.

• Protein deposits: prepare a 1% pepsin solution in 0.1M of HCL. Place the electrode in the solution for 10 minutes. Rinse the electrod with distilled water.

Rectivating the pH Electrode:

If stored and cleaned properly, the electrode should be ready for immediate use. However, a dehydrated sensitive membrane may cause sluggish response. To rehydrate the sensitive membrane, immerse the electrode in a pH4.01 buffer solution for 10 to 30 minutes. If this fails, the electrode requires activation.

- 1. Soak the electrode in 0.1M HCl for 5 minutes.
- 2. Remove and rinse with deionized water, then place in 0.1M NaOH for 5 minutes.
- 3. Remove and rinse again, and soak in 3M KCL solution for at least 30 minutes.

ORP electrode

Ensure that the ORP electrode is thoroughly washed with distilled water after use.

• In aggressive chemicals, dirty or viscous solutions, and solutions with heavy metals or proteins, take readings quickly and rinse electrode

immediately.

If you do not use the electrode for long periods, store the electrode with 4M KCL solution.

Cleaning the Electrode:

Contamination of the sensing element often results in slow response and inaccurate readings. If necessary, clean the element by one of the following

procedures.

Inorganic Deposits:

- 1.1 Soak the ORP electrode in 0.1M HCl for 10 minutes.
- 1.2 Remove and rinse with distilled water, then place in alcohol for 5 minutes.
- 1.3 Remove and rinse again, and soak in pH4.01 buffer solution for 15 minutes.

Oil and Grease Films:

- 2.1 Wash the electrode gently in some detergents and water.
- 2.2 Dip the electrode in the 4M KCL solution for at least 30 minutes.

Troubleshooting

LCD DISPLAY	CAUSE	CORRECTIVE ACTION			
	Electrode dried out	Soak the pH electrode in 3M KCL solution at least 30 minutes.			
	Measured value is out of range	Check the electrode whether clogged, dirty or broken.			
5	Incorrect calibration solutions	Using the fresh calibration solutions for calibration.			
Err	pH electrode has expired	Replace the pH electrode.			

Specifications

	Model	METRIA M92				
рН	Range	-2.000~20.000pH				
	Accuracy	±0.002pH				
	Resolution	0.1, 0.01, 0.001pH				
	Calibration Points	1 to 5 points				
		USA (pH1.68/4.01/7.00/10.01/12.45)				
	pH Buffer Options	NIST (pH1.68/4.01/6.86/9.18/12.45)				
		DIN (pH1.09/4.65/6.79/9.23/12.75)				
	Temperature Compensation	0~100°C, 32~212°F, Manual or Automatic				
	Range	-1999.9~1999.9mV				
mV	Accuracy	±0.2mV				
	Resolution	0.1, 1mV				
	Calibration Points	1 point (Only for relative mV mode)				
	Range	0~105°C, 32~221°F				
T	Accuracy	±0.5°C				
Temperature	Resolution	0.1°C				
	Calibration Points	1 point				
General	Memory	Stores up to 500 data sets				
	Output	USB communication interface				
	Connector	BNC				
	Display	LCD				
	Operating Temperature	0~60°C				
	Relative Humidity	< 80%				
	Power Requirements	DC5V, using AC adapters, 220VAC/50Hz				
	Dimensions	210 (L) × 188 (W) × 60 (H)mm				
	Weight	1.5kg				

Addendum 1: pH Electrode Selection Guide

The METRIA M92 benchtop pH/ORP meter comes with a general purpose pH electrode that is used to measure the pH of the liquids. If this electrode can not meet your measurement requirements, please refer to the table below to select an applicable probe.

SAMPLE TYPE	P11	P12	P13	P15	P16	P18	P19	P21	E201	E202
Agar										•
Beer	•	•	•					•	•	•
Blood Products	•	•	•					•		•
Bread, Dough						•	•			
Cement	•									
Cosmetics	•	•	•					•	•	•
Dairy Products	•	•	•				•			•
Education	•								•	•
Fats/Cream							•			
Field Use						•			•	•
Fish Products							•			•
Lab Flasks		•								
Low Ionic	•			•				•		
Meat, Cheese							•			•
Micro Samples			•							
Paint		•	•							•
Photographic										
Soil						•	•			
Surface										•
Test Tubes		•			•					
Tris Buffer					•					
Viscose Samples										•

Addendum 2: ORP Electrode Selection Guide

ORDER CODE	APPLICATION
501	Suitable for the sample with strong redox potential, plastic body, temperature range: 0~80°C
502	Suitable for the sample with weak redox potential, plastic body, temperature range: 0~80°C
504	Suitable for the high temperature samples, glass body, temperature range: 0~100°C

Addendum 3: Preparation of pH Buffer Solutions

• Open the pH7.00 buffer packet, place the reagent into a 250ml volumetric flask. Pour the distilled water 250ml to scale line, mix the solution until the reagent is completely dissolved.

• Preparation of pH4.01 and 10.01 standard buffer solutions are the same as above. Prepared standard buffer solutions should be stored in hermetically sealed glass containers.



Addendum 4: Preparation of ORP Standard Solutions

Add 3 grams of quinhydrone to 500ml buffer pH4.01 and stir for 15 minutes. Un-dissolved quinhydrone powder must be present.
 Potential @ 25°C =+263mV (±10mV)

• Add 3 grams of quinhydrone to 500ml buffer pH7.00 and stir for 15 minutes. There must be an excess of undissolved quinhydrone powder.

Potential @ 25°C =+87mV (±10mV)

Hazardous Substance Statement

METRIA Instruments is committed to the reduction and eventual elimination of all hazardous substances in both the manufacturing process and finished products we supply. We have an active manufacturing and procurement program to minimize and eliminate the use of harmful heavy metals such as cadmium, lead, mercury and the like. New technologies and design parameters are also promoting these efforts and we expect to have little or no such materials in our product in the coming years. We welcome our customer suggestions on how to speed up these efforts.



Warranty

The warranty period for meter is one year from the date of shipment. Above warranty does not cover the sensor and calibration solutions. Out of warranty products will be repaired on a charged basis. The warranty on your meter shall not apply to defects resulting from:

- Improper or inadequate maintenance by customer.
- Unauthorized modification or misuse.
- Operation outside of the environment specifications of the products.

For more information, please contact the nearest authorized distributor.