



Exported to
106 countries



3 years warranty




PERISTALTIC PUMP

Innofluid Co.,Ltd

APPLICATION



Laboratory chemical dosing

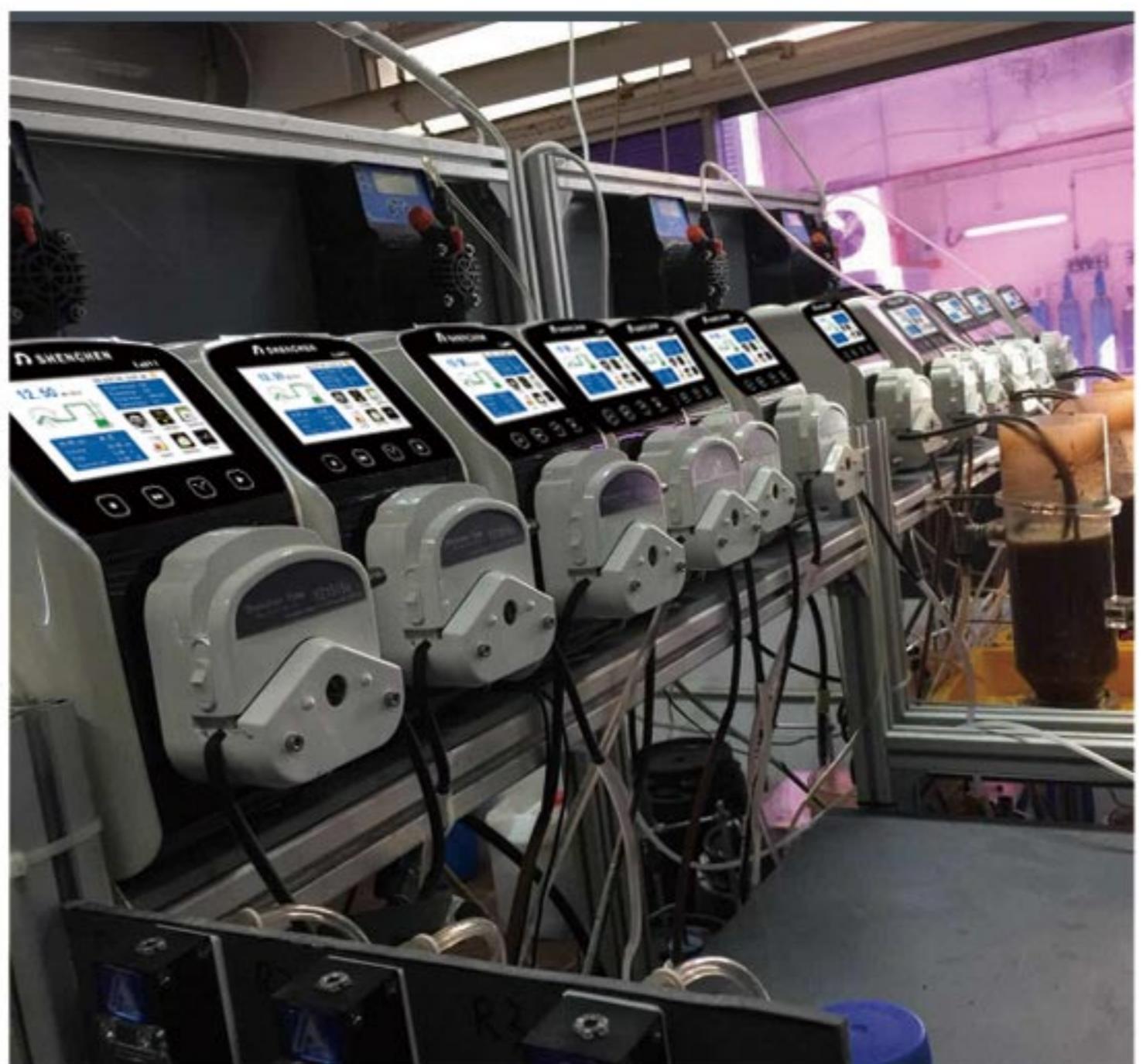

Analytical instrument sampling



Waste water treatment



Chromatography



Cosmetic filling



Pharmaceutical filling



Food & beverage filling



Diagnostic reagents filling



SHENCHEN

Pump Head Flow Rate

	0.1~150rpm	0.0002~65.17mL/min
AMC Series		
	0.1~350rpm	0.004~149.23mL/min
MicroPump		
	0.1~300rpm	0.0024~190mL/min
MiniPump		
	0.1~300rpm	0.0033~365.69mL/min
HandyPump		
	0.1~600rpm	0.0033~560.04mL/min
KT15		
	0.1~350rpm	0.08~930mL/min
UD15		
	0.1~600rpm	0.024~1500mL/min
SN Series		
	0.1~600rpm	0.0053~3100mL/min
EasyPump Series		
	0.1~600rpm	0.211~3600mL/min
DZ25-3L		
	0.1~350rpm	0.01~4340mL/min
DY Series		
	0.1~600rpm	0.3~6000mL/min
DZ25 Series		
	0.1~600rpm	0.0875~6663mL/min
UC25		
	0.1~600rpm	0.69~12000mL/min
YZ35		
	30~350rpm	1.0~28.15L/min
DZ45		



Features:

The elastic positioning mechanism enables users remove and install cartridge with one hand.
Stepless adjustment of the tube pressure gap, effectively improving the flow rate accuracy between channels.

NEW New Generation


Multichannel Type Pump Head




Features:

The tube clamp linkage mechanism makes it more convenient to install the tube.
The rubbing wheel adjustment mechanism can easily fix different size tube.
The lever assist mechanism makes the operation more labor-saving.
The tube pressure gap fine-tuning mechanism can adjust the pressure, extend the life of the tube and improve the dispensing accuracy.
Tube tubes can be installed at the same time to realize single pump head with two channels.

NEW New Generation
Quick Load Type
Pump Head

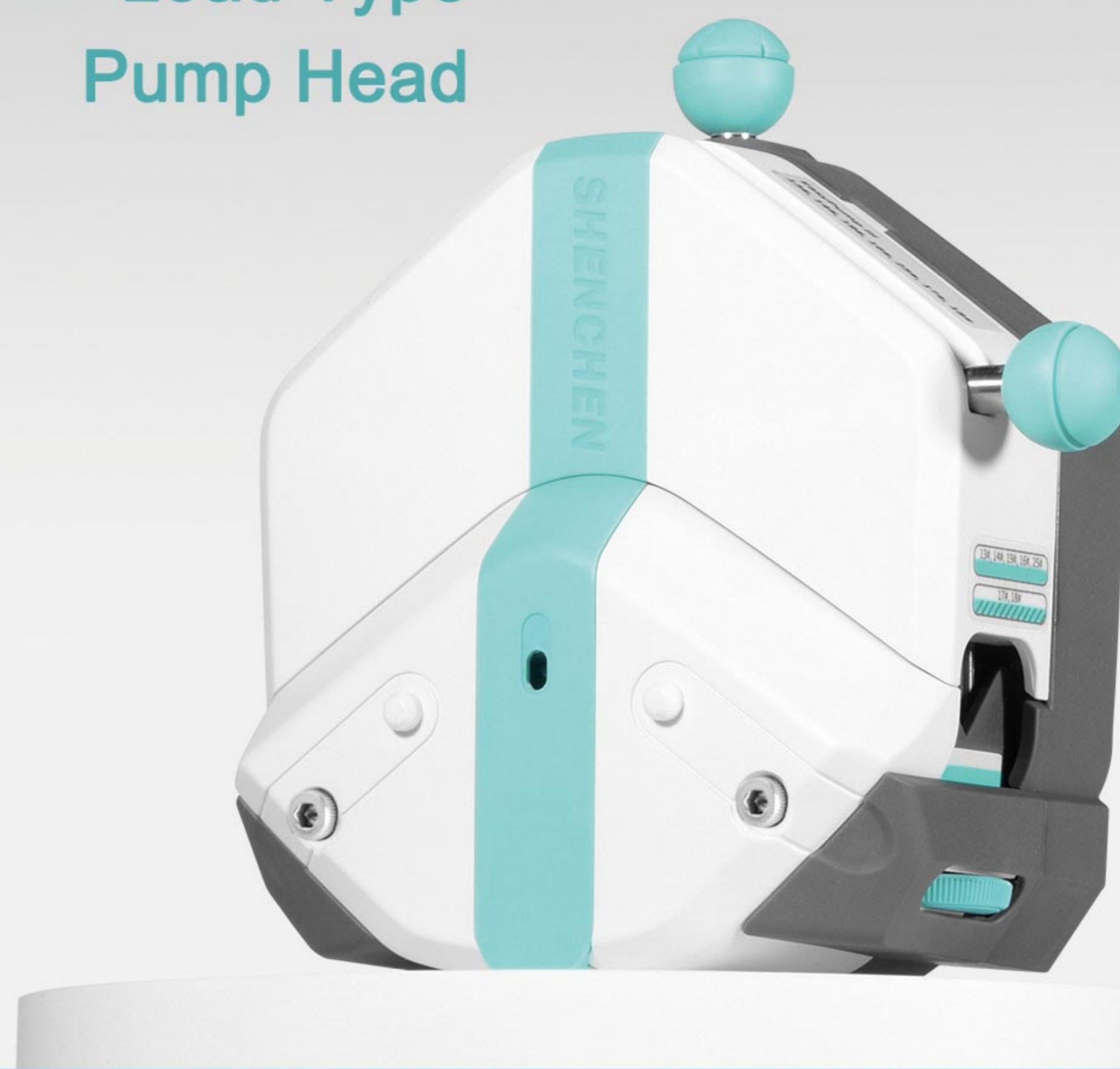

Dual channel

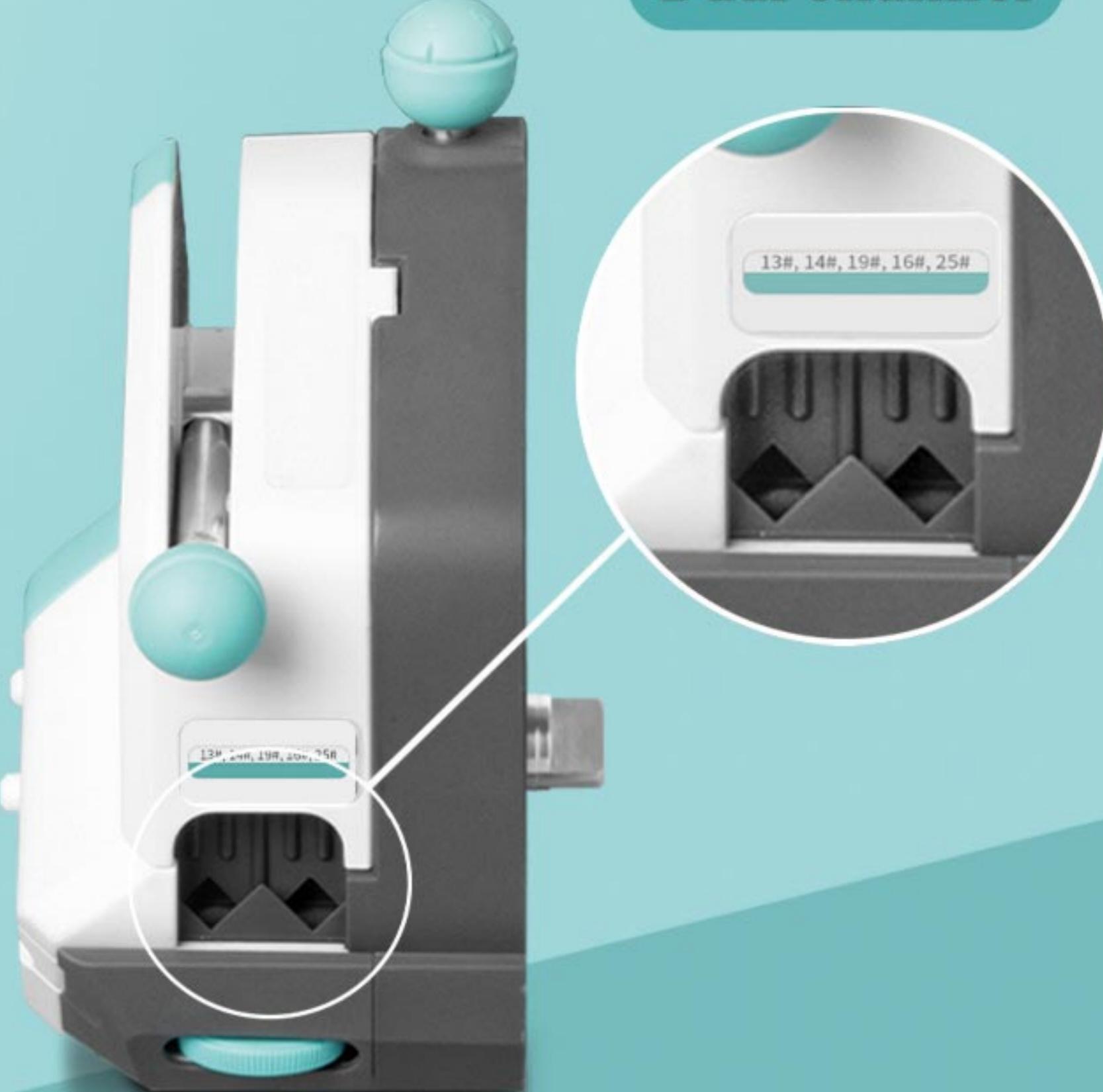

Single channel





DIRECTORY



COMPANY
APPLICATION
NEW PRODUCT
PRODUCTS

◆ OCM Fluid Solution.....	01
◆ Compact Peristaltic Pump.....	06
LabQ.....	06
LabK1.....	07
LabT.....	08
ST-HandyPump.....	09
SP-MiniPump.....	10
◆ Flow Rates Peristaltic Pump.....	11
LabV/V Series.....	11
V6 Series	13
Lab-IV Series.....	15
LabN Series.....	16
N6 Series	18
◆ Dispensing Peristaltic Pump.....	19
LabF Series/F Series.....	19
F6 Series.....	22
IF3 Low Pulsation.....	25
◆ Filling system.....	26
KF300 Series.....	26
Split Type CF600 PlusII.....	28
Split Type CF350 Plus.....	29
Integrated Type DF600 PlusII.....	30
Integrated Type DF600 PlusIV/DF350 Plus.....	31
◆ Basic Peristaltic Pump.....	35
LabM Series.....	35
M6 Series.....	36
BT-N Series.....	37
◆ Planetary Gear Industrial Pump.....	38
J025.....	38
◆ Explosion Proof Peristaltic Pump.....	39
EXP600.....	39
QD600.....	40
◆ Peristaltic Pump Tubing and Accessory.....	43

What is ODM


ODM (Original Design Manufacturer)

- Peristaltic pump core technology
- ODM System Solutions
- User-defined work mode



PROBLEMS

Conventional products can't be solved

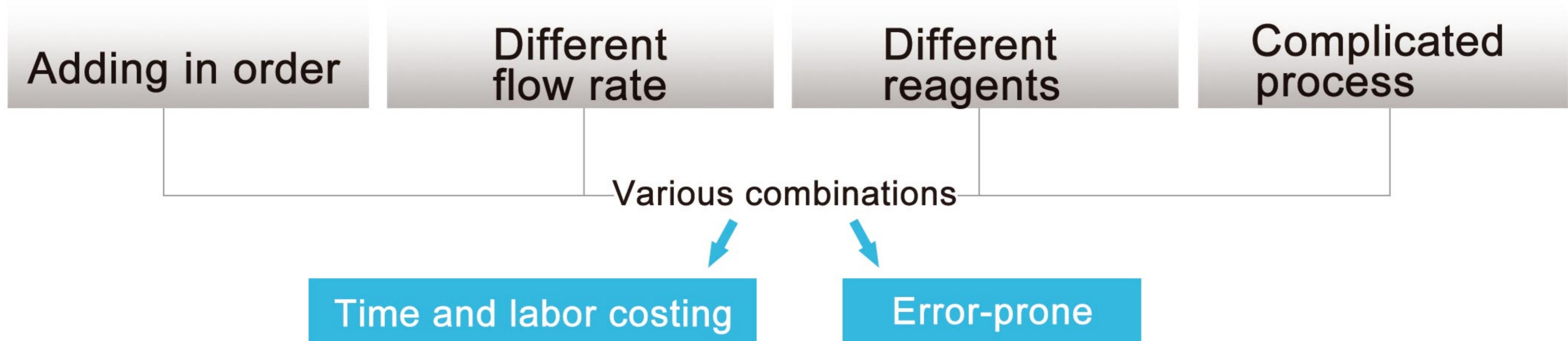

How to solve the problems?

ODM provides system solutions for you.

1 Analyze problem 2 Design solution 3 Design ODM product

Our advantages

CUSTOMIZED Professional engineer offer one to one customized service.

EXPERIENCE 16 years of fluid transmission.

TECHNOLOGY Peristaltic pump technology research center.



ODM Application



01 Printing and Dyeing Industry

Pigment accurate allocation, use several different pump heads to allocate different colors of pigment.



02 Pharmaceutical Industry

Filling one bottle with several times during filling process to prevent foaming.


03 Reagent Dispensing

It can fill a variety of different liquids in same time. About 5000-6000 pcs 96 deep well plates can be filled every day.



04 Fermentor Supporting

Peristaltic pump used in fermentor supporting, 4 Handypump head used for adding acid, alkali, antifoaming agents and nutritional agents. The Lab series pump used for waste discharge.

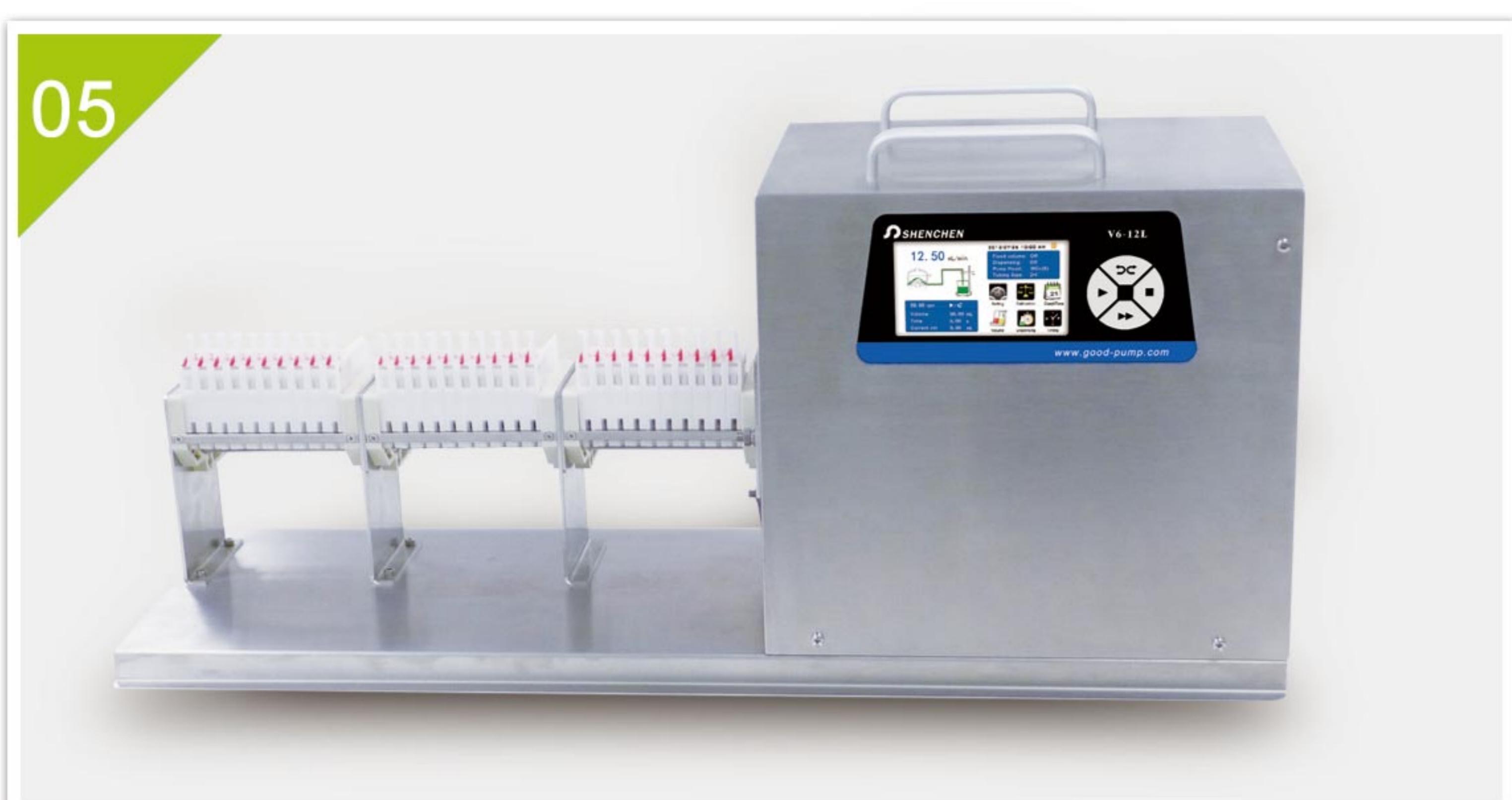

05 Ink Cartridge Filling

30 channels peristaltic pump filling 30 ink cartridges in same time, high accuracy and high efficiency.

Nucleic acid detection reagent dispensing system



INTRODUCTION

This product is an intelligent platform for automatic separation of nucleic acid detection reagents. It has achieved high precision automatic batch production of reagents.

FEATURES

8 different reagents are packed in unequal or **equal amounts**.

1

Up to 0.5% packing accuracy.

2

7 inch industrial true color LCD touch screen operation, man-machine interface is friendly, can store more than a group of commonly used work mode, simple and convenient.


3

The utility model has **two-way quick liquid filling function**, can be used for bidirectional operation and is convenient for cleaning the hose, and can quickly fill the reagent in the preparation stage or the packing line before and after the packing.

4

Mechanical arm operation program can be edited to apply different size, different shape of reagent box, accurate positioning, production efficiency.

5

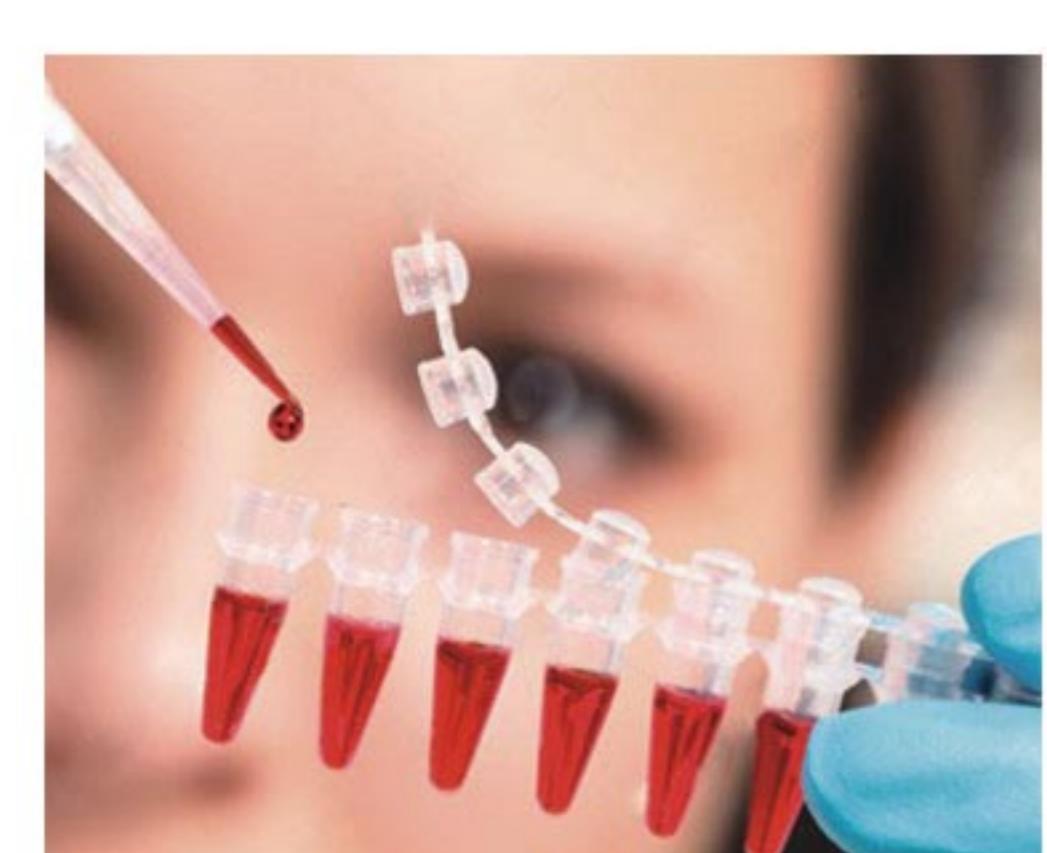

Gene testing reagent packing



Microporous plate packing



Reagent box packing



Micro reagent packing



Biological reagent packing

APPLICATION

This product is used for gene testing reagent packing, microporous plate packing, reagent box packing, micro reagent packing, biological reagent packing and so on.

Programmable Dispensing System



PDS

Introduction

- | One controller can control maximum 8 pump units.
- | Each pump unit with same or different pump head.
- | The OCM controller and pump unit, can be integrated type or split type.

Four working modes


- | Logic working mode: PDSI
- | Independent working mode: PDSII
- | Independent working mode and steps filling mode: PDSIII
- | Comprehensive working mode: PDSIV

Suitable pump head

- | AMC series, MC series, MiniPump, UD15, HandyPump, KT15
- | EasyPump series, YZ1515x, YZ2515x, YZ35
- | DZ25-3L, DZ25-6L, DY15, DY25

Technical Specifications


Speed range	0.1~600rpm, also depend on pump head	Power supply	AC 220V±10%, AC 110V±10%
Speed resolution	0.01rpm	Power consumption	15W
Back suction angle	0-360°	Controller dimension	240*221*111mm
Display	7 inch- industrial grade- true color LCD TFT screen	Controller weight	2.1kg
Control	Touch screen	Memory function	Storage the running parameters when power off
Start/stop, direction signal	Active switch signal 5V	Pump housing material	Mirror stainless steel

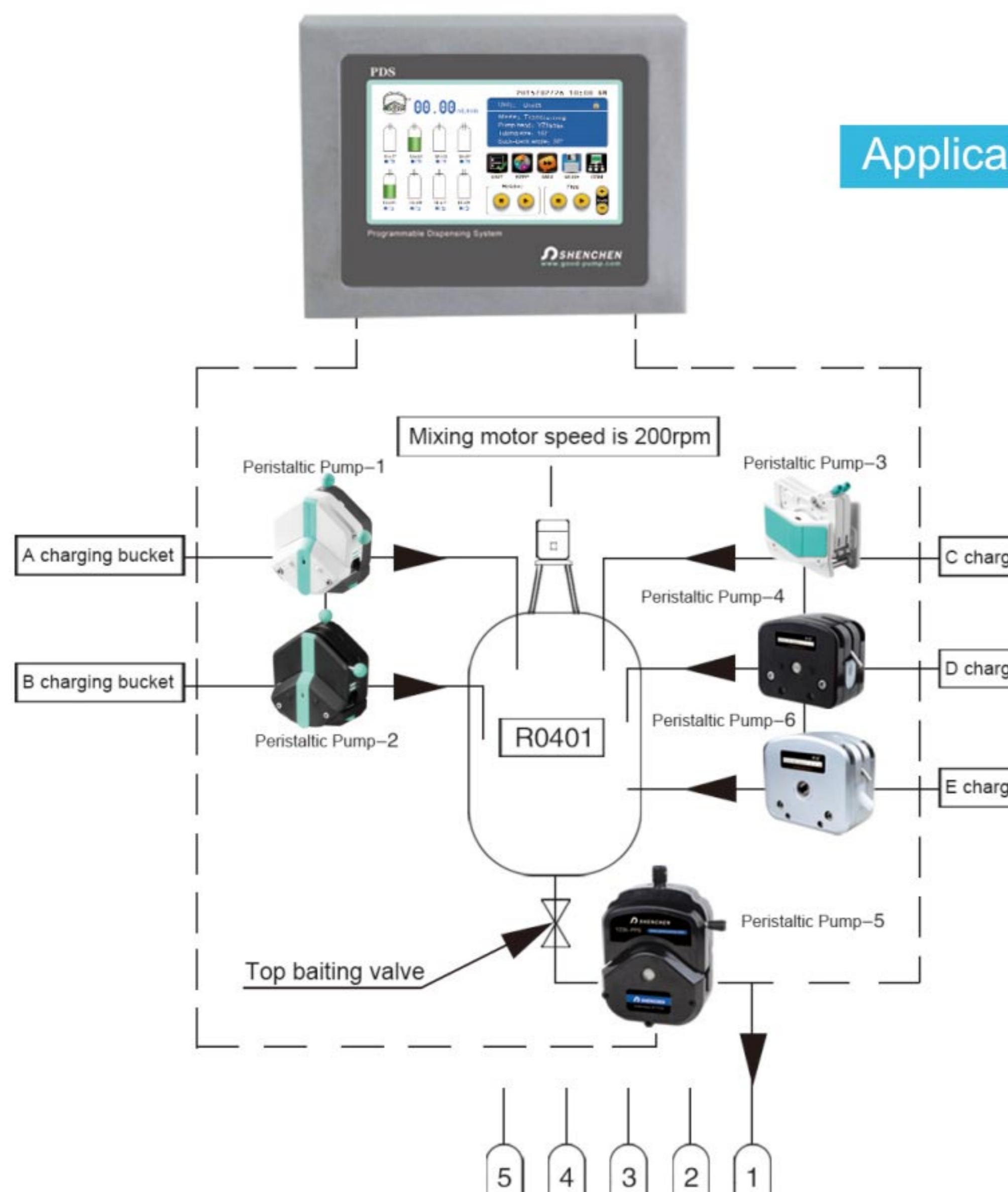

Application1




Application2: PDSIV, one controller with 7 pump units, 4 minipump and 3 YZ1515x pump heads, for different liquid filling with different time sequence.

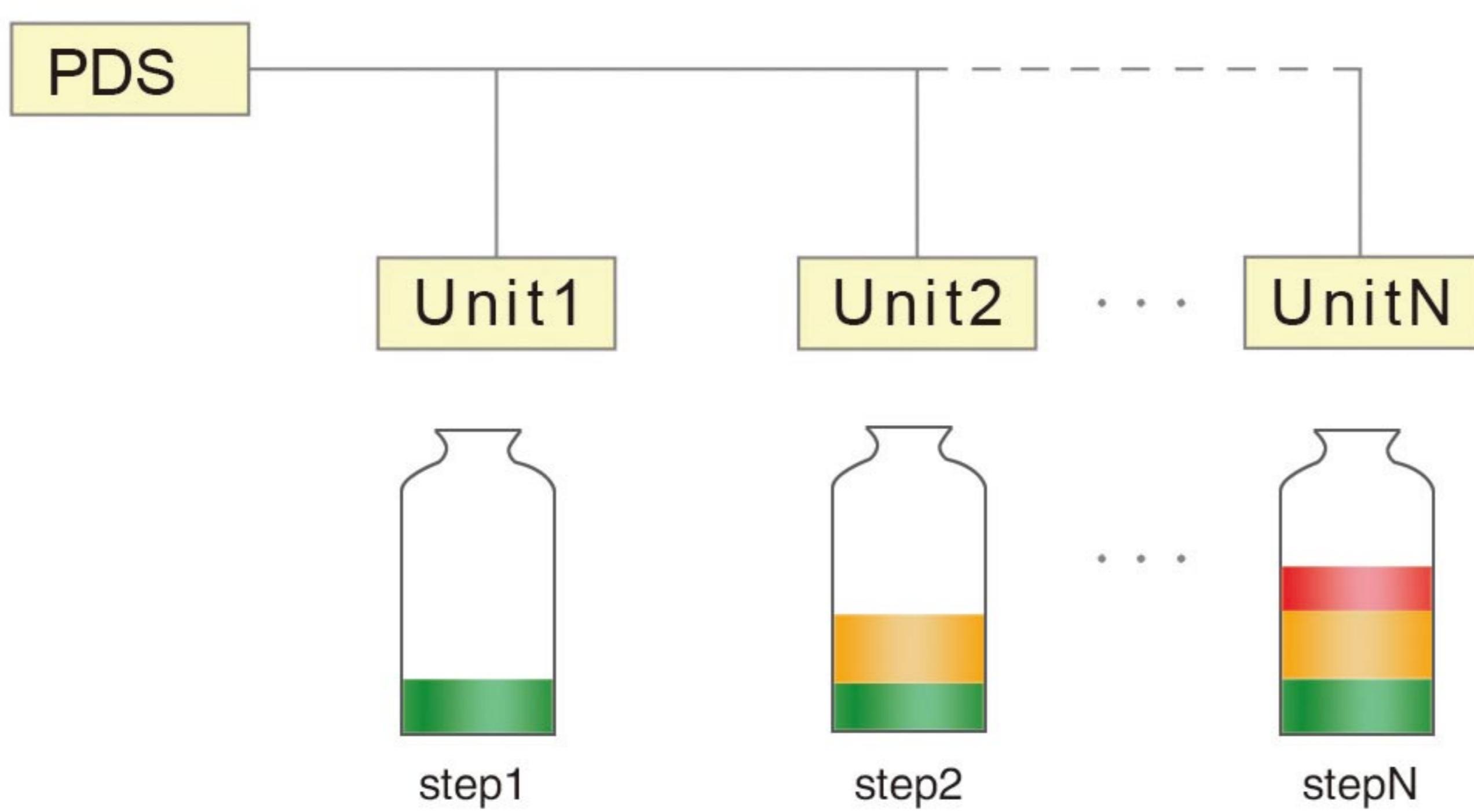


Logic working mode

In logic working mode, each unit working cooperatively. Filling or transferring with orders, also with different proportion adding, finish automation mixing work.
 Maximum steps: 50 steps
 Common mode: Can save maximum 5 groups data.
 Working time: 0.1sec ~9999hour
 Pause time: 0.1sec ~9999 hour
 Steps trigger way: Time trigger or external trigger.
 Calibration: Can calibrate each step separately, online micro adjusting function.

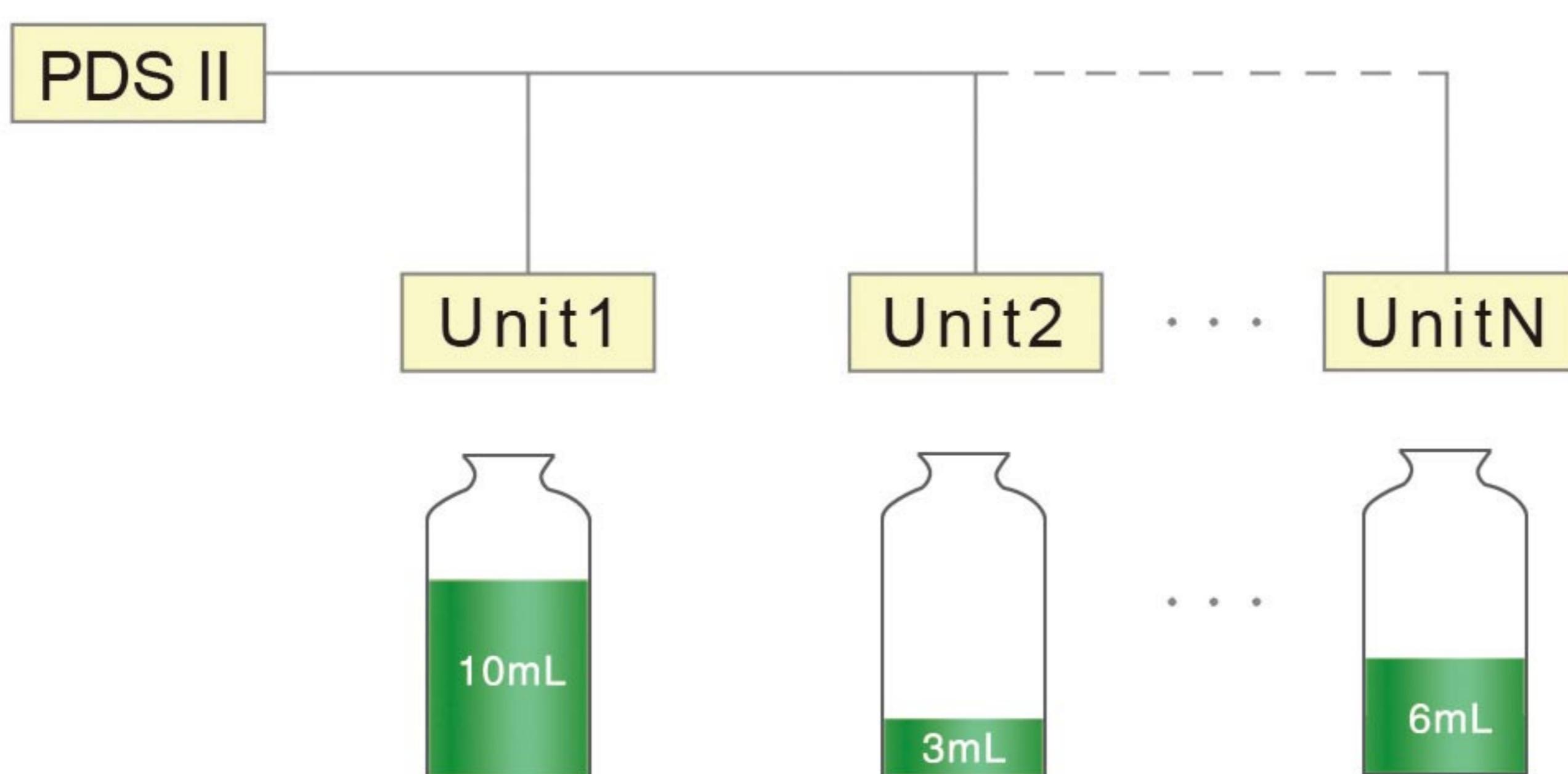

Independent working mode

In independent working mode, each unit working independently, can fill different volume.
 Working mode: transferring or dispensing
 Common mode: In dispensing mode can save 5 groups data.
 Calibration: Each pump unit can calibrate separately, online micro adjusting function.
 Dispensing time: 0.1~9999sec
 Pause time: 0.1~9999sec
 Repeat number: 1~9999 times, set '0' for unlimited times
 Communication: RS232/RS485, Modbus protocol (RTU mode)

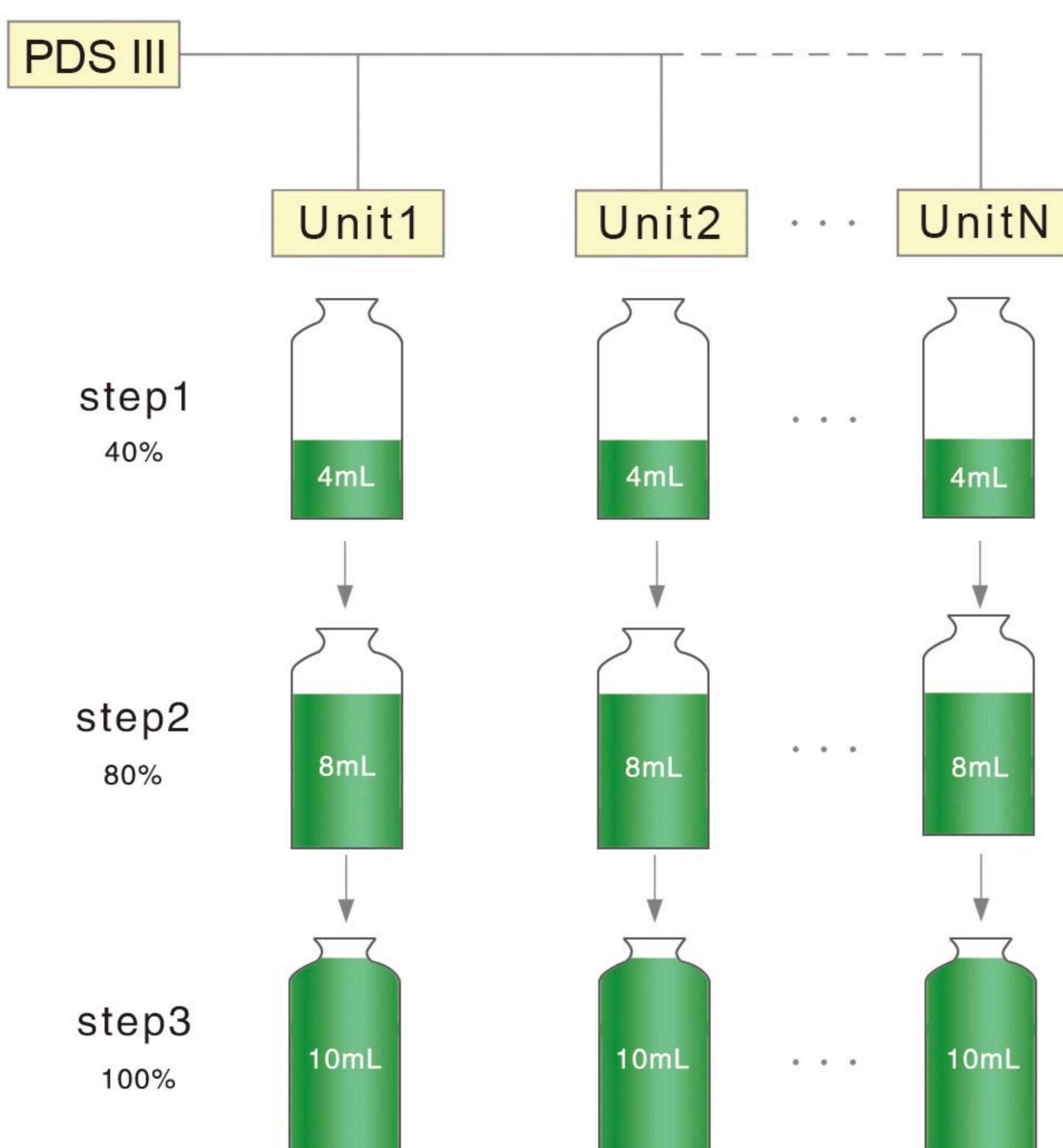

Independent and steps filling working mode

In independent working mode, each unit can finish filling process with multiple different filling volume. This function suitable for prevent fluid splashing and foaming. This process can also repeat many times, achieving complicated dispensing function.
 Working mode: Transferring and dispensing mode
 Transferring total volume: Can record the total liquid volume transferred by each unit in transfer mode.
 Dispensing step: It can be dispensed in three steps, and different parameters can be set for each step.
 Calibration: Can calibrate each step separately, online micro adjusting function.
 Common mode: In dispensing mode can save 5 groups data.
 Dispensing time: 0.1~9999.99sec
 Pause time: 0.1~9999.99 sec
 Repeat number: 1~9999times, set '0' for unlimited times.
 Communication: RS232/RS485, Modbus protocol (RTU mode)

PDS IV

Comprehensive working mode

Working mode: Include logic working mode, independent working mode, independent and steps filling mode.




Compact Peristaltic Pump

LabQ 3 years warranty



Dimension Drawing (Unit: mm)

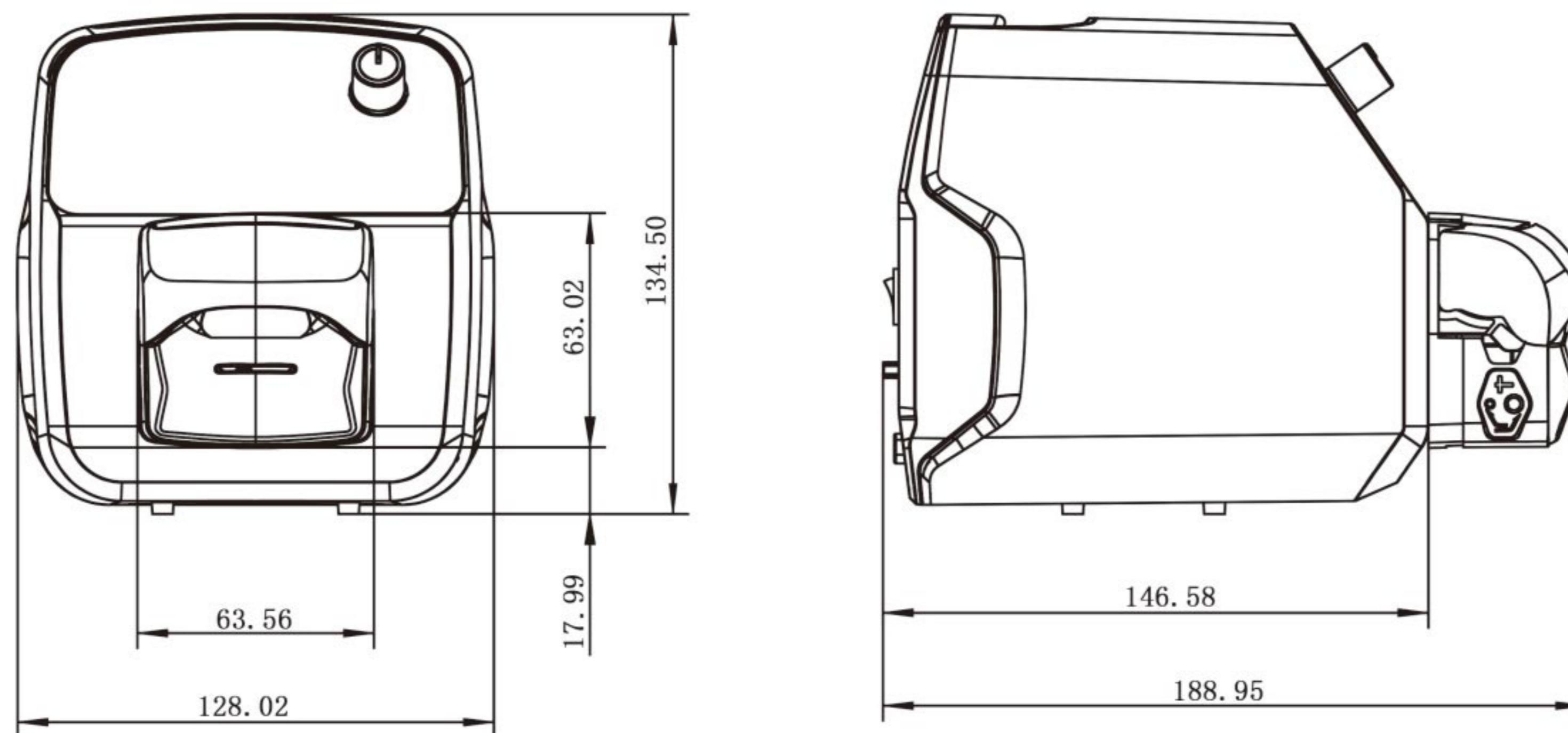

Technical Specifications

Speed range	0.1-350rpm	Start/stop, reversing signal	Switch signal(The default is passive signal, active signal is optional)
Speed resolution	0.1rpm	Communication interface	RS485, (MODBUS protocol, RTU mode)
Control method	Mechanical keypad and digital Knob	Dimension	188*128*135mm(L×W×H)
Display	2.4 inch LCD screen	Weight	1.1kg
External speed control signal	0-5V, 0-10V, 4-20mA	Power consumption	<30W
Output interface	Open-Collector output	Temperature	0-40°C
Power supply	Output: (24V/1.25A)	Relative humidity	<80%
	Input: AC100V-240V, 50Hz/60Hz		

Product Composition and Flow Rate Range

Model	Channel number	Tubing	ID×Wall thickness(mm)	mL / r	Speed(rpm)	Flow Rate(mL/min)	Weight(kg)
LabQ/KT15	Single channel	13#	0.8×1.6	0.033	0.1~350	0.0033~11.55	1.1
		14#	1.6×1.6	0.156		0.0156~54.60	
		19#	2.4×1.6	0.286		0.0286~100.10	
		16#	3.1×1.6	0.477		0.0477~166.95	
		25#	4.8×1.6	0.933		0.0933~326.55	





Compact Peristaltic Pump

3 years warranty

LabK



Dimension Drawing (Unit: mm)

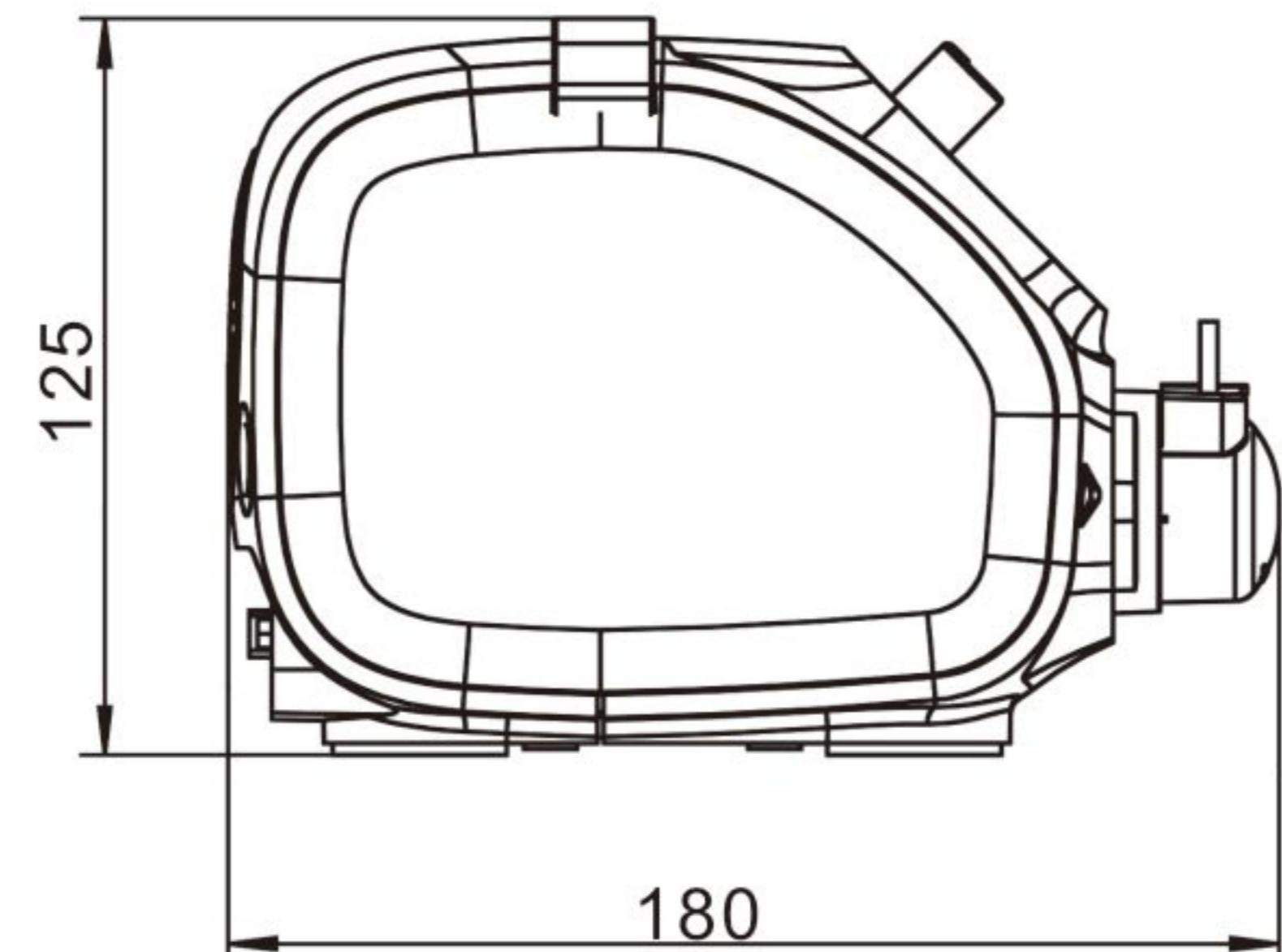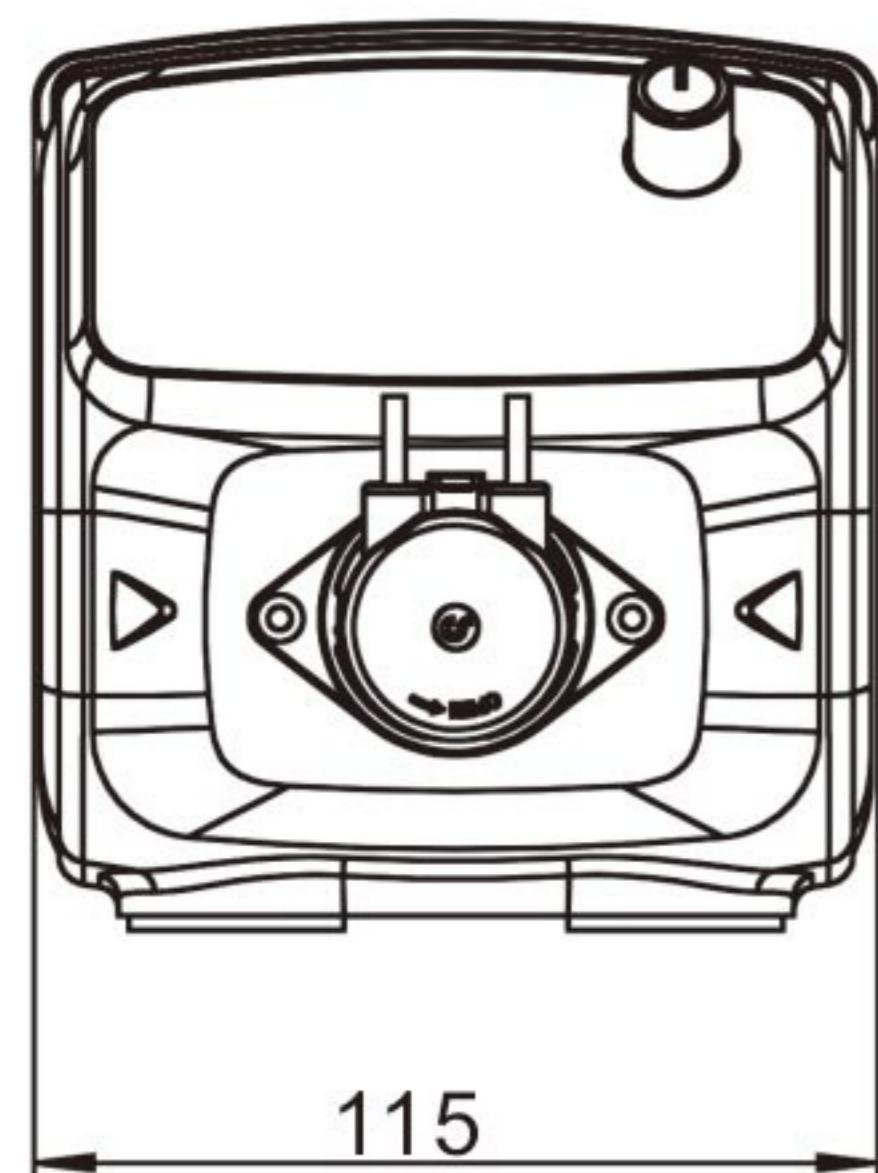

Color Selection



Technical Specifications

Flow rate range	0.004-63.96mL/min	External control	Start/stop direction control (switch signal), 0-5V, 4-20mA (standard), 0-10V (optional)
Speed range	0.1-150rpm reversible	Power adapter	Output: (12V/1A); Input: AC100V-240V, 50Hz/60Hz
Speed resolution	0.1rpm	Dimension	180*115*125mm(L×W×H)
Speed control	Digital knob	Weight	0.8kg
Control method	Digital knob control and mechanical keypad	Power consumption	<12W
Keypad lifetime	300,000 times	Condition temperature	0-40°C
Display	0.96" OLED display	Relative humidity	<80%
Communication interface	USB connector, RS485 interface (MODBUS protocol, RTU mode)		

Product Composition and Flow Rate Range


Model	Pump Head	Speed(rpm)	Tubing Size (ID×Wall Thickness(mm))	Flow Rate Range (mL/min)
LabK	MicroPump	0.1-150	1*1	0.004-6.38
			2*1	0.014-21.45
			3*1	0.031-47.26
			4*1	0.042-63.96



Compact Peristaltic Pump



Dimension Drawing (Unit: mm)

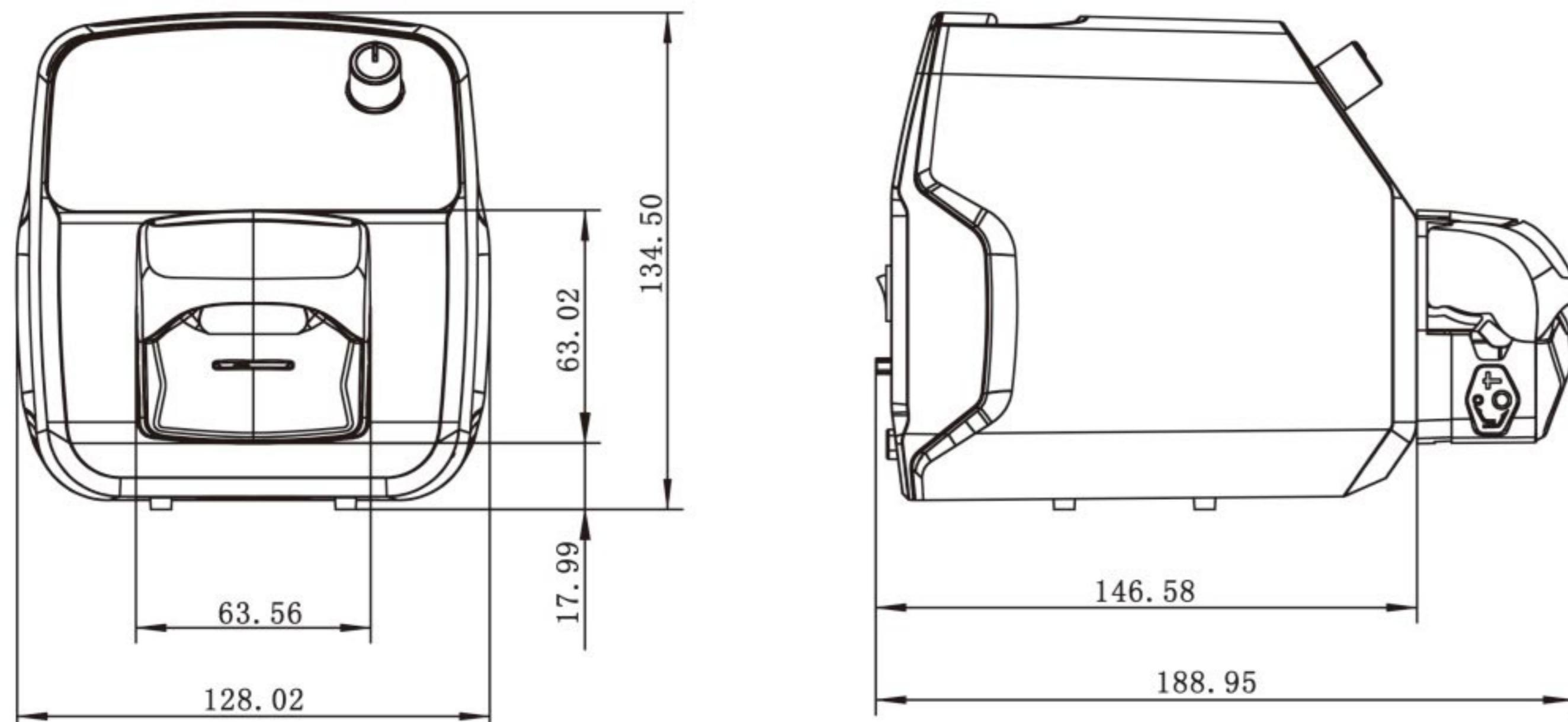

Product Introduction

- ABS engineering plastic housing, anti-corrosion, anti-static.
- OLED screen display motor speed, digital knob control speed.
- Compact design, various external control.
- Easy to observe the pump head and tubing working situation.

Product Features

- Low power consumption, mute working.
- Stable flow rate and suitable for continuous dosing applications.
- Easy to replace long life PharMed tubing.
- Digital knob control speed, memory back up, user setting saved if power lost.

Pump Head



LabT/UD15

LabT/Minipump

Technical Specifications

Model	LabT/KT15	Display	0.96" OLED display
	LabT/MiniPump	Power adapter	Output: 24V/1.25A; Input: AC100V-240V, 50Hz/60Hz
	LabT/UD15		
Speed range	LabT/KT15: 0.1-350rpm	Communication interface	USB connector, RS485 interface (MODBUS protocol, RTU mode)
	LabT/MiniPump: 0.1-300rpm		
	LabT/UD15: 0.1-350rpm	Dimension	189*128*135mm (L×W×H)
Speed resolution	0.1rpm	Weight	1100g
Control method	Digital knob control and mechanical keypad	Power consumption	< 15W
External control	Start/stop direction control (switch signal) , 0-5V, 4-20mA (standard), 0-10V (optional)	Condition temperature	0-40°C
		Relative humidity	<80%

Product Composition and Flow Rate Range

Model	Pump Head	Speed(rpm)	Tubing Size	Flow Rate Range(mL/min)
LabT	KT15	0.1-350	13"	0.0033~11.55
			14"	0.0156~54.6
			19"	0.0286~100.10
			16"	0.0477~166.95
			25"	0.0933~326.55
	MiniPump01	0.1-300	13", 14", 19", 16", 25", 1x1, 2x1, 2.5x1, 3x1	0.0024~190
			1x1, 2x1, 2.5x1, 3x1	0.005~108.39
	UD15	0.1-350	16", 17", 25"	0.08~930



Compact Peristaltic Pump

3 years warranty
ST-HandyPump




Product Features

- 304 stainless steel shell, resist corrosion. Can be stacked to use.
- Closed-loop stepper motor, compact structure, low noise, low vibration.
- OLED screen display working speed and status. Digital knob control speed, easy operation.
- Support RS232 and RS485 communication.
- Various external control interface.
- Wide voltage design, high versatility.

Model number | ST-HandyPump

Dimension Drawing (Unit: mm)

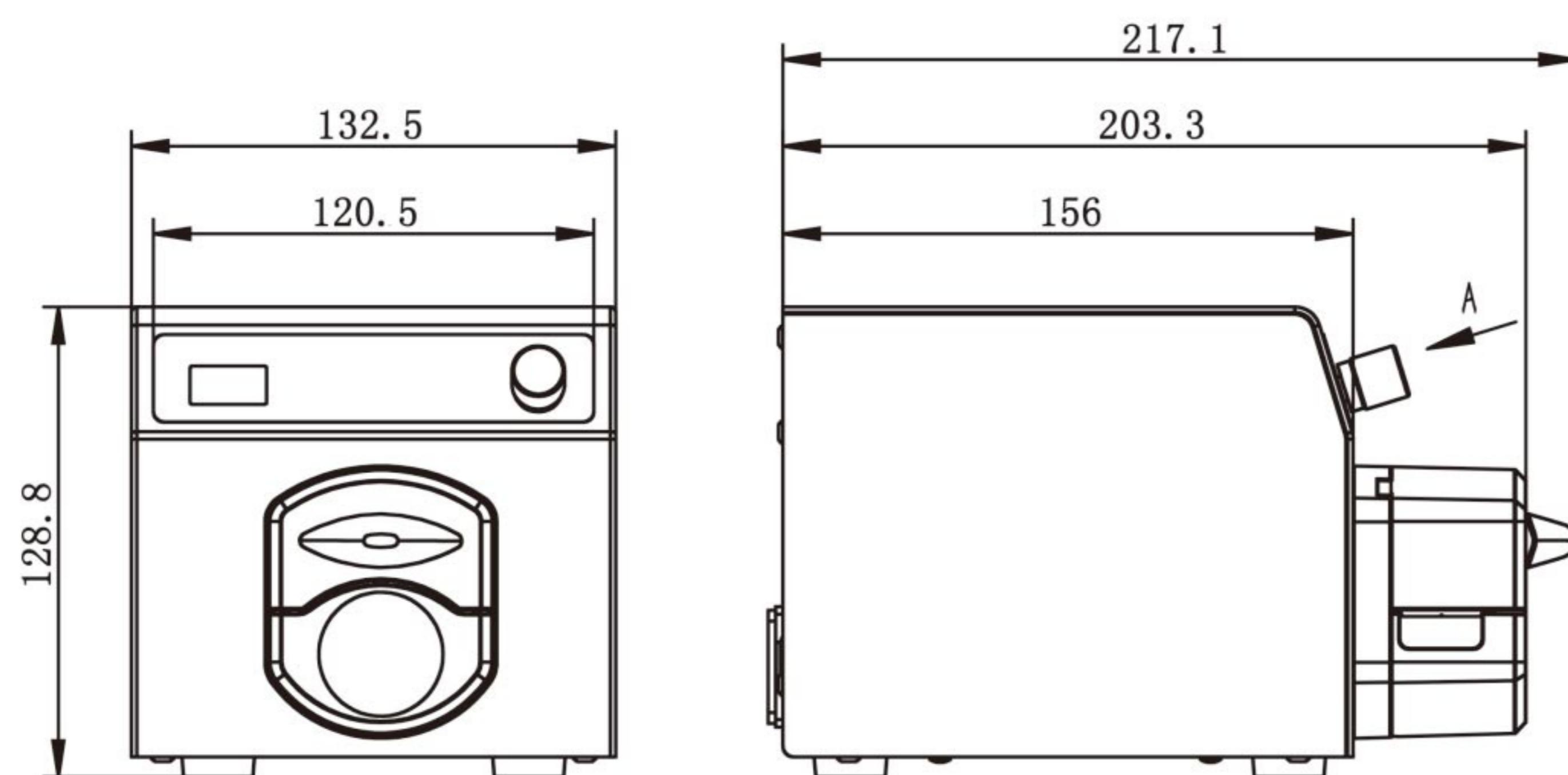

Pump Head



HandyPump01



HandyPump02

Technical Specifications

Flow rate range	0.0033~365.69mL/min	Power supply	Output: (24V/1.25A); Input: AC100V-240V, 50Hz/60Hz
Speed range	0.1~300rpm	Communication interface	RS485, RS232
Speed resolution	0.1rpm	Output interface	Output motor working status
Control method	Digital knob control and mechanical keypad	Drive dimension	217.1*132.5*128.5mm(L*W*H)
Motor type	57 closed-loop stepper motor	Weight	2.92Kg
Keypad lifetime	300,000 times	Power consumption	≤75W
Display	OLED display(0.96")	Condition temperature	0-40°C
External control	Start/stop control(switch signal) Speed: 0-5V, 4-20mA(standard), 0-10V (optional)	Relative humidity	<80%

Product Composition and Flow Rate Range

Model	Pump Head	Speed(rpm)	Tubing Size	Flow Rate Range(mL/min)
ST-HandyPump	HandyPump01	0.1-300	13#, 14#, 19#, 16#, 25#	0.0033-365.69
	HandyPump02		13#, 14#, 19#, 16#	0.0033-190.76






Compact Peristaltic Pump

3 years warranty

SP-MiniPump



Dimension Drawing (Unit: mm)

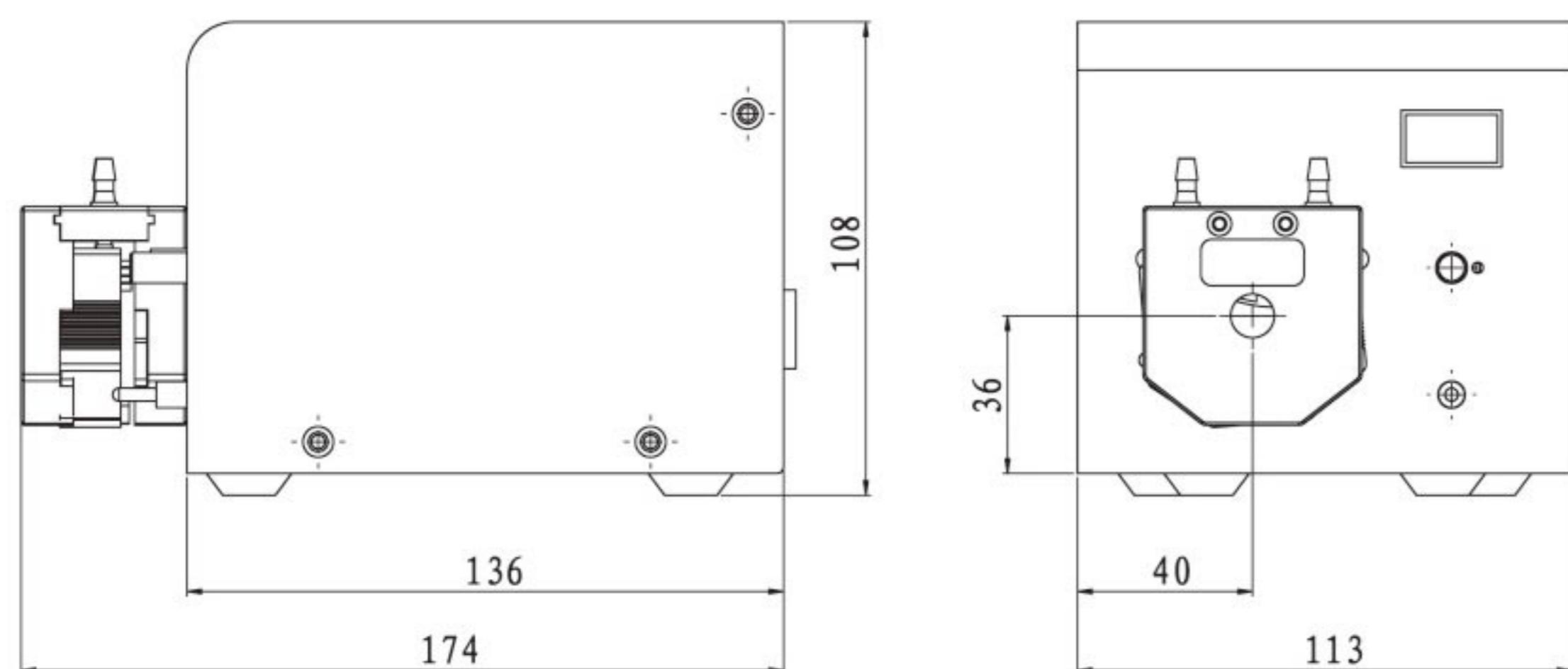

Features

- 304 stainless steel housing, resisting corrosion, space-saved.
- OLED screen, display the current motor speed and working status.
- Digital knob control speed, toggle switch control direction and start/stop.
- Various external control functions, support RS485 standard MODBUS protocol.

Model Number

SP-MiniPump01, SP-MiniPump02

External Control Schematic Diagram



Technical Specifications

Flow rate range	0.0024~190 mL/min	Communication interface	RS232, RS485 support MODBUS protocol (RTU mode)
Speed range	0.1~300 rpm	Power supply	24V/1.25A DC power
Speed resolution	0.1rpm	Drive dimension	136×113×108 mm (L×W×H)
Speed control	Digital knob control speed	Weight	1.64kg
Motor type	Stepper motor	Power consumption	< 30W
Display	OLED display	Condition temperature	0-40°C
External control	Start/Stop control (switch signal)	Relative humidity	< 80%
	0-5V, 4-20mA(standard), 0-10V(optional)		

SP-MiniPump Speed/Flow Rate Reference

Model	Pump Head	Speed(rpm)	Tubing Size	Per Channel Flow Rate Range(mL/min)
SP-MiniPump	MiniPump01	0.1-300	13"	0.0024-8.28
			14"	0.0112-33.88
			19"	0.0252-77.23
			16"	0.0394-114.31
			25"	0.0652-190.00
	MiniPump01/ MiniPump02	0.1-300	1*1	0.005-15.01
			2*1	0.018-54.63
			2.5*1	0.0256-76.84
			3*1	0.0356-108.39

Flow Rates Peristaltic Pump

LabV/V Series

3 years warranty



- | Laboratory
- | Industrial equipment supporting
- Industrial grade 4.3" true color LCD screen, touch screen control.
- Dynamic display transferring status. Flow rate data, setting parameters and system configuration display in the same screen.
- 3 Kinds of working mode: fixed volume metering, fixed time and volume, timer start and stop, meet different transferring and dispensing request.
- Intelligent calibration function and online micro adjusting function.

Suitable Pump Head

Flow Rates Peristaltic Pump

V1, V3, V6

3 years warranty



EasyPump Series
(Pressure Adjustable)



EasyPump Series
(Fixed Pressure)




AMC Series



YZ Series



External Control Schematic Diagram

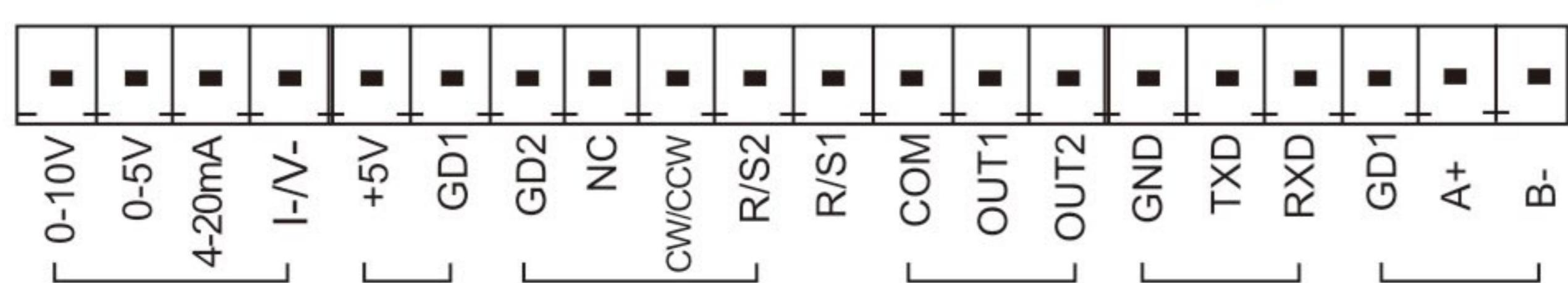

V1, V3, V6

Laboratory
Industrial equipment




Technical Specifications

Flow rate range	LabV1-III: 0.0053~775 mL/min	External speed control signal	0-5V, 0-10V, 4-20mA for option
	LabV3-III: 0.0053~1808 mL/min		
	LabV6-III: 0.0053~3100 mL/min	Start/stop, direction signal	Passive switch signal, such as foot pedal switch
	LabV1/V1: 0.0002~570 mL/min		Active switch signal: 5-24V universal
	LabV3/V3: 0.0002~1330 mL/min	Communication interface	RS232, RS485 support MODBUS protocol (RTU mode)
	LabV6/V6: 0.0002~2280 mL/min		
Speed range	0.1-600 rpm	Output interface	Output motor working status (Open-Collector output)
Speed resolution	0.01 rpm		
Flow rate resolution	0.01 mL	Power supply	AC 220V±10% 50Hz/60Hz (standard)
Flow rate accuracy	<±0.5%		AC 110V±10% 50Hz/60Hz (optional)
Back suction angle	0-360°	Drive dimension (L×W×H)	LabV Series: 261.4×157.3×236.9 mm
Motor type	LabV-III series: Closed-loop stepper motor LabV/V series: Stepper motor		V Series: 252×152×243 mm
Display	Industrial grade 4.3" LCD color display	Drive weight	LabV Series: 4.40 kg
Control method	Touch screen and Mechanical keypad		V Series: 4.20 kg
Keypad lifetime	300,000 times	Power consumption	<80W
		Condition temperature	0-40°C
		Relative humidity	< 80%

Product Composition and Flow Rate Range

Flow Rates Peristaltic Pump		Pump Head & Flow Rate (mL/min)		
		New Generation Quick Load Type Pump Head		
Drive&speed	Tubing	EasyPumpI/III 13#, 14#, 19#, 16#, 25#, 17#, 18#	EasyPumpII/IV 15#, 24#, 35#, 36#	EasyPumpV/VI(dual channel) 13#, 14#, 19#, 16#, 25#
LabV1-III	0.1-150 rpm	0.0053~645	0.18~775	0.0053~295
LabV3-III	0.1-350 rpm	0.0053~1505	0.18~1808	0.0053~688
LabV6-III	0.1-600 rpm	0.0053~2580	0.18~3100	0.0053~1180
Drive&speed	Tubing	YZ1515x 13#, 14#, 19#, 16# 25#, 17#, 18#	YZ2515x 15#, 24#	AMC1~AMC12(10) Inner diameter: 0.13-3mm Wall thickness: 0.8-1mm
LabV1/V1	0.1-150 rpm	0.007~570	0.17~435	0.0002-48(working speed≤150rpm)
LabV3/V3	0.1-350 rpm	0.007~1330	0.17~1015	0.0002-65(working speed≤150rpm)
LabV6/V6	0.1-600 rpm	0.007~2280	0.17~1740	

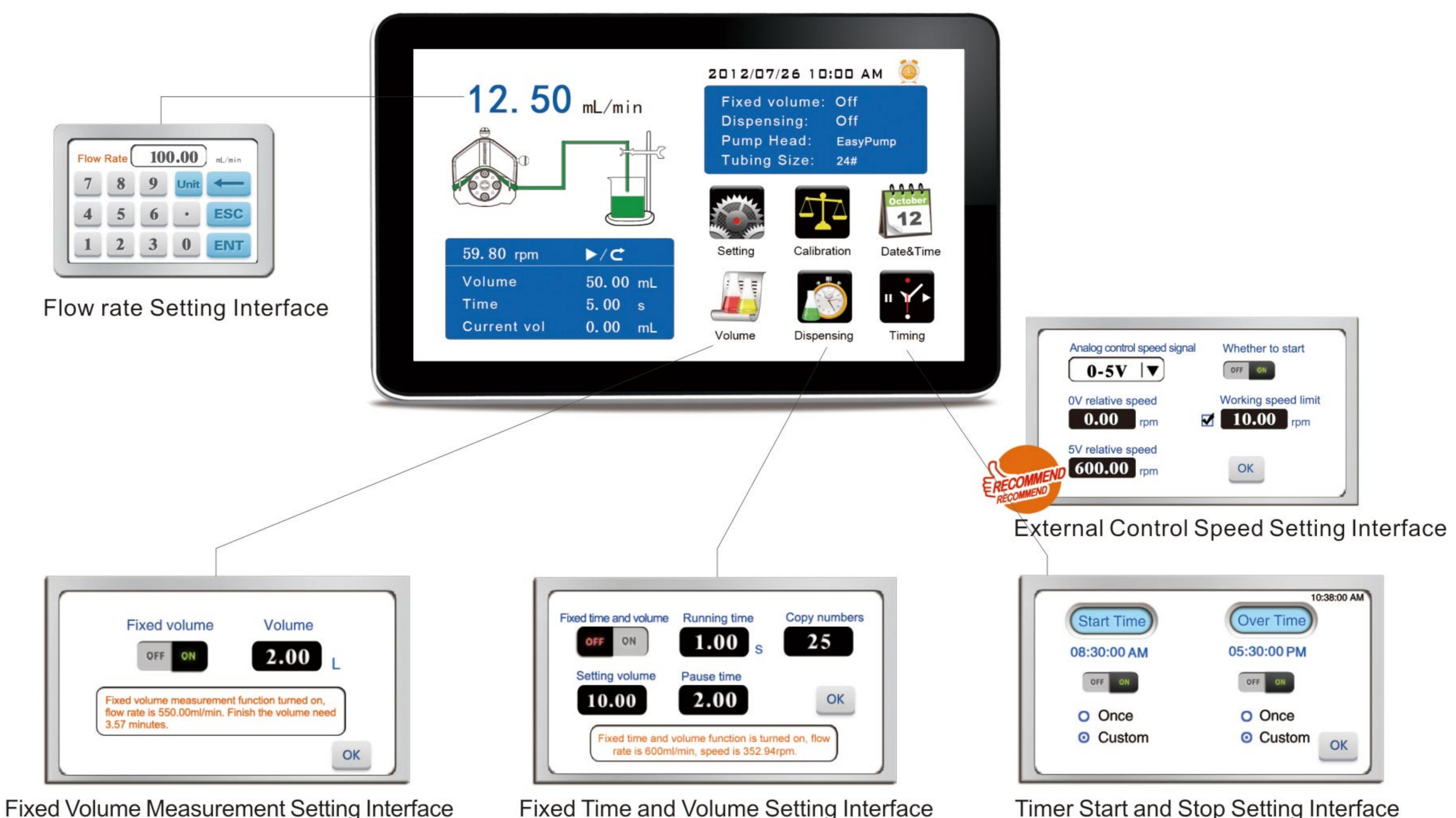

Flow Rates Peristaltic Pump

V6-12L



V6-6L



V6-3L



V6-3L/EasyPump

3 years warranty

304
SS Housing

Features

Large flow rate, high precision, intelligent control of liquid transferring.

Closed-loop stepper motor drive, accurate control, strong driving force.

304 stainless steel housing, the first choice for high level industrial sites.

Model Number

- V6-3L/EasyPump
- V6-3L/DZ25-3L
- V6-6L/DZ25-6L
- V6-12L/YZ35

Technical Specifications

Flow rate range	V6-3L: 0.0053~3600 mL/min V6-6L: 0.3~6000 mL/min V6-12L: 0.69~12000 mL/min	Start/stop,direction signal	Passive switch signal, such as foot pedal switch Active switch signal: 5-24V universal
Speed range	0.1-600 rpm	Communication interface	RS232, RS485 support MODBUS protocol (RTU mode)
Speed resolution	0.01 rpm	Output interface	Output motor working status (Open-Collector output)
Flow rate resolution	0.01 mL	Power supply	AC 220V±10% 50Hz/60Hz (standard) AC 110V±10% 50Hz/60Hz (optional)
Flow rate accuracy	<±0.5%	Drive dimension	V6-3L: 223×152×231mm V6-6L: 283×192×274mm V6-12L: 302×222×331mm
Back suction angle	0-360°	Drive weight	V6-3L: 5.02kg; V6-6L: 7.85kg; V6-12L: 13.14kg
Motor type	Closed-loop stepper motor	Power consumption	V6-3L: <80W ; V6-6L: <180W ; V6-12L: <300W
Display	Industrial grade 4.3" color LCD display	Environment temperature	0-40°C
Control method	Touch screen and Mechanical keypad	Relative humidity	< 80%
Keypad lifetime	300,000 times		
External speed control signal	0-5V,0-10V,4-20mA		

Product Composition and Flow Rate Range

Drive	Motor Type	Pump Head	Tubing Size	Speed Range(rpm)	Flow Rate(mL/min)
V6-3L	Closed-loop stepper motor	EasyPump	13#, 14#, 19#, 16#, 25#, 17# 18#/ 15#, 24#, 35#, 36#	0.1-600	0.0053~3100
		DZ25-3L	15#, 24#, 35#, 36#		0.211~3600
V6-6L		DZ25-6L	15#, 24#, 35#, 36#		0.3~6000
V6-12L		YZ35	26#, 73#, 82#		0. 69~12000



V Series peristaltic pump has various external control interface, can meet different equipment supporting requirements.

| External control interface definition

1. Analog signal input port: Choose the external control speed signal from the external control setting interface. Turn on the external control speed function, control the motor speed range through the analog signal.

2. Internal isolation 5VDC output.

3. External control start/stop,direction signal input port:


active signal input.

4. R/S1 external control start/stop signal input port: passive signal input.

5. Motor running status output port: output current running status of the motor.

6. RS232 communication interface: choose RS232 in communication setting interface, this port is effective.

7. RS485 communication interface: choose RS485 in communication setting interface, this port is effective.

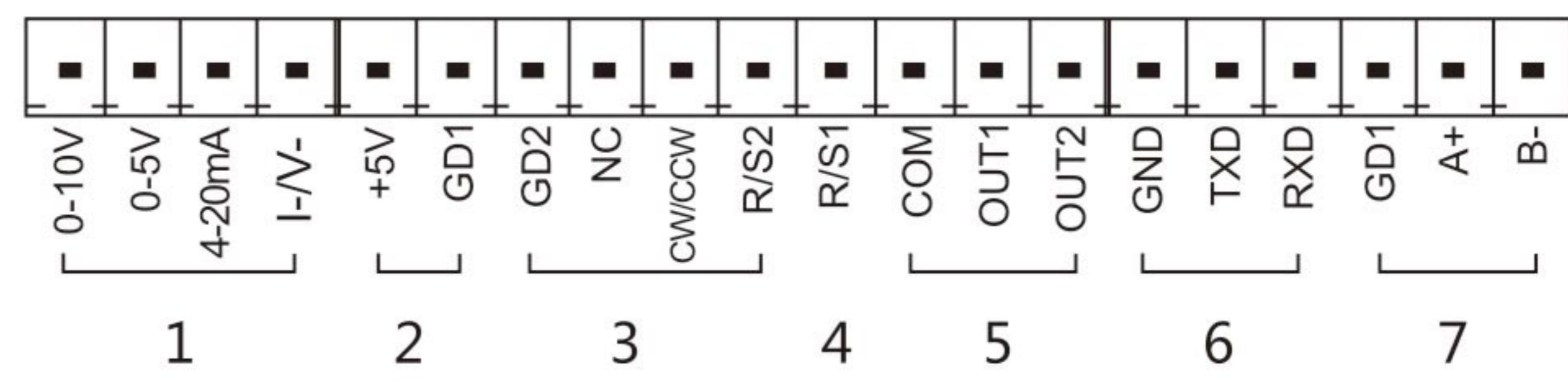

V Series peristaltic pump external control setting interface.


| Choose external control start/stop,direction signal

1. Can respectively set external control start/stop or reversing switch whether effective or not.

2. Can choose control mode according to requirement: pulse mode or level mode.

3. Can set the foot switch whether effective or not.

4. Can choose switch value signal 's high level values according to customers' main control equipment: 5-24V universal.

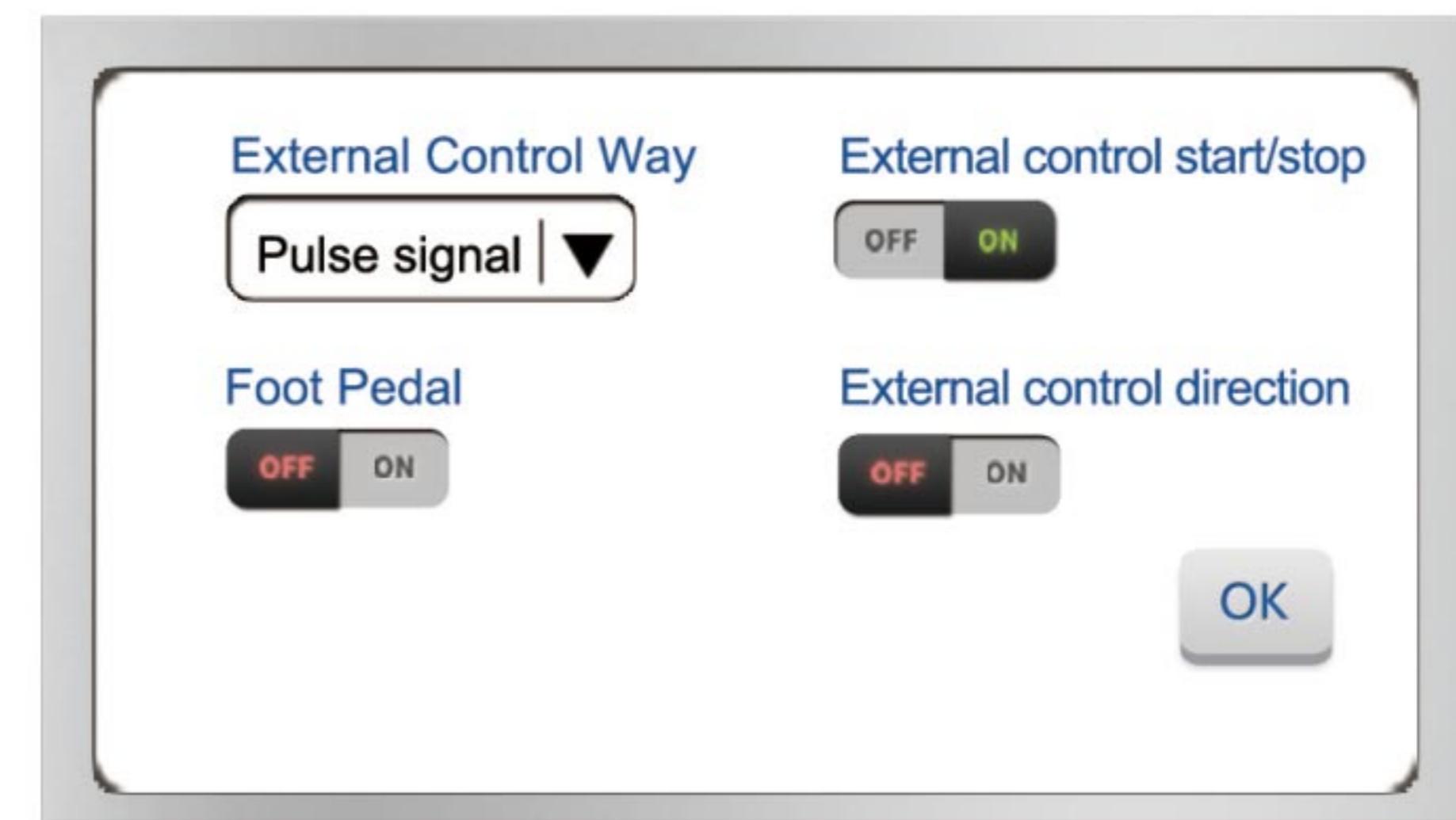


| Programmable external control speed setting

1. Can choose analog signal source according to requirements: 0-5V, 0-10V, 4-20mA.

2. When the signal source is chosen, can choose maximum speed and minimum speed which corresponding to the signal source's maximum value and minimum value according to actual demand, to reach userdefined rotate speed range purpose.

3. Working speed limit is at the situation that the linear relation of analog quantity signal and rotating speed keeps invariant, set peristaltic pump's maximum working rotate speed. This setting can avoid production accident caused by sudden change of transfer fluid amount with the sudden change of external analog signal.

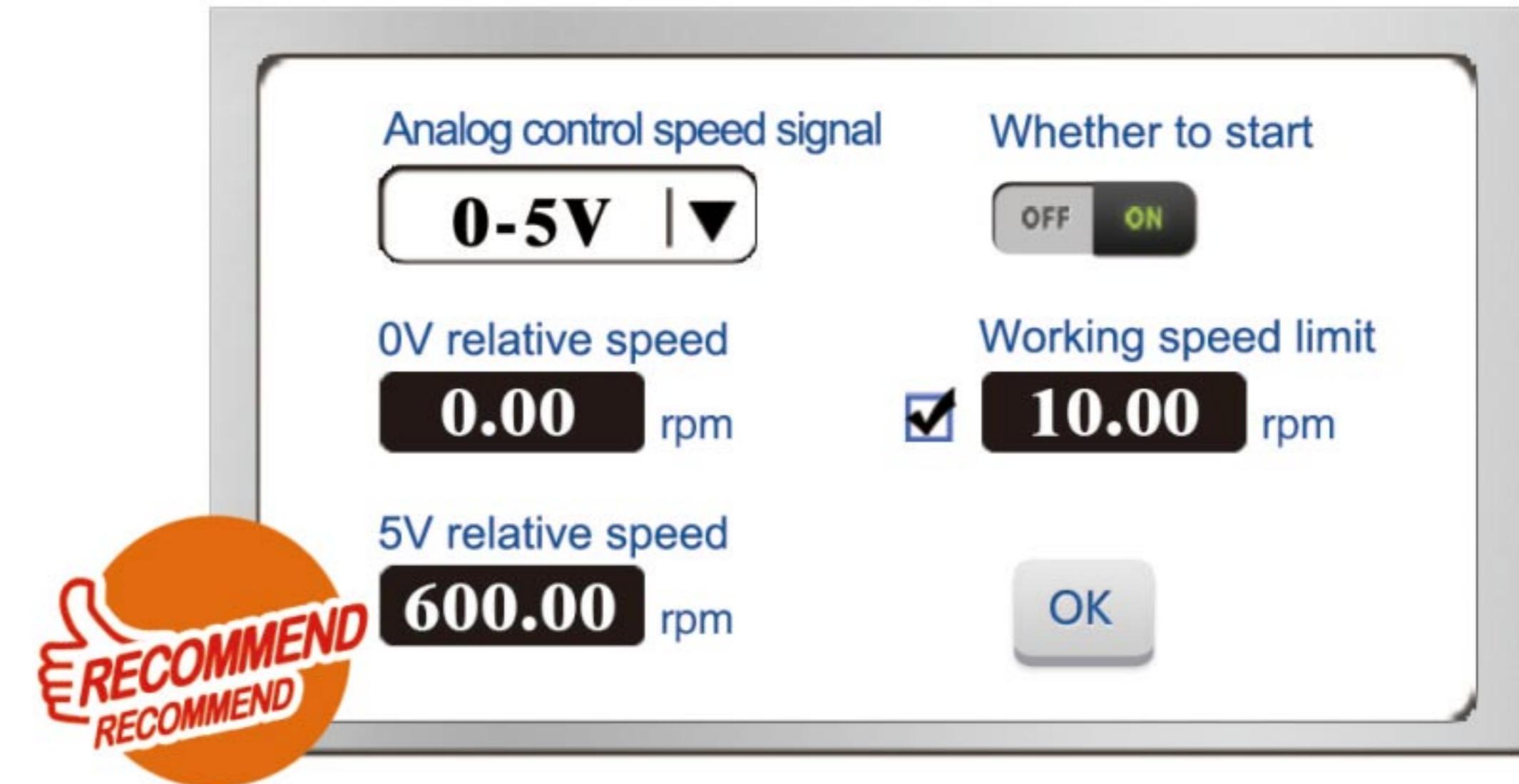


| Communication setting


1. Support RS232 and RS485 interface, can be chosen in software interface.

2. Support various Baud rate : 2400bps, 4800bps, 9600bps, 19200bps.


3. Can choose standard modbus communication protocol or Shenzhen communication protocol to control, modbus communication protocol is suitable to match with industrial site HMI, PLC or other upper computer which support standard modbus communication protocol; Shenzhen communication protocol is suitable to match with single chip micro computer or upper computer on computer programme.

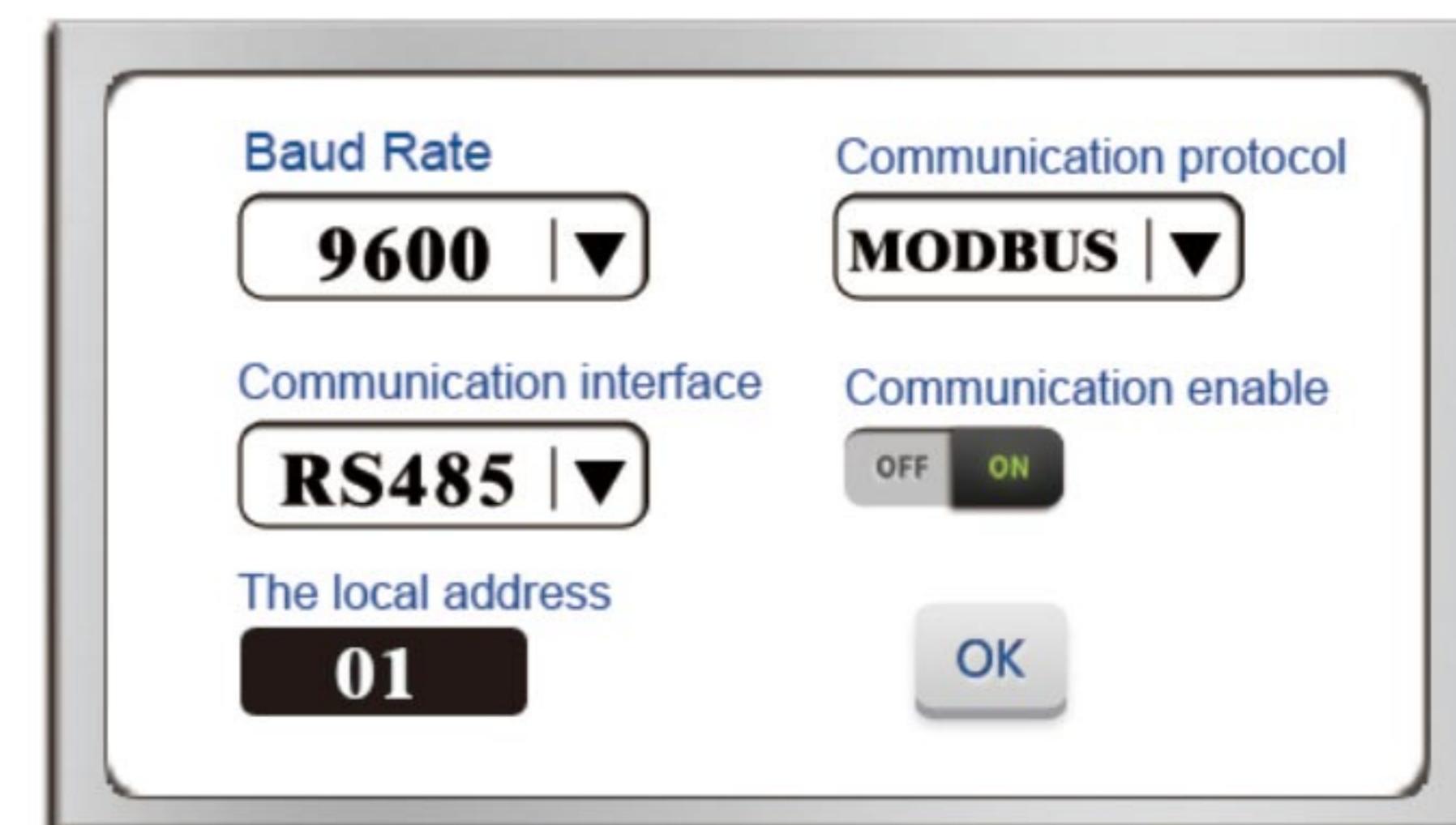


Dimension Drawing(Unit: mm)

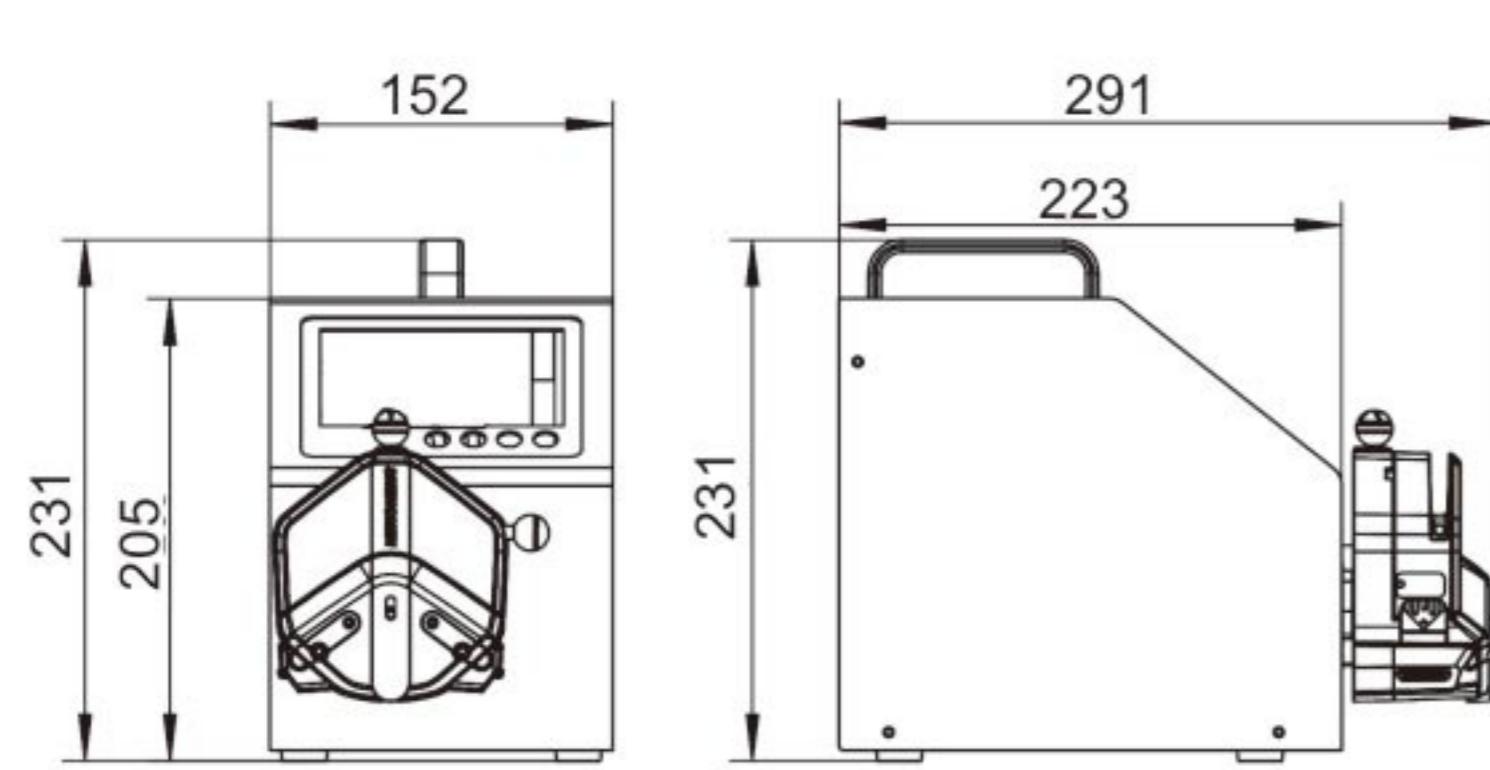

V6-3L/EasyPump

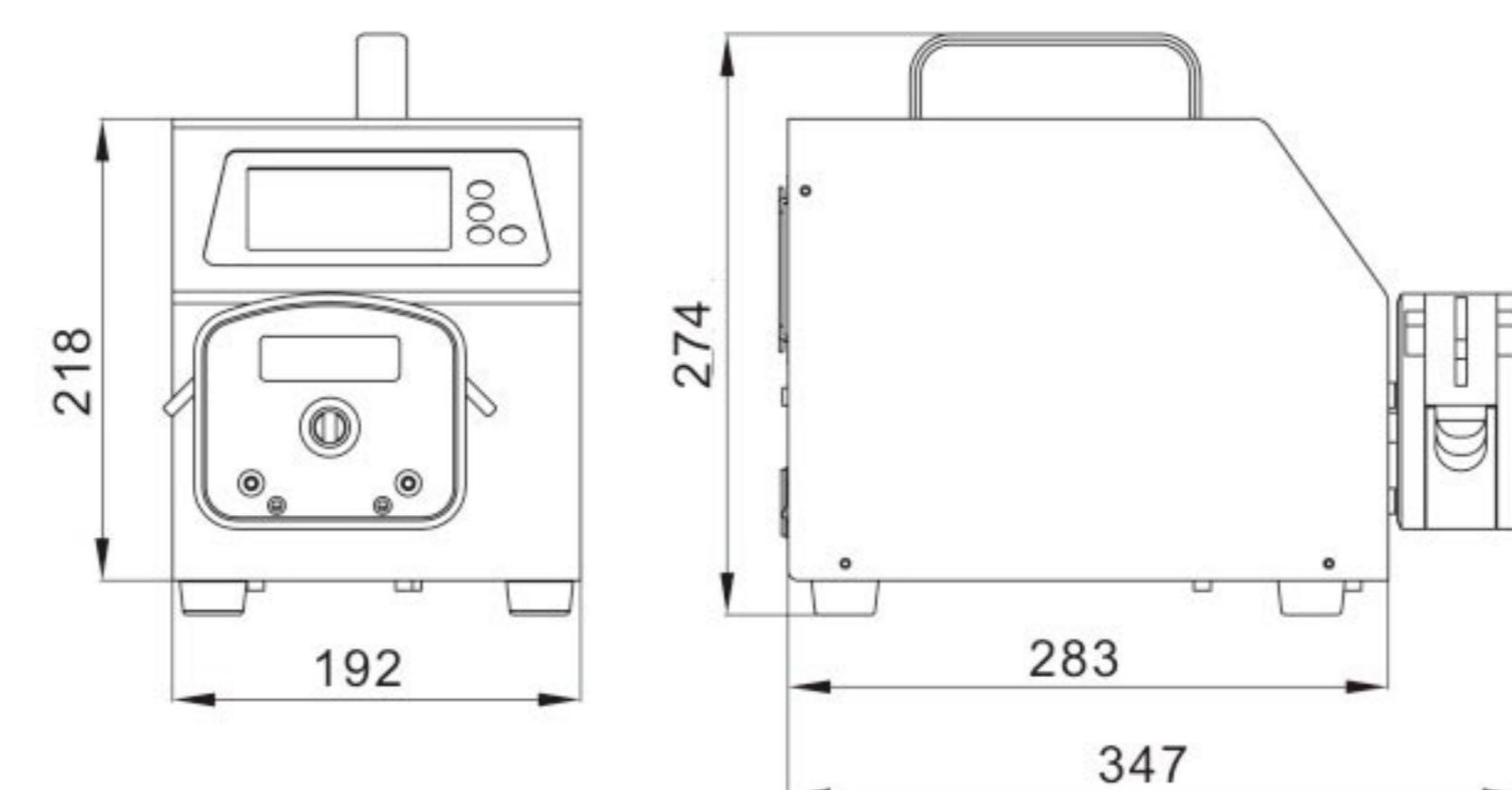

V6-6L

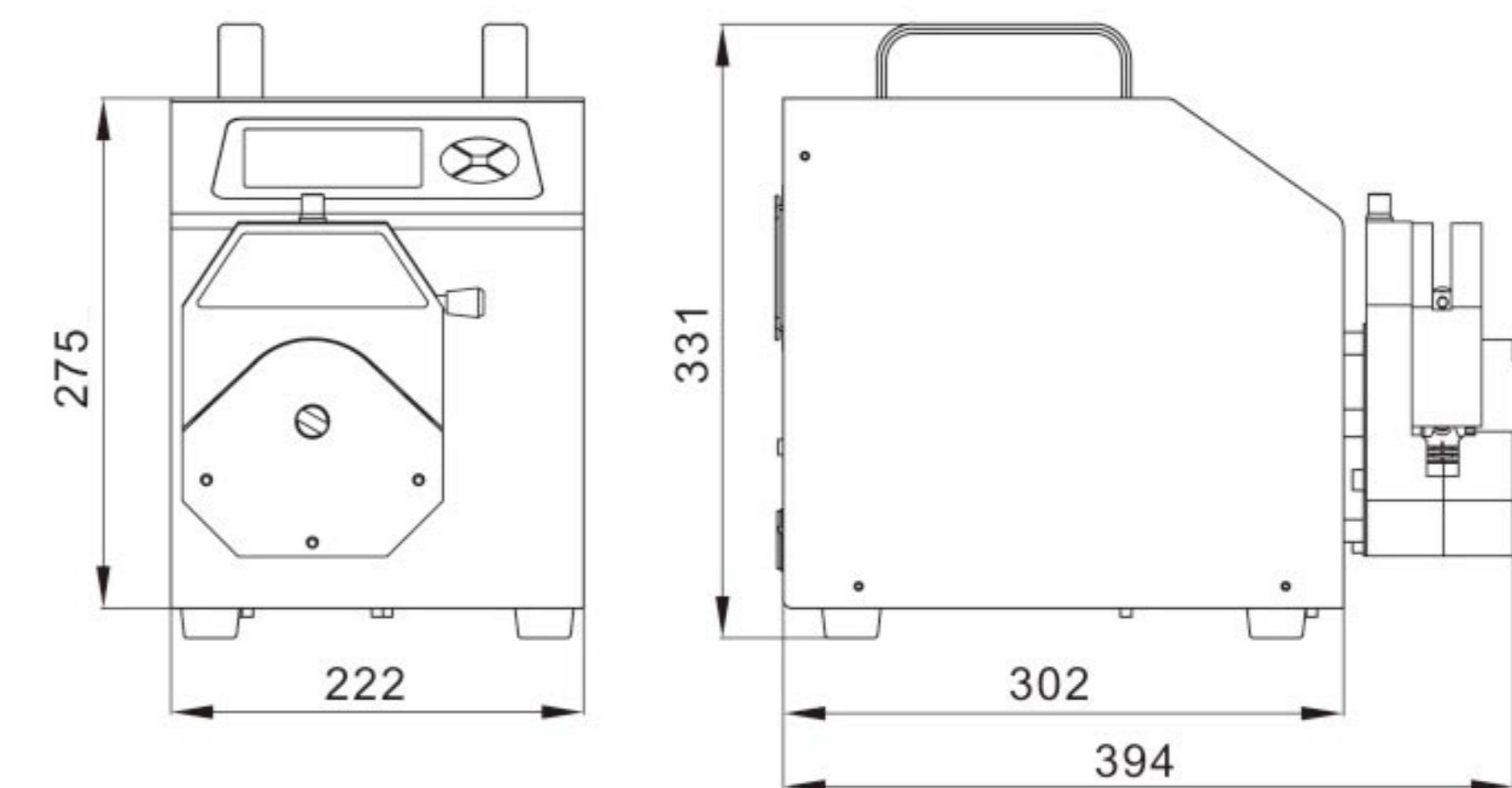

V6-12L



Flow Rates Peristaltic Pump

LabN Series

3 years warranty



Suitable Pump Head



EasyPump Series
(Pressure Adjustable)



EasyPump Series
(Fixed Pressure)



AMC Series




YZ Series

Typical Application

- Special for university laboratory and research institute.
- Ion chromatography and titrator
- Pilot scale and industry production

LabN Series External Control Schematic Diagram

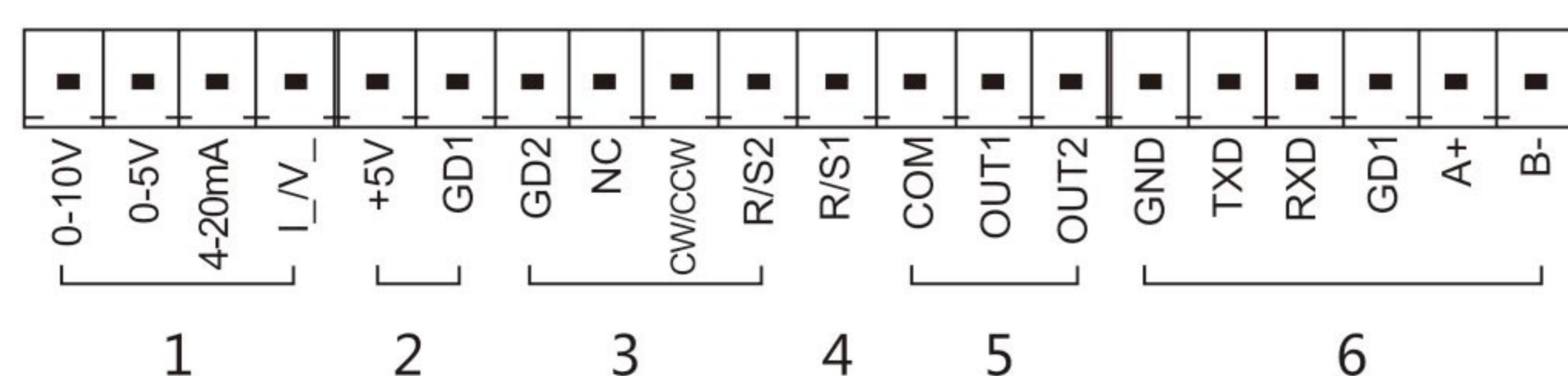

Features


- | 3.2" color LCD screen display.
- | Flow rate and motor speed display in same screen.
- | Timing function, time range 0.1s-9999 hours, can be used for simple dispensing function.

Technical Specifications

Flow rate range	LabN1-III: 0.0053~775 mL/min	Start/stop, reversing signal	Passive switch signal, such as foot pedal
	LabN3-III: 0.0053~1808 mL/min		Active switch signal: 5V-24V universal
	LabN6-III: 0.0053~3100 mL/min	Communication interface	RS232, RS485 support Modbus protocol(RTU mode)
	LabN1: 0.0002~570 mL/min		
	LabN3: 0.0002~1330 mL/min	Output interface	
	LabN6: 0.0002~2280 mL/min		Output motor working status (Open-Collector output)
Speed resolution	0.1rpm	Power supply	AC 220V±10% 50Hz/60Hz (standard)
Back suction angle	0-360°		AC 110V±10% 50Hz/60Hz (optional)
Testing time range	0.1 s-9999 h	Drive dimension	
Motor type	LabN-III series: Closed-loop stepper motor		261.4×157.3×236.9mm
	LabN series: Stepper motor		
Display	3.2" high definition LCD screen	Drive weight	4.40 kg
Control method	Digital knob and Mechanical keypad		
Keypad lifetime	300,000 times	Power consumption	<80W
External speed control signal	0-5V, 0-10V, 4-20mA for option		
		Condition temperature	0-40°C
		Relative humidity	< 80%




Product Composition and Flow Rate Range						
Flow Rates Peristaltic Pump		Pump Head & Flow Rate (mL/min)				
		New Generation Quick Load Type Pump Head				
Tubing		EasyPumpI/III		EasyPumpII/IV		EasyPumpV/VI(dual channel)
Drive&speed		13#, 14#, 19#, 16#, 25#, 17#, 18#		15#, 24#, 35#, 36#		13#, 14#, 19#, 16#, 25#
LabN1-III	0.1-150 rpm	0.0053~645		0.18~775		0.0053~295
LabN3-III	0.1-350 rpm	0.0053~1505		0.18~1808		0.0053~688
LabN6-III	0.1-600 rpm	0.0053~2580		0.18~3100		0.0053~1180
Tubing		YZ1515x	YZ2515x	AMC1~AMC12(10)	AMC1~AMC12(6)	
Drive&speed		13#, 14#, 19#, 16# 25#, 17#, 18#	15#, 24#	Inner diameter: 0.13-3mm Wall thickness: 0.8-1mm		
LabN1	0.1-150 rpm	0.007~570	0.17~435	0.0002-48(working speed≤150rpm)	0.0002-65(working speed≤150rpm)	
LabN3	0.1-350 rpm	0.007~1330	0.17~1015			
LabN6	0.1-600 rpm	0.007~2280	0.17~1740			

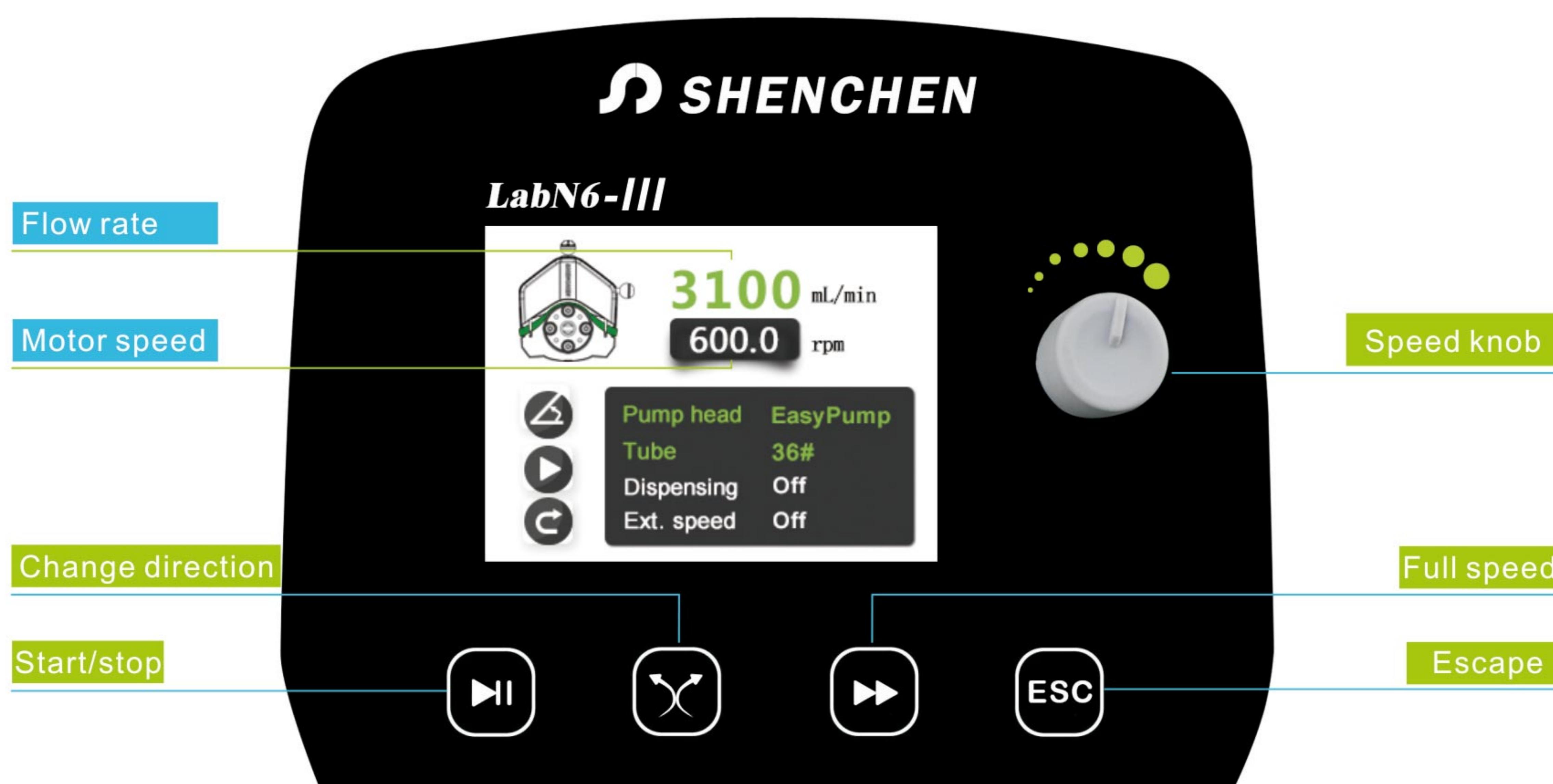

LabN Series Interface and Keypad

Application

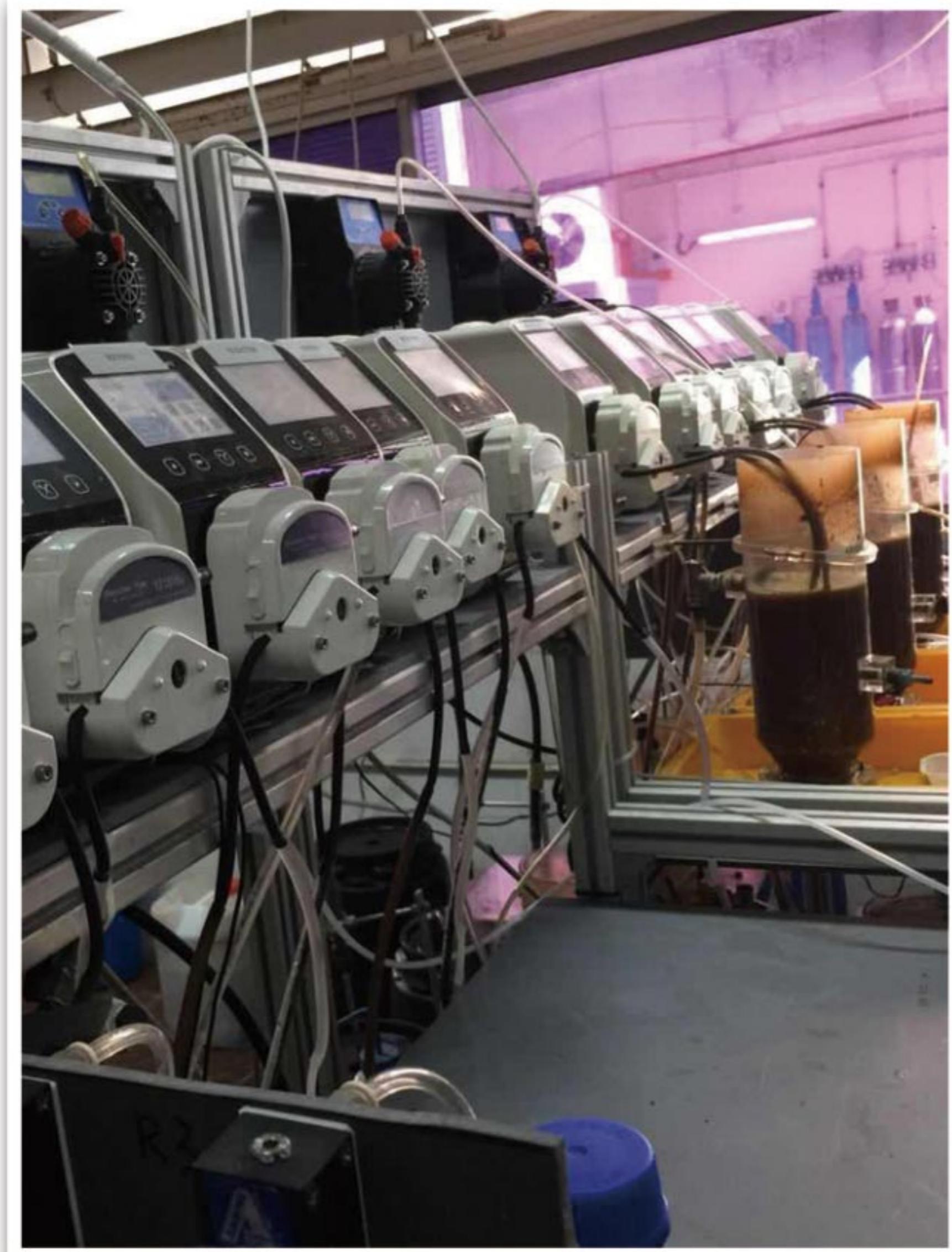

Flow Rates Peristaltic Pump

N6-12L

3 years warranty

upgrade

N6-6L

N6-3L N6-3L/EasyPump



Features

- | 3.2" LCD display.
- | Ultra-quiet drive setting, precise control, small vibration and low noise.
- | Imported button control, menu interface, convenient for users to set various parameters at any time.
- | With timing dispensing function, various external control functions.

Model Number

304
SS Housing

- | N6-3L/EasyPump
- | N6-3L/DZ25-3L
- | N6-6L/DZ25-6L
- | N6-12L/YZ35

Technical Specifications

Flow rate range	N6-3L: 0.0053~3600 mL/min N6-6L: 0.3~6000 mL/min N6-12L: 0.69~12000 mL/min	Output interface	Output motor working status (Open-Collector output)
Speed range	0.1~600 rpm	Power supply	AC 220V±10% 50Hz/60Hz (standard) AC 110V±10% 50Hz/60Hz (optional)
Flow rate accuracy	<±0.5%	Power consumption	N6-3L: <80W ; N6-6L: <180W ; N6-12L: <300W
Speed resolution	0.1rpm	Communication interface	RS232, RS485 support MODBUS protocol (RTU mode)
Fixed time dispensing Function	0.1s~9999h	Motor type	Closed-loop stepper motor
Back suction angle	0~360°	Copy numbers	1~9999 times, setting '0' means unlimited
Control method	Mechanical keypad+Digital knob	Drive dimension (L×W×H)	N6-3L: 223×152×230mm N6-6L: 283×192×264mm N6-12L: 302×222×321mm
Display	3.2" high definition LCD screen	Drive weight	N6-3L: 5.06kg; N6-6L: 7.88kg; N6-12L: 13.01kg
Start/stop, direction signal	Passive switch signal, such as foot pedal switch Active switch signal: 5V, 12V and 24V for option	Condition temperature	0~40°C
External speed control signal	0~5V, 0~10V, 4~20mA for option	Relative humidity	<80%

Product Composition and Flow Rate Range

Drive	Motor Type	Pump Head	Tubing Size	Speed Range(rpm)	Flow Rate(mL/min)
N6-3L	Closed-loop stepper motor	EasyPump	13#, 14#, 19#, 16#, 25#, 17# 18#/ 15#, 24#, 35#, 36#	0.1~600	0.0053~3100
		DZ25-3L	15#, 24#, 35#, 36#		0.211~3600
N6-6L		DZ25-6L	15#, 24#, 35#, 36#		0.3~6000
N6-12L		YZ35	26#, 73#, 82#		0.69~12000



Dispensing Peristaltic Pump

3 years warranty LabF/F Serirs



Features

- | Imported 4.3" industrial grade color LCD screen display, with touch screen control.
- | Can preset dispensing volume, dispensing time, pause time and copy numbers.
- | With intelligent calibration function and online micro adjusting function.
- | The pump can store 60 commonly used filling modes.
- | Back suction angle setting, avoid liquid drop off when the pump stops working.
- | Two working mode: Volume dispensing and speed dispensing (special for viscous liquid)
- | Can communicate with balance, closed-loop control.

Typical Application

- Medicine and chemical dispensing, such as oral liquid, diagnostic reagents.
- Cosmetic dispensing, such as perfume, essential oil.


Suitable Pump Head



EasyPump Series
(Pressure Adjustable)



EasyPump Series
(Fixed Pressure)

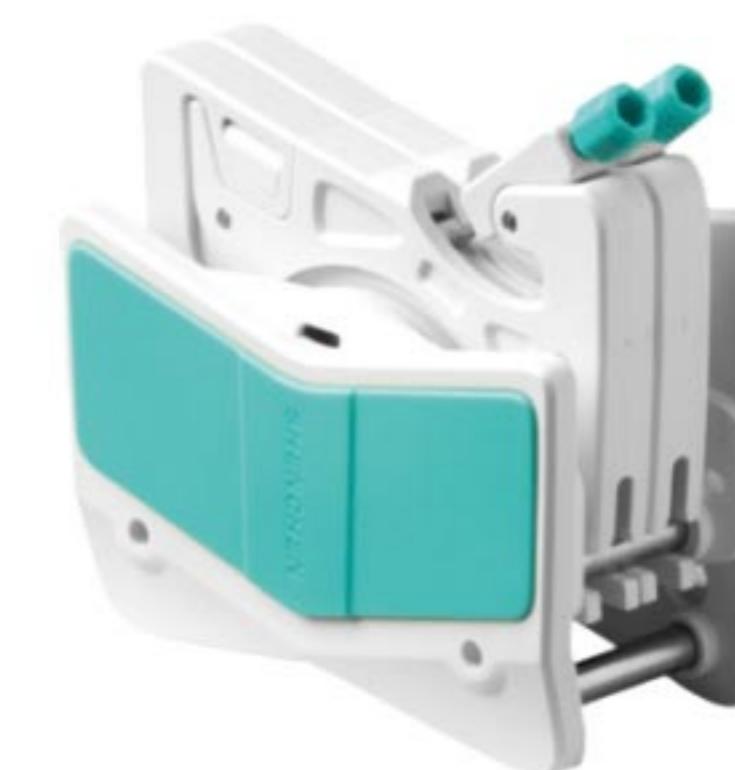

AMC Series



YZ Series



Accessories



① Filling Countersunk

② Filling Nozzle

③ Foot Pedal Switch

Dispensing Peristaltic Pump

F1, F3, F6

3 years warranty



**304
SS Housing**




Technical Specifications

Flow rate range	LabF1-III: 0.0053~775 mL/min	Display	Industrial grade 4.3" color LCD display
	LabF3-III: 0.0053~1808 mL/min		
	LabF6-III: 0.0053~3100 mL/min	Control method	Touch screen and Mechanical keypad
	LabF1/F1: 0.0002~570 mL/min		Keypad lifetime 300,000 times
	LabF3/F3: 0.0002~1330 mL/min	Start/stop, direction signal	Passive switch signal, such as foot pedal switch; Active switch signal: 5-24V universal
	LabF6/F6: 0.0002~2280 mL/min		
Speed range	0.1-600rpm	Output interface	Output motor working status (Open-Collector output)
Speed resolution	0.01rpm		
Dispensing volume range	0.1-9999.99mL	Communication interface	RS232, RS485 Modbus protocol (RTU mode)
Dispensing volume resolution	0.01mL		
Dispensing time	0.1-9999. 99s	Power supply	AC 220V±10% 50Hz/60Hz (Standard)
Pause time	0.1-9999. 99s		AC 110V±10% 50Hz/60Hz (Optional)
Time resolution	0.01s	Drive dimension (L×W×H)	LabF Series: 261.4×157.3×236.9 mm
Copy numbers	1-9999 times, setting'0' means unlimited		F Series: 252×152×243 mm
Back suction angle	0-360°	Drive weight	LabF Series: 4.40 kg
Dispensing accuracy	<±0.5%		F Series: 4.20 kg
Motor type	LabF-III series: Closed-loop stepper motor LabF/F series: Stepper motor	Power consumption	< 80W
		Condition temperature	0-40°C
		Relative humidity	< 80%

Product Composition and Flow Rate Range

Flow Rates Peristaltic Pump		Pump Head & Flow Rate (mL/min)		
		New Generation Quick Load Type Pump Head		
Tubing	EasyPumpI/III	EasyPumpII/IV	EasyPumpV/VI(dual channel)	
Drive&speed	13#, 14#, 19#, 16#, 25#, 17#, 18#	15#, 24#, 35#, 36#	13#, 14#, 19#, 16#, 25#	
LabF1-III	0.1-150 rpm	0.0053~645	0.18~775	0.0053~295
LabF3-III	0.1-350 rpm	0.0053~1505	0.18~1808	0.0053~688
LabF6-III	0.1-600 rpm	0.0053~2580	0.18~3100	0.0053~1180
Tubing	YZ1515x	YZ2515x	AMC1~AMC12(10)	AMC1~AMC12(6)
Drive&speed	13#, 14#, 19#, 16# 25#, 17#, 18#	15#, 24#	Inner diameter: 0.13-3mm Wall thickness: 0.8-1mm	
LabF1/F1	0.1-150 rpm	0.007~570	0.17~435	
LabF3/F3	0.1-350 rpm	0.007~1330	0.17~1015	0.0002-48(working speed≤150rpm)
LabF6/F6	0.1-600 rpm	0.007~2280	0.17~1740	0.0002-65(working speed≤150rpm)



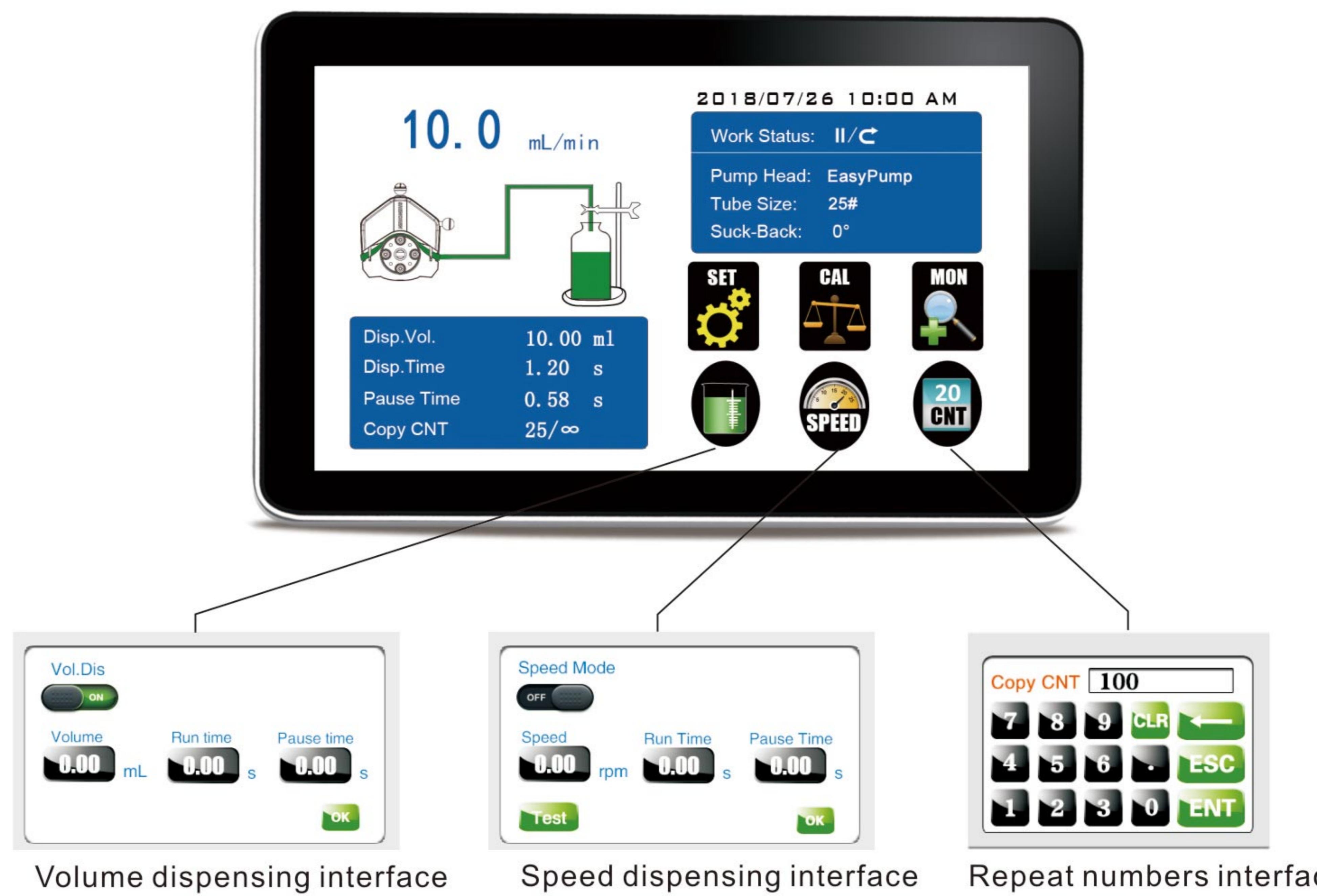

Experimental conditions: standard atmospheric pressure, room temperature at 20°C, the liquid is pure water, no pressure, no suction and lift.

Note: Actually, it is affected by many factors such as transmission medium, inlet and outlet pressure, hose material and error, working environment, etc. This data is for reference only.

LabF6-III/LabF6/F6 Filling Volume Reference Parameter(Media is water)							
Drive	Pump Head	Tubing	Filling Volume (mL)	Filling Time(s)	Filling Accuracy (±%)	Output(pcs/min)	Motor Speed(rpm)
LabF6-III	EasyPump	13#	0.4	1.2	0.8	27	377.36
		13#	1	2.5	0.5	17	452.83
		14#	2	1	0.5	30	444.44
		19#	5	1.2	0.5	27	454.55
		16#	7	1	0.5	30	450.16
		25#	10	0.8	0.8	33	381.29
		25#	15	1	0.5	30	457.55
		25#	20	1.5	0.5	24	406.71
		17#	30	1.2	0.5	27	450.05
		18#	50	1.5	0.5	24	465.12
		15#	15	1.2	0.5	27	416.67
		24#	20	1.2	0.5	27	365.90
		35#	30	1.2	0.5	27	391.34
		36#	50	1.5	0.5	24	387.07
		13#	0.5	0.8	0.8	33	353.77
2*EasyPump	2*EasyPump	13#	1	1.2	0.5	27	471.70
		14#	2	0.6	1	38	370.37
		14#	5	1.5	0.5	24	370.37
		19#	10	1.2	0.5	27	454.55
		16#	15	1.2	0.5	27	401.93
		16#	20	1.5	0.5	24	428.72
		25#	30	1	0.5	30	457.55
		17#	50	1	0.5	30	450.05
		18#	70	1	0.5	30	488.37
		15#	30	1.2	0.5	27	416.67
		24#	50	1.2	0.5	27	457.37
		35#	70	1.2	0.5	27	456.56
		36#	100	1.5	0.5	24	387.07

Dispensing Peristaltic Pump



Features

- Suitable for large flow rate, high efficiency, high precision filling.
- Closed-loop stepper motor drive, accurate control, strong driving force.
- 304 stainless steel drive housing can support the filling line.
- Two working mode: Voleme dispensing and speed dispensing (special for viscous liquid filling)**

Model Number

- F6-3L/EasyPump
- F6-3L/DZ25-3L
- F6-6L/DZ25-6L
- F6-12L/YZ35

Technical Specifications

Flow rate range	F6-3L: 0.0053~3600 mL/min F6-6L: 0.3~6000 mL/min F6-12L: 0.69~12000 mL/min	Keypad lifetime	300,000 times
Speed range	0.1-600 rpm	Start/stop, direction signal	Passive switch signal, such as foot pedal switch
Speed resolution	0.01 rpm	Output interface	Active switch signal: 5-24V universal
Dispensing volume range	0.1-9999.99 mL	Communication interface	Output motor working status (Open-Collector output)
Dispensing volume resolution	0.01 mL	Power supply	Rs232, RS485 support Modbus protocol (RTU mode)
Dispensing time	0.1-9999.99 s		AC 220V±10% 50Hz/60Hz (Standard)
Pause time	0.1-9999.99 s		AC 110V±10% 50Hz/60Hz (Optional)
Time resolution	0.01 s	Drive dimension	F6-3L: 223×152×231mm
Copy numbers	1-9999 times, setting '0' means unlimited		F6-6L: 283×192×274mm
Back suction angle	0-360°		F6-12L: 302×222×331mm
Dispensing accuracy	<±0.5%	Drive weight	F6-3L: 5.02kg; F6-6L: 7.85kg; F6-12L:13.14kg;
Motor type	Closed-loop stepper motor	Power consumption	F6-3L:<80W ; F6-6L:<180W ; F6-12L:<300W
Display	Industrial grade 4.3" color LCD display	Condition temperature	0-40°C
Control method	Touch screen and Mechanical keypad	Relative humidity	<80%

Product Composition and Flow Rate Range

Drive	Motor Type	Pump Head	Tubing Size	Speed Range(rpm)	Flow Rate(mL/min)
F6-3L	Closed-loop stepper motor	EasyPump	13#, 14#, 19#, 16#, 25#, 17# 18#/ 15#, 24#, 35#, 36#	0.1-600	0.0053~3100
		DZ25-3L	15#, 24#, 35#, 36#		0.211~3600
F6-6L		DZ25-6L	15#, 24#, 35#, 36#		0.3~6000
F6-12L		YZ35	26#, 73#, 82#		0.69~12000



Experimental conditions: standard atmospheric pressure, room temperature at 20°C, the liquid is pure water, no pressure, no suction and lift.
 Note: Actually, it is affected by many factors such as transmission medium, inlet and outlet pressure, hose material and error, working environment, etc. This data is for reference only.

F6-3L/6L/12L Filling Volume Reference Parameter(Media is water)							
Drive	Pump Head	Tubing	Filling Volume (mL)	Filling Time(s)	Filling Accuracy (±%)	Output(pcs/min)	Motor Speed(rpm)
F6-3L	EasyPump	13#	0.4	1.2	0.8	27	377.36
		13#	1	2.5	0.5	17	452.83
		14#	2	1	0.5	30	444.44
		19#	5	1.2	0.5	27	454.55
		16#	7	1	0.5	30	450.16
		25#	10	0.8	0.8	33	381.29
		25#	15	1	0.5	30	457.55
		25#	20	1.5	0.5	24	406.71
		17#	30	1.2	0.5	27	450.05
		18#	50	1.5	0.5	24	465.12
		15#	15	1.2	0.5	27	416.67
		24#	20	1.2	0.5	27	365.90
		35#	30	1.2	0.5	27	391.34
		36#	50	1.5	0.5	24	387.07
	2*EasyPump	13#	0.5	0.8	0.8	33	353.77
		13#	1	1.2	0.5	27	471.70
		14#	2	0.6	1	38	370.37
		14#	5	1.5	0.5	24	370.37
		19#	10	1.2	0.5	27	454.55
		16#	15	1.2	0.5	27	401.93
		16#	20	1.5	0.5	24	428.72
		25#	30	1	0.5	30	457.55
		17#	50	1	0.5	30	450.05
		18#	70	1	0.5	30	488.37
		15#	30	1.2	0.5	27	416.67
		24#	50	1.2	0.5	27	457.37
		35#	70	1.2	0.5	27	456.56
		36#	100	1.5	0.5	24	387.07
	DZ25-3L	15#	16	1	0.5	30	443.54
		24#	30	1.2	1	27	454.88
		35#	150	4	0.6	12	447.94
		36#	200	4	0.6	12	481.80
F6-6L	DZ25-6L	15#	30	1.2	0.5	27	500.00
		24#	50	1.2	0.5	27	454.55
		35#	70	1.2	0.5	27	437.50
		36#	100	1.5	0.5	24	400.00
F6-12L	YZ35	26#	50	1	0.5	30	434.78
		26#	100	2	0.5	20	434.78
		73#	100	1.2	0.5	27	406.50
		73#	150	2	0.5	20	365.85
		82#	200	1.5	0.5	24	400.00
	2*YZ35	73#	100	1.2	0.5	27	203.25
		73#	200	1.5	0.5	24	325.20
		82#	500	2	0.5	20	375.00



F6 Series Filling Accuracy Reference Parameter(Media is water)

Drive	Pump Head	Tubing	Filling Volume	Filling Time	Actual Filling Volume(mL)												Filling Accuracy	
F6-3L	EasyPump I/III		13#	0.1mL	0.5s	0.1015	0.1013	0.1003	0.1013	0.1018	0.1004	0.1009	0.1014	0.1012	0.1008	-0.10%	1.80%	
			13#	0.5mL	1.2s	0.5000	0.4996	0.5022	0.5011	0.4977	0.5018	0.5007	0.5005	0.4980	0.5017	-0.46%	0.44%	
			13#	1mL	2.5s	1.0027	0.9996	0.9984	1.0044	1.0001	0.9980	1.0020	1.0025	0.9998	1.0000	-0.20%	0.44%	
			14#	2mL	1s	2.0038	2.0023	1.9993	1.9979	1.9966	1.9960	1.9958	1.9957	1.9952	1.9948	-0.30%	0.19%	
			19#	5mL	1.2s	5.0057	5.0045	5.0027	4.9994	4.9953	4.9896	5.0007	5.0100	5.0114	5.0090	-0.21%	0.27%	
			16#	7mL	1s	7.00	7.02	7.02	7.03	7.00	7.01	6.99	7.00	6.99	6.99	-0.29%	0.43%	
			25#	10mL	1s	10.02	10.03	10.05	10.04	10.02	10.04	10.01	10.01	10.02	10.01	-0.10%	0.50%	
			25#	20mL	1.2s	20.04	20.06	19.98	19.96	19.98	20.07	20.08	20.10	20.02	19.96	-0.30%	0.50%	
			17#	30mL	1.2s	30.15	30.07	29.96	30.19	30.16	30.04	29.89	30.19	30.15	30.00	-0.37%	0.70%	
			17#	50mL	2s	50.14	50.04	50.28	49.98	50.22	49.88	50.23	50.14	49.98	50.23	-0.24%	0.56%	
			18#	80mL	2.5s	80.05	80.27	80.32	80.05	80.37	79.99	80.31	80.13	80.42	80.08	-0.37%	0.53%	
			18#	100mL	3s	100.22	100.14	100.05	100.03	99.90	99.85	100.03	100.07	99.99	100.05	-0.28%	0.22%	
			15#	16mL	1s	16.05	15.98	16.03	16.02	15.96	16.08	15.92	16.12	15.86	16.13	-0.69%	0.81%	
			15#	20mL	1.2s	19.91	19.92	20.00	20.02	20.02	20.02	20.00	20.02	19.98	19.97	-0.45%	0.10%	
F6-3L	EasyPump II/IV		24#	30mL	1.2s	30.05	30.03	30.04	30.06	30.06	30.05	29.92	29.97	30.08	30.06	-0.27%	0.27%	
			24#	50mL	2s	50.09	50.04	49.84	50.15	50.08	50.07	49.78	50.18	50.14	49.97	-0.44%	0.36%	
			35#	150mL	4s	149.79	150.15	149.75	150.21	149.77	150.08	149.92	150.06	150.07	149.85	-0.36%	0.14%	
			36#	200mL	4s	200.22	200.07	200.15	200.26	200.20	200.20	199.94	200.13	200.02	200.22	-0.16%	0.13%	
			15#	16mL	1s	15.97	15.97	16.00	16.06	15.97	15.92	15.96	16.02	16.01	15.97	-0.50%	0.50%	
			24#	30mL	1.2s	29.98	30.12	30.20	30.29	30.06	29.94	30.11	30.18	30.02	29.99	-0.20%	1.00%	
F6-3L	DZ25-3L		35#	150mL	4s	149.50	149.50	150.10	150.00	149.80	149.00	149.70	149.50	150.10	149.60	-0.15%	0.60%	
			36#	200mL	4s	199.70	200.30	200.10	200.30	200.00	200.00	200.70	200.80	200.70	200.40	-0.67%	0.07%	
			15#	80mL	4s	80.2	80.0	80.3	80.2	80.1	80.1	80.2	80.1	80.1	80.1	-0.125%	0.375%	
			24#	150mL	4s	80.1	79.9	80.1	79.9	80.0	80.1	19.9	80.1	80.1	80.0	-0.067%	0.200%	
F6-6L	DZ25-6L		35#	200mL	3.5s	149.9	150.2	150.0	150.1	150.0	150.3	150.0	150.1	150.3	150.0	-0.400%	0.100%	
			36#	300mL	3.5s	199.7	199.9	200.2	200.2	200.2	199.6	199.8	199.8	200.1	199.9	-0.167%	0.200%	
			15#	80mL	4s	300.4	300.1	300.3	300.3	300.4	300.1	300.3	300.4	300.6	300.2	-0.47%	0.33%	
			24#	150mL	4s	300.2	300.0	300.1	300.2	300.2	299.9	299.5	300.5	300.5	300.6	-0.07%	0.27%	
F6-12L	YZ35		26#	150mL	3s	499.5	498.3	501.5	501.9	500.0	501.8	501.1	500.5	501.6	499.6	-0.44%	0.38%	
			73#	300mL	3s	500.0	497.8	499.7	498.3	498.4	499.5	498.2	499.7	499.2	498.4	-0.21%	0.27%	
			82#	500mL	3s	500.8	500.6	500.5	500.5	500.6	500.7	500.5	500.6	500.6	500.6	-0.21%	0.27%	

Low Pulsation Dispensing Peristaltic Pump




IF3

Model Number | IF3

Dimension Drawing (Unit: mm)

Features

- | Closed-loop stepper motor drive, high precision, low pulsation dispensing peristaltic pump.
- | The low pulsation pump head is special for high precision filling. Through the phase difference between the two sets of rollers, make the fluid peaks and valleys complementary, then reduce the pulsation of the fluid.
- | Adaptive pressure tubing space, extend the tubing lifetime effectively.
- | Achieving high precision dispensing of micro flow rate.
- | New flow rate mode, can be used for continuous transferring.
- | Can communicate with balance, closed-loop control.

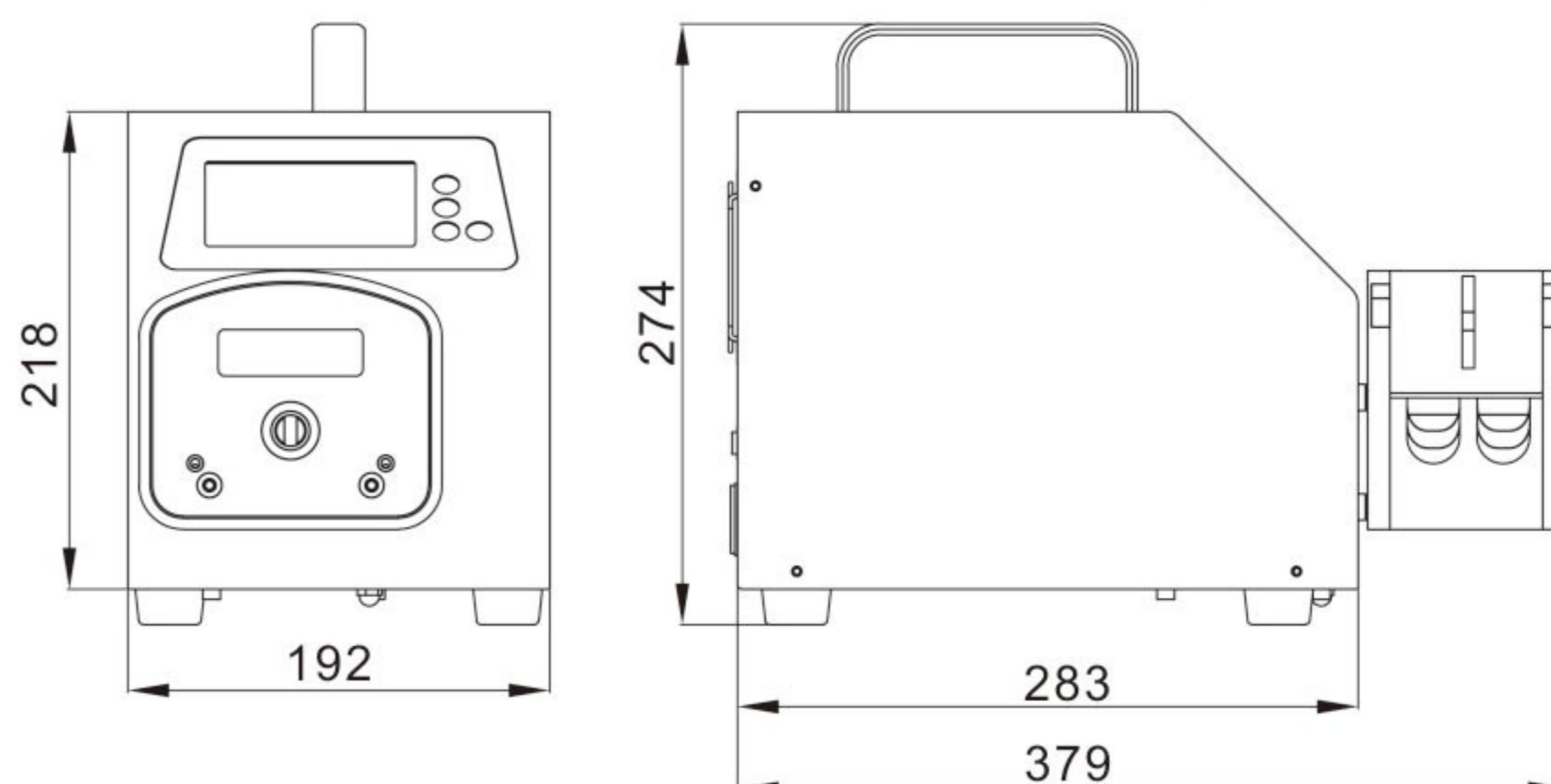

Product Composition and Flow Rate Range

Drive	Motor Type	Pump Head	Tubing	Speed Range(rpm)	Flow Rate(mL/min)
IF3	Closed-loop stepper motor	DY15	13#, 14#, 19#, 16#, 25#, 17#, 18#	0.1-350	0.01~3337
		DY25	15#, 24#, 35#, 36#		0.42~4340

Experimental conditions: standard atmospheric pressure, room temperature at 20°C, the liquid is pure water, no pressure, no suction and lift.

Note: Actually, it is affected by many factors such as transmission medium, inlet and outlet pressure, hose material and error, working environment, etc. This data is for reference only.


Filling Volume Reference Parameter(Media is water)

Tubing	Filling Volume	Filling Time	Actual Filling Volume(mL)												Filling Accuracy	
			0.5003	0.5010	0.4982	0.5015	0.5015	0.5010	0.5014	0.5005	0.4981	0.4998	0.4977	0.5008		
14#	0.5mL	0.5s	0.5025	0.4983	0.5024	0.5021	0.5020	0.5004	0.5005	0.4997	0.4977	0.5008	-0.37%	0.32%		
			1.0037	1.0000	1.0026	1.0003	1.0031	1.0015	1.0000	1.0007	1.0023	1.0035	-0.44%	0.42%		
14#	1mL	1s	1.0002	1.0020	0.9968	1.0010	0.9986	1.0026	1.0030	1.0029	1.0023	1.0017	-0.58%	0.60%		
			2.0020	2.0060	2.0039	2.0037	2.0047	1.9944	2.0062	2.0040	1.9951	2.0060	-0.33%	0.31%		
14#	2mL	1s	2.0036	1.9940	2.0024	2.0086	1.9946	2.0005	1.9951	2.0068	1.9971	1.9965	-0.47%	0.13%		
			3.0058	2.9935	3.0092	3.0008	2.9927	3.0072	3.0046	2.9914	3.0040	3.0000	-0.57%	0.20%		
19#	3mL	1s	2.9902	3.0048	3.0008	2.9919	2.9984	3.0012	3.0084	2.9915	2.9990	3.0051	-0.39%	0.50%		
			5.0215	4.9982	5.0145	5.0038	4.9864	5.0183	4.9962	5.0227	5.0158	4.9806	-0.20%	0.40%		
16#	5mL	1s	5.0248	5.0070	4.9861	5.0108	4.9995	5.0080	5.0044	4.9868	5.0231	4.9977	-0.47%	0.13%		
			10.02	10.04	10.00	9.98	10.02	10.03	10.04	10.04	10.01	9.99	-0.57%	0.20%		
25#	10mL	1s	10.00	10.02	10.02	10.04	9.99	9.98	10.00	10.04	10.02	10.03	-0.60%	0.10%		
			15.02	14.97	14.97	14.94	15.00	15.00	14.94	14.95	14.96	14.99	-0.67%	0.33%		
25#	15mL	1s	14.96	14.95	14.96	14.99	14.98	14.94	14.93	14.93	14.98	14.96	-0.43%	0.43%		
			29.92	30.01	29.99	29.91	29.83	29.86	29.91	29.91	29.91	29.89	-0.20%	0.30%		
17#	30mL	1s	29.88	29.96	30.03	30.06	30.02	30.06	29.96	29.96	29.96	29.83	-0.20%	0.30%		
			99.80	99.90	99.80	99.90	100.00	100.10	100.10	100.00	100.10	100.10	-0.20%	0.10%		
18#	100mL	2.5s	99.80	100.00	100.10	100.20	99.90	100.10	100.00	100.30	100.10	100.10	-0.47%	0.13%		
			9.97	10.01	9.96	9.98	10.01	9.96	9.97	9.99	9.96	9.96	-0.43%	0.43%		
15#	10mL	1s	10.00	9.94	9.97	10.00	9.96	9.97	10.00	10.00	10.01	9.97	-0.60%	0.10%		
			29.90	29.80	30.10	29.90	29.90	30.00	29.90	29.90	30.00	29.90	-0.67%	0.33%		
24#	30mL	1s	69.80	69.90	69.90	69.90	69.70	70.10	70.00	69.70	69.90	69.90	-0.43%	0.43%		
			69.90	70.10	70.00	69.90	70.30	69.70	70.10	69.70	69.90	69.90	-0.20%	0.10%		
36#	100mL	2s	99.80	99.90	99.80	99.80	100.00	99.90	99.90	99.90	99.90	99.90	-0.20%	0.10%		
			99.80	99.90	99.90	99.90	100.00	99.80	100.10	99.90	99.90	99.90	-0.20%	0.10%		



Desktop Filling System

3 years warranty

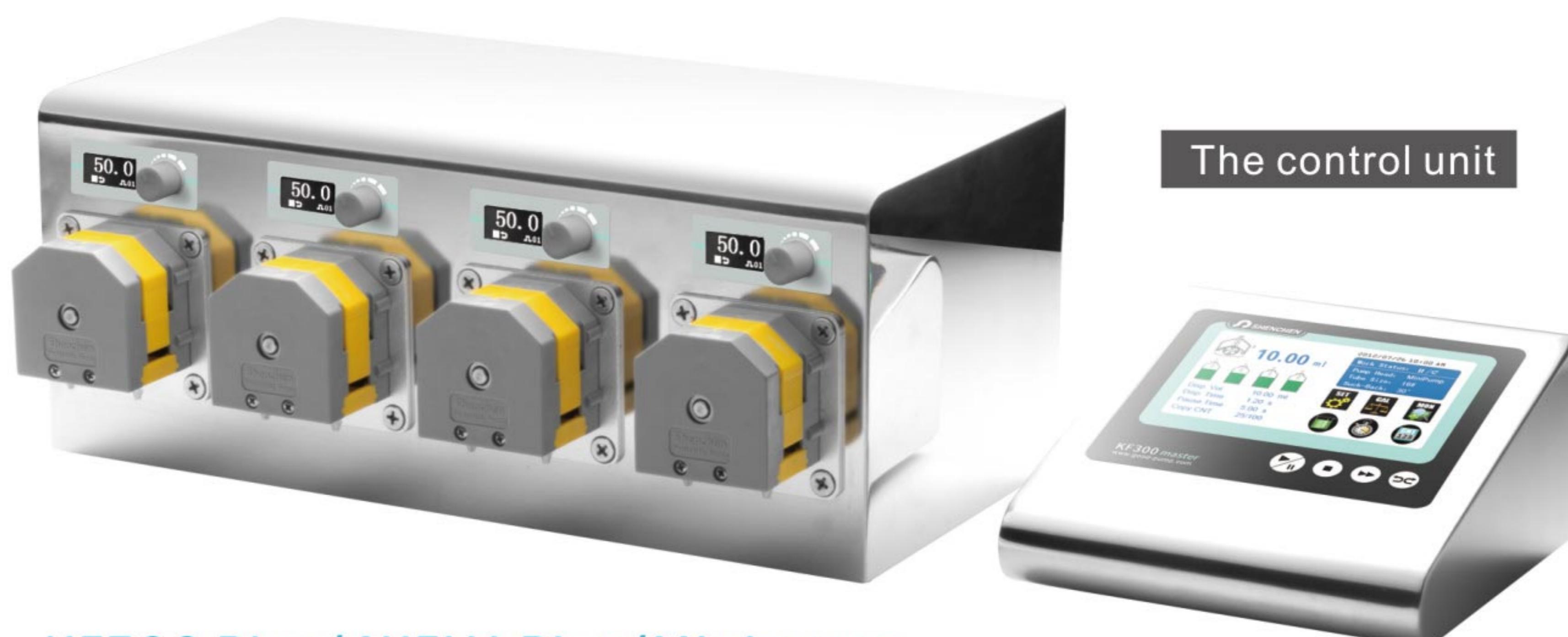

KF300 Plus/4KFU1 Plus/Minipump

The control unit

3 years warranty



KF300 Plus/4KFU3 Plus/Handypump

The control unit

3 years warranty



KF300 Plus/4KFU2 Plus/KT15

The control unit

Features


- Integrated design, a controller can control 1-32 units.
- Compact structure, compact size, beautiful appearance and saving space
- Suitable for high precision micro liquids filling in laboratory.

Model Number | KF300 Plus Series

Typical Application

- Laboratory liquid distribution
- Diagnostic reagent components
- Medium dispensing

External Control Schematic Diagram

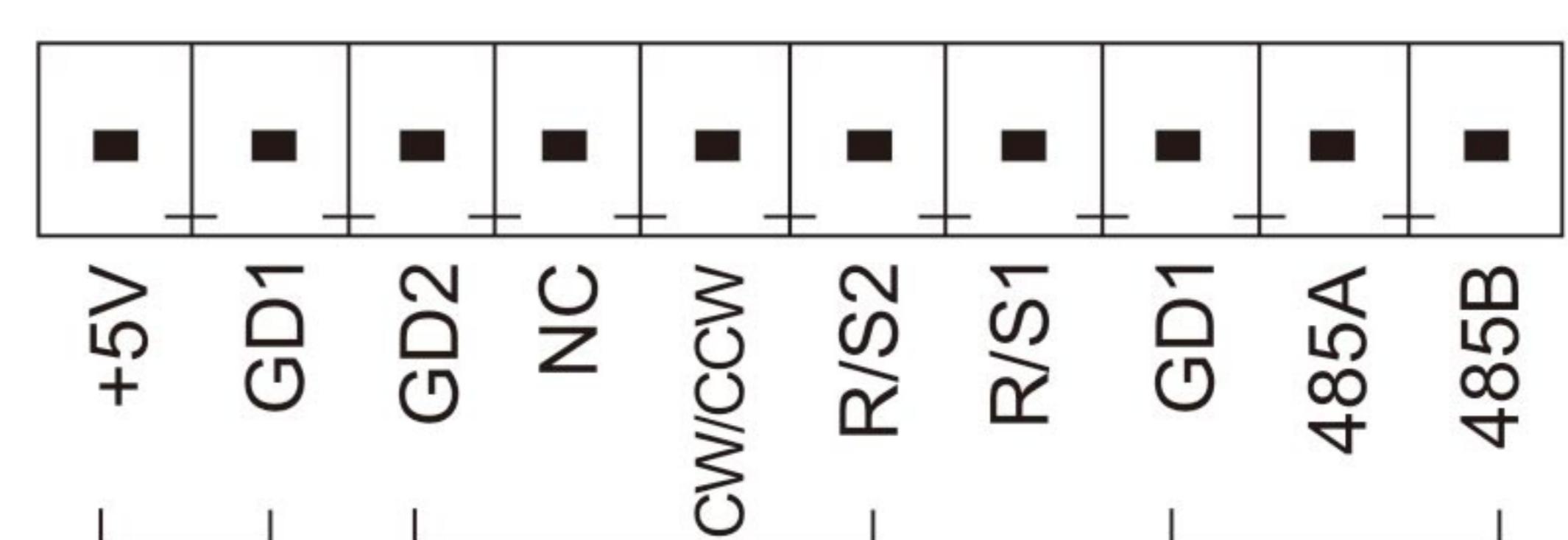



Technical Specifications

	KF300 Plus/4KFU1 Plus/Minipump:0.1-300 rpm	Keypad lifetime	300,000 times
Speed range	KF300 Plus/4KFU3 Plus/Handypump:0.1-300 rpm	Control mode	Touch screen and Mechanical keypad
	KF300 Plus/4KFU2 Plus/KT15:0.1-350 rpm		Passive switch signal, such as :foot
Filling volume	0.01-9999.99 ml	Start/stop, direction signal	pedal switch; Active switch signal, 5-24V universal
Dispensing time	0.1-9999.99 s		
Interval time	0.1-9999.99 s	Output interface	Output motor working status (Open-collector output)
Volume resolution	0.01 ml	Communication interface	RS485 support Modbus protocol (RTU mode)
Time resolution	0.01 s		
Copy numbers	1-9999 times,0 represent unlimited		
Back suction angle	0-3600°	Power supply	AC 100V-240V(50Hz/60Hz)
Filling accuracy	<±0.5%		
Filling units	1-32 units	Condition Temperature	0-40°C
Display	4.3"Industrial grade –color LCD display	Relative humidity	<80%

Experimental conditions:standard atmospheric pressure, room temperature at 20°C,the liquid is pure water, no pressure, no suction and lift.

Note: Actually, it is affected by many factors such as transmission medium, inlet and outlet pressure, hose material and error, working environment, etc. This data is for reference only.

KF300 Filling Volume Reference Parameter(Media is water)

Drive	Pump Head	Tubing	mL/r	Filling Volume (mL)	Filling Time(s)	Filling Accuracy (±%)	Output(pcs/min)	Motor Speed(rpm)
KF300 Plus	MiniPump	13#	0.027	0.06	0.8	3	33	166.67
		13#	0.027	0.1	1	2	30	222.22
		13#	0.027	0.2	2	0.8	20	222.22
		14#	0.114	0.8	1.5	0.8	24	280.70
		14#	0.114	1	2	0.5	20	263.16
		19#	0.248	2	2	0.5	20	241.94
		16#	0.377	3	2	0.5	20	238.73
	KT15	14#	0.156	1	1.5	1	24	256.41
		14#	0.156	2	3	0.5	15	314.69
		19#	0.286	3	2	0.8	20	314.47
		16#	0.477	5	2	0.8	20	225.08
		25#	0.933	7	2	1	20	241.94
	HandyPump (6 rollers)	13#	0.031	0.05	0.4	3	43	241.94
		13#	0.031	0.1	0.8	2	33	258.06
		13#	0.031	0.2	1.5	1	24	353.77
		13#	0.031	0.5	3	0.5	15	322.58
		14#	0.126	1	1.5	0.8	24	317.46
		19#	0.215	2	2.5	0.5	17	223.26
		16#	0.345	3	2	0.5	20	260.87
	HandyPump (3 rollers)	13#	0.033	0.2	1.5	1	24	242.42
		14#	0.187	1.5	2	0.8	20	240.64
		19#	0.371	3	2	0.5	20	242.59
		16#	0.636	5	2	0.8	20	235.85
		25#	1.219	10	2	0.5	20	246.10



Split Type Filling System



3 years warranty

CF600 PlusII

**304
SS Housing**

Features

- | Split design, adding or deleting filling units freely.
- | One controller can control 1-32 filling units.
- | Each filling unit can receive stop filling signal when bottle absent.
- | Controller screen display the working status of the filling units.


Model Number

| CF600 Plus, CF600 PlusII




Typical Application

- | Laboratory dispensing
- | Industry filling

External Control Schematic Diagram

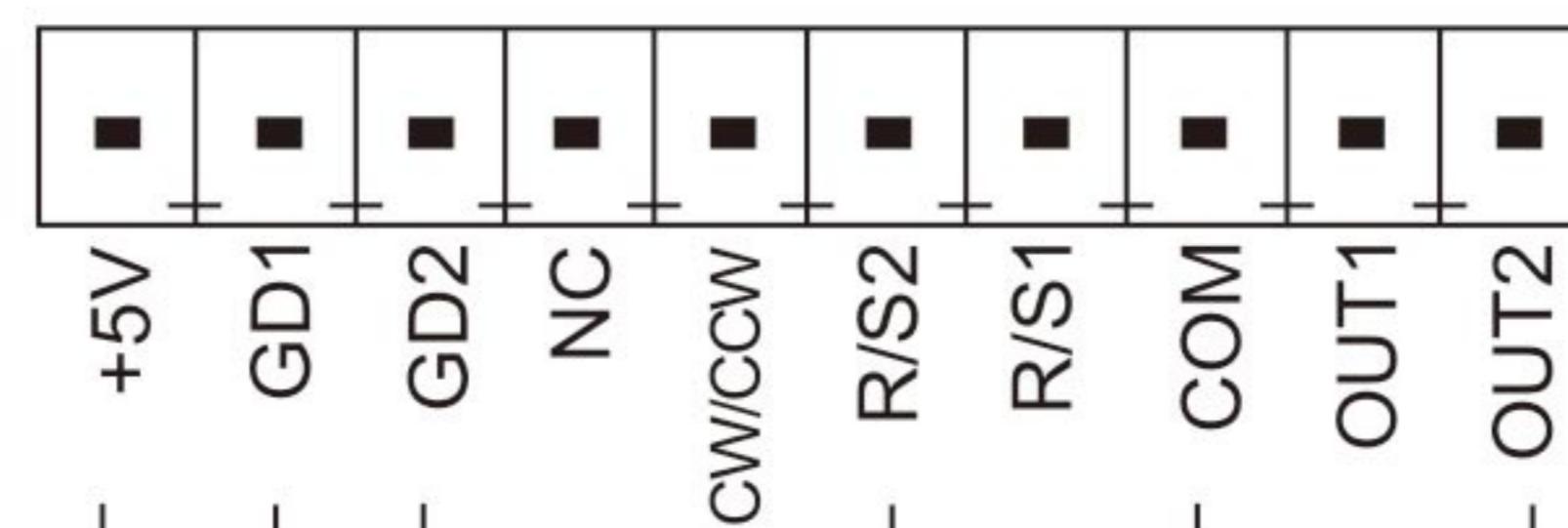

Dimension Drawing(Unit: mm)

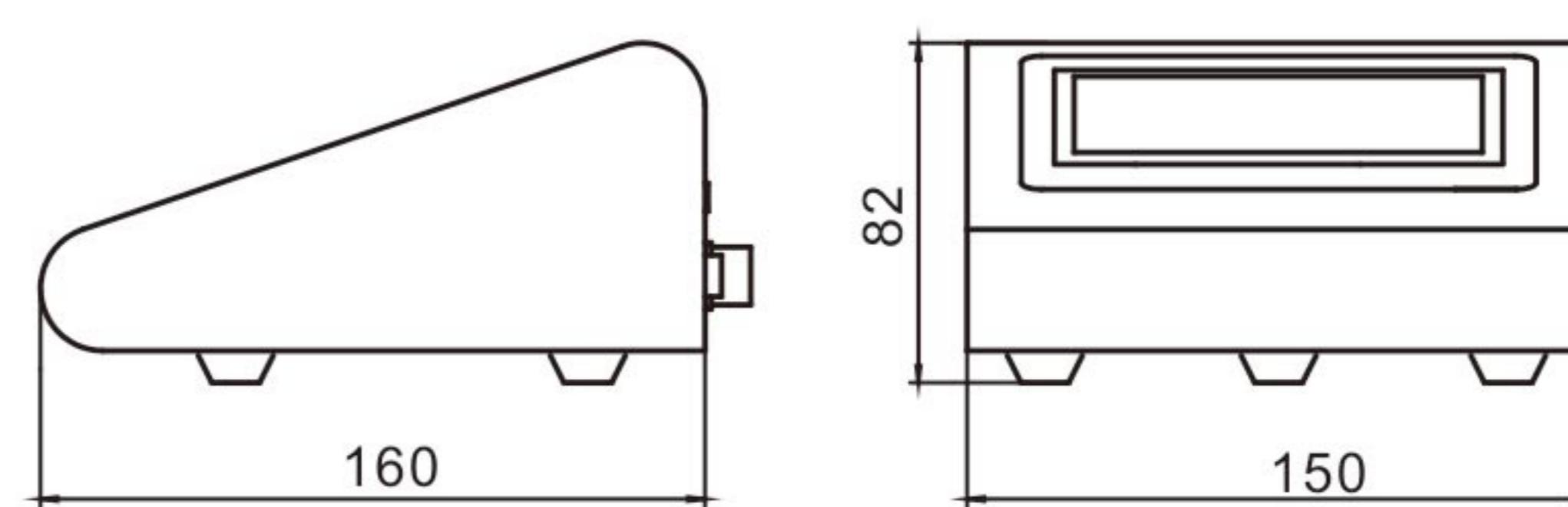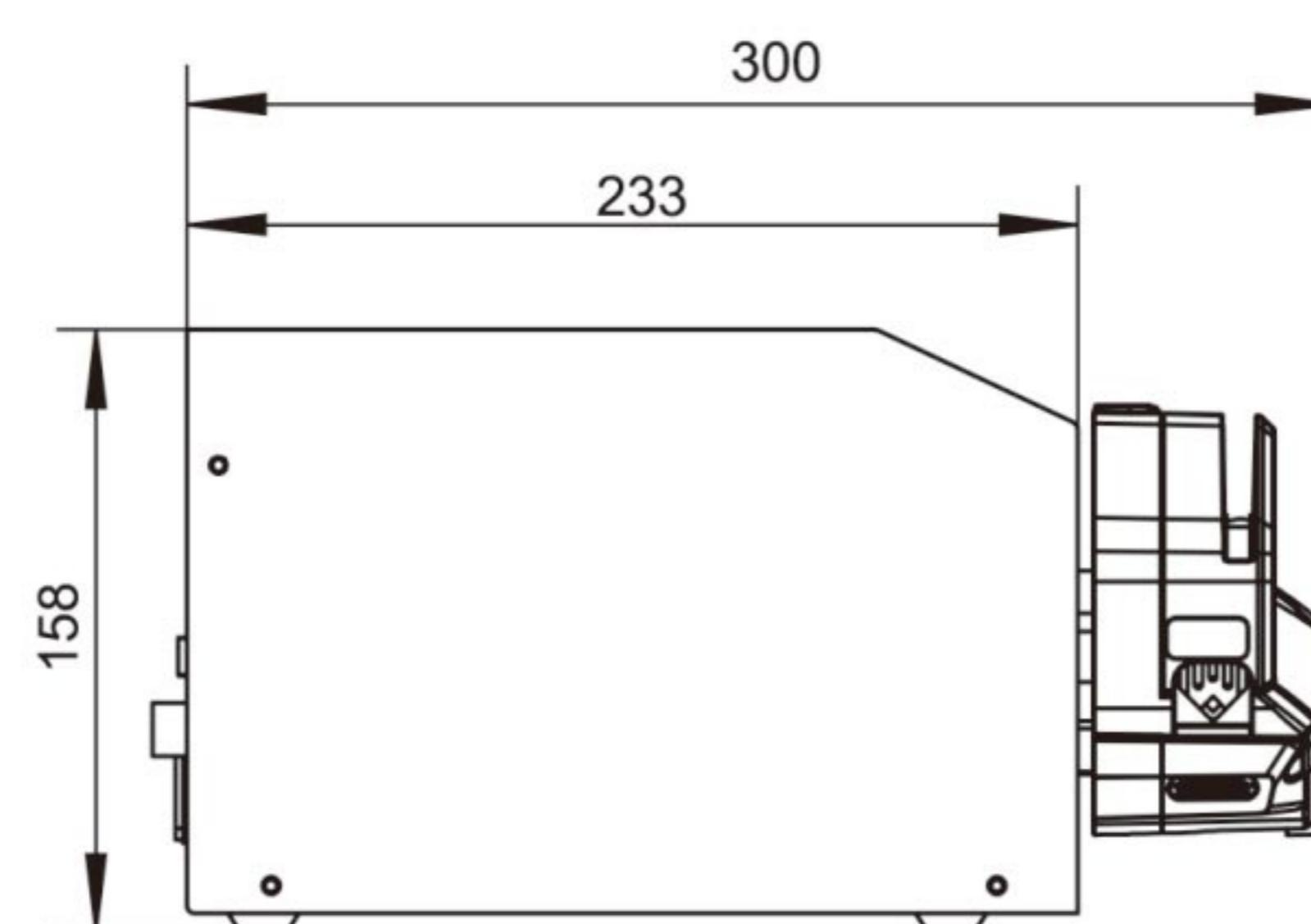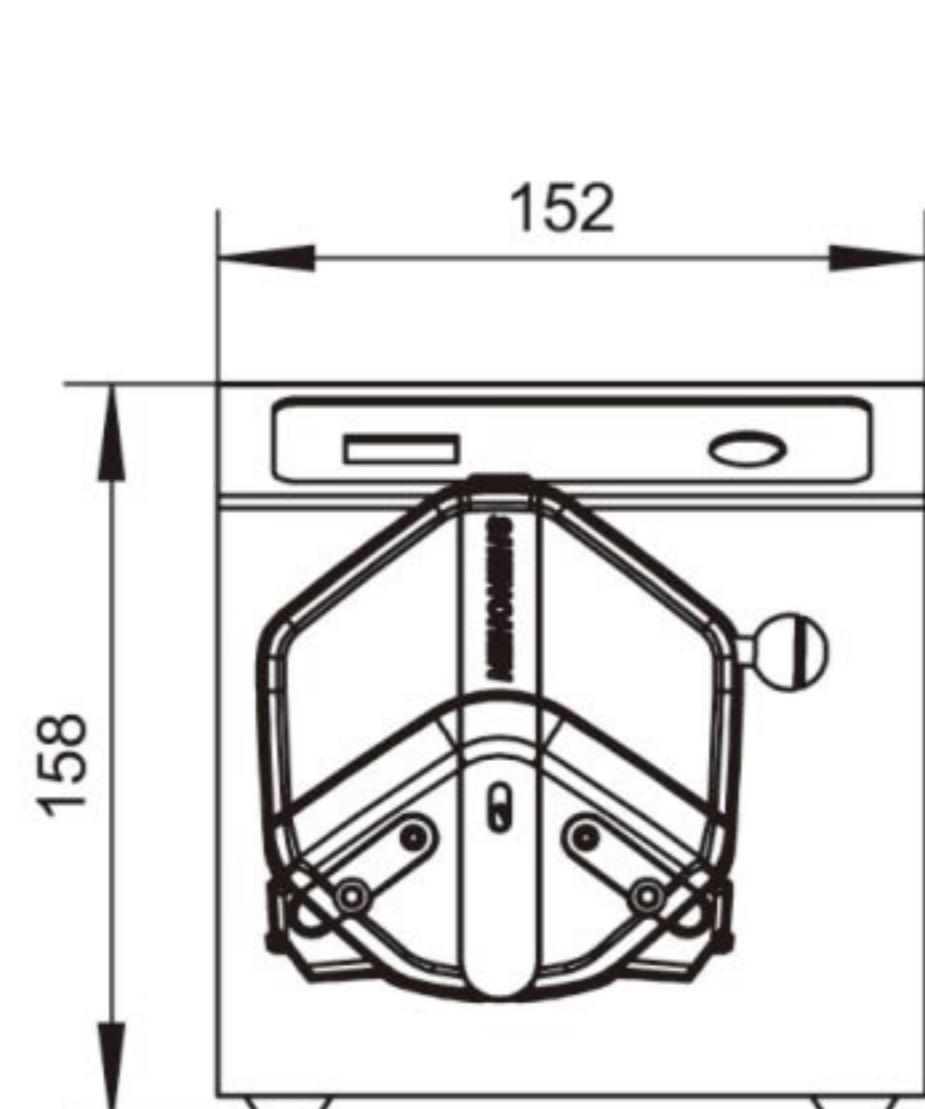

Drive Dimension Drawing

Controller Dimension Drawing



Filling Unit ①

Filling Unit ②

• • •

Filling Unit ⑯

Split Type Filling System

3 years warranty

CF350 Plus

304
SS Housing



Features


- | Closed-loop stepper motor drive, precision control.
- | Split type design, space-saving, easy installation.
- | Use low pulsation pump head to reduce the pulsation of the fluid effectively.
- | Suitable for micro volume high precision filling.

Model Number



- | CF350 Plus
- | CF600 PlusIII, CF600 PlusIV




Dimension Drawing (Unit: mm)

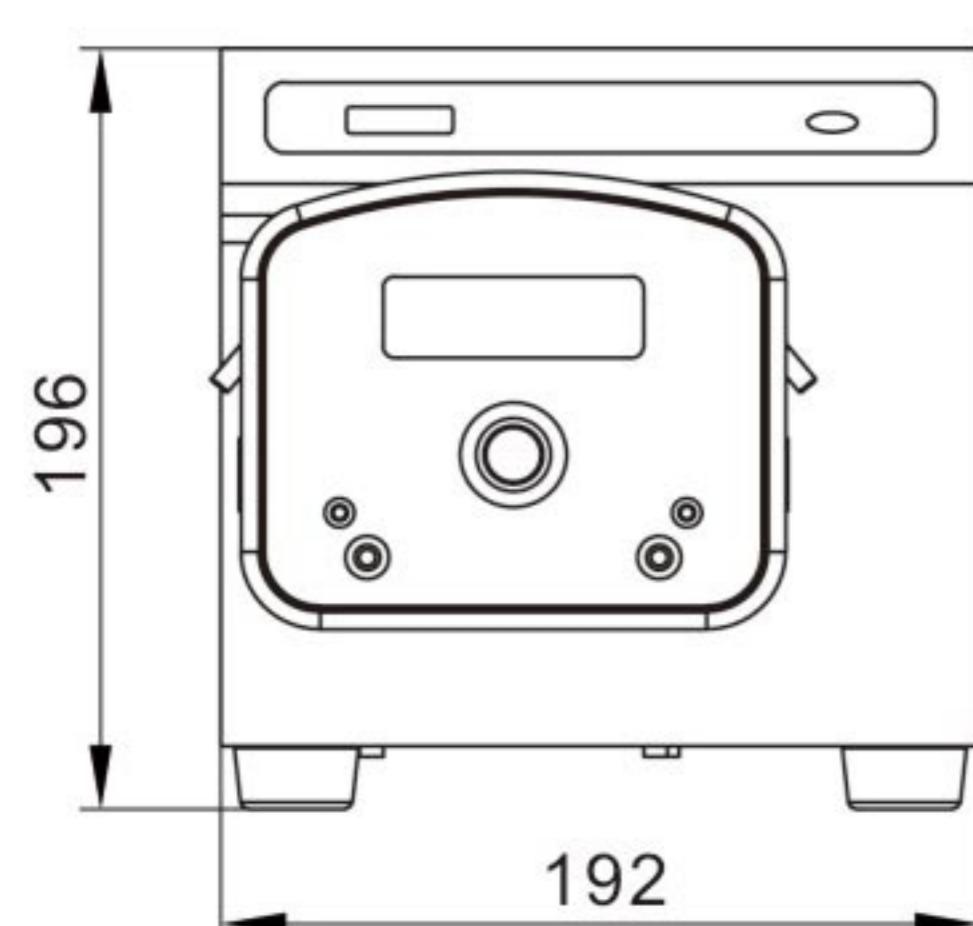

Drive Dimension Drawing

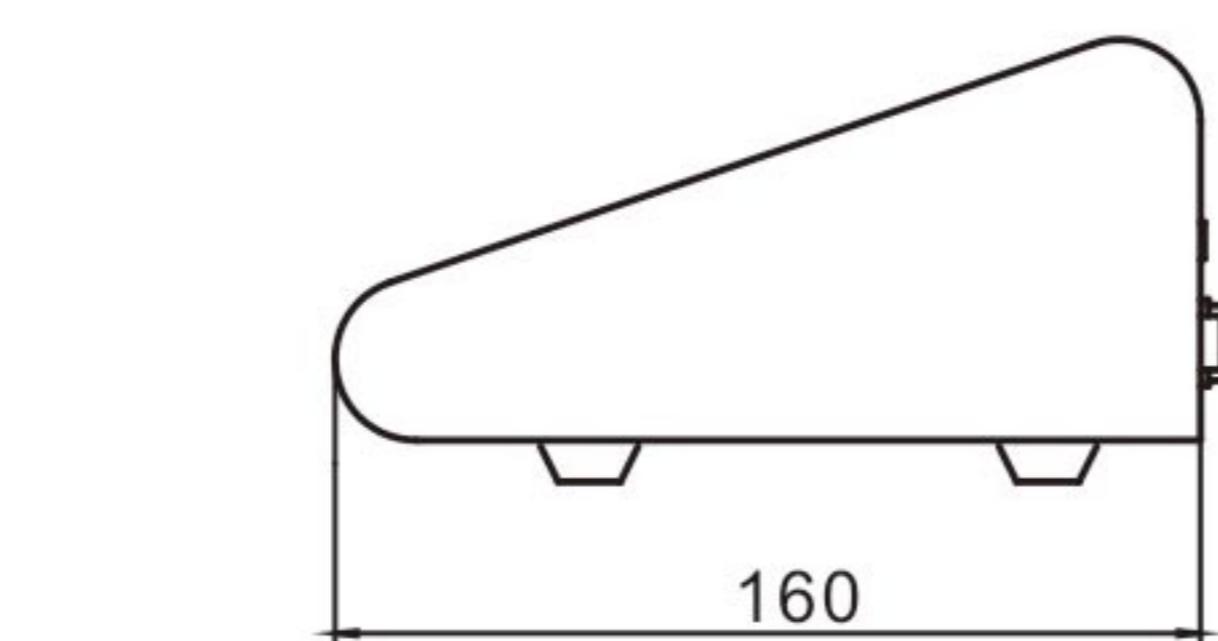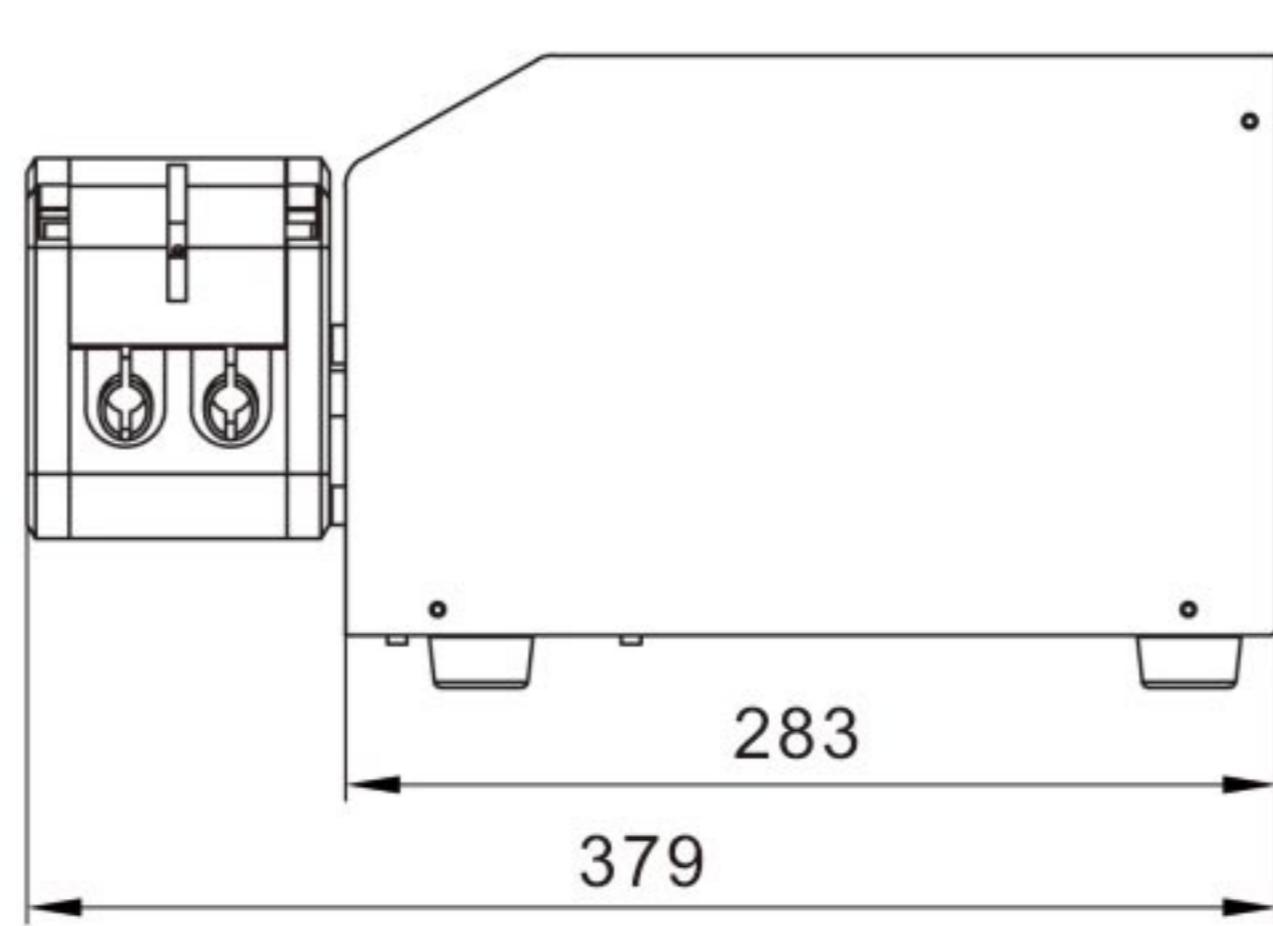

Controller Dimension Drawing

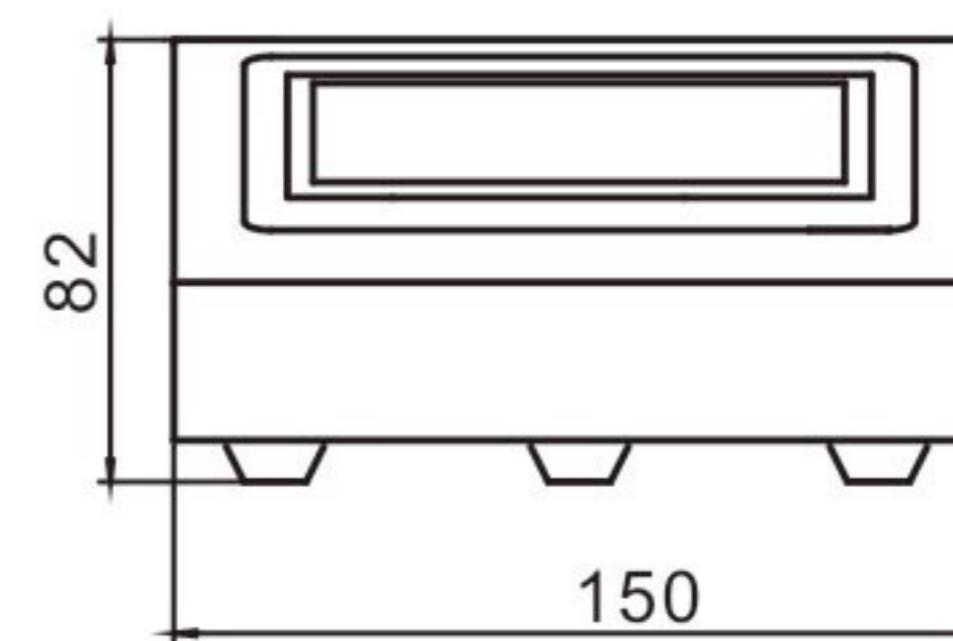

Product Composition and Flow Rate Range						
Drive	Pump Head	Flow Rate (mL/min)	Motor Type	Drive Dimension (L×W×H)	Power Consumption	Drive Weight
CF600 Plus	YZ1515x	0.007-2280	Stepper motor	237.4×152×158mm	Each unit<50W	4.2kg
	YZ2515x	0.17-1740				
CF600 PlusII	EasyPump	0.0053-3100		237.4×152×158mm	Each unit<80W	4.95kg
	DZ25-3L	0.211-3600				
CF600 PlusIII	DZ25-6L	0.3-6000	Closed-loop stepper motor	283×192×196mm	Each unit<180W	7.8kg
CF350 Plus	DY15	0.01-3337		283×192×196mm	Each unit<180W	7.8kg
	DY25	0.42-4340		310×228×248mm		
CF600 PlusIV	YZ35	0.69-12000		310×228×248mm	Each unit<300W	11.9kg



Integrated Filling System

DF600 PlusII

3 years warranty

304
SS Housing



Features


- Integrated design, one controller can control 8 groups (total 32) filling units.
- It can support the production filling line, store 30 commonly used filling modes.
- Each filling unit can independently receive stop filling signal when bottle absent.

Model Number

- DF600 Plus, DF600 PlusII

Typical Application

- Pharmaceutical, health product filling, diagnostic reagent dispensing.
- Food, beverage filling.
- Cosmetic filling.

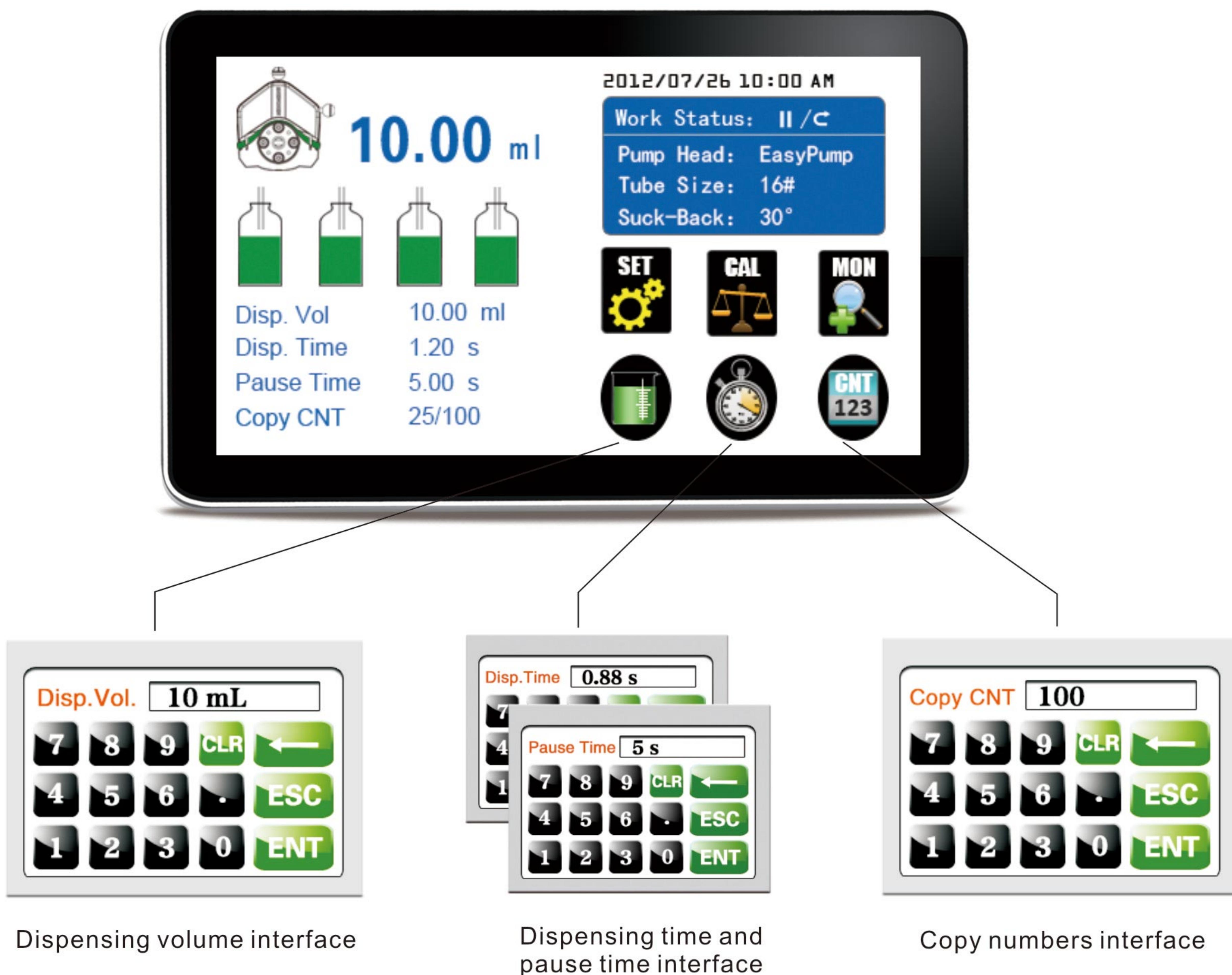

Dispensing volume interface

Dispensing time and pause time interface

Copy numbers interface

Integrated Filling System

DF600 PlusIV

3 years warranty



DF350 Plus

3 years warranty



Product Introduction

DF600 Plus series closed-loop stepper motor control series is integrated type intelligent filling system with high precision and low pulsation driven by closed-loop stepper motor. It is composed by control unit and integrated type filling units, each group have 8 filling units, total 32 filling units. This system use imported 4.3" industrial grade true color display with touch screen control, can preset filling volume, filling time, pause time, copy numbers and back-suction angle. Dynamics display working status, filling data, setting parameter, system configuration display at the same screen; with intelligent calibration and online micro adjust function. Can connect foot pedal and receive switch signal, realize long-distance control. With motor working status output signal, can let other equipment realtime monitor the current filling status of peristaltic pump. This filling system use closed-loop stepper motor, and high precision, low pulsation pump head to make the filling precision higher, up to 0.1%-0.6%, make the micro volume and big volume high precision filling come true.

Model Number

- | DF350 Plus
- | DF600 PlusIII
- | DF600 PlusIV

Typical Application

- | Pharmaceutical, health product filling, diagnostic reagent dispensing.
- | Food, beverage filling.
- | Cosmetic filling.





Integrated Filling System

Product Composition and Flow Rate Range						
Drive	Pump Head	Flow Rate (mL/min)	Motor Type	Drive Dimension (L×W×H)	Power Consumption	Drive Weight
DF600 Plus	YZ1515x	0.007~2280	Stepper motor	706×250×186mm	Each Unit<50W	21.8kg
	YZ2515x	0.17~1740		706×250×186mm	Each Unit<80W	21.8kg
DF600 PlusII	EasyPump	0.0053-3100	Closed-loop stepper motor	910×290×212mm	Each Unit<180W	41.45kg
	DZ25-3L	0.211-3600		910×290×212mm	Each Unit<180W	41.45kg
DF600 PlusIII	DZ25-6L	0.3-6000	Closed-loop stepper motor	1056×310×248m	Each Unit<300W	58.6kg
DF350 Plus	DY15	0.01-3337				
	DY25	0.42-4340				
DF600 PlusIV	YZ35	0.69-12000				

Technical Specifications						
Speed range	CF350 Plus/DF350 Plus	0.1-350 rpm	Control method	Touch screen and mechanical keypad		
	CF600 Plus/DF600 Plus/CF600 PlusII/			Keypad lifetime		
	DF600 PlusII/CF600 PlusIII/DF600 PlusIII/			300,000 times		
Filling volume range	CF PlusIV/DF600 PlusIV	0.1-600 rpm	Start/stop, direction signal	Passive switch signal, such as foot pedal switch; Active switch signal: 5-24V universal		
	0.1-9999.99 mL					
Filling time range	DF600 Plus	0.1-9999.99 s	Communication interface	RS485 communication		
Pause time range	DF600 Plus	0.5-9999.99 s		Modbus protocol(RTU mode)		
Filling volume resolution	0.01 mL		Output interface	Output motor working status (Open-Collector output)		
Time resolution	0.01 s					
Copy numbers	1-9999, '0' means unlimited		Power supply	AC 220V±10% 50Hz/60Hz (standard)		
Back suction angle	0-3600°			AC 110V±10% 50Hz/60Hz (optional)		
Filling accuracy	<±0.5%		Condition temperature	0-40°C		
Display	4.3" industrial grade true color LCD screen					
			Relative humidity	< 80%		

Dimension Drawing(Unit: mm)

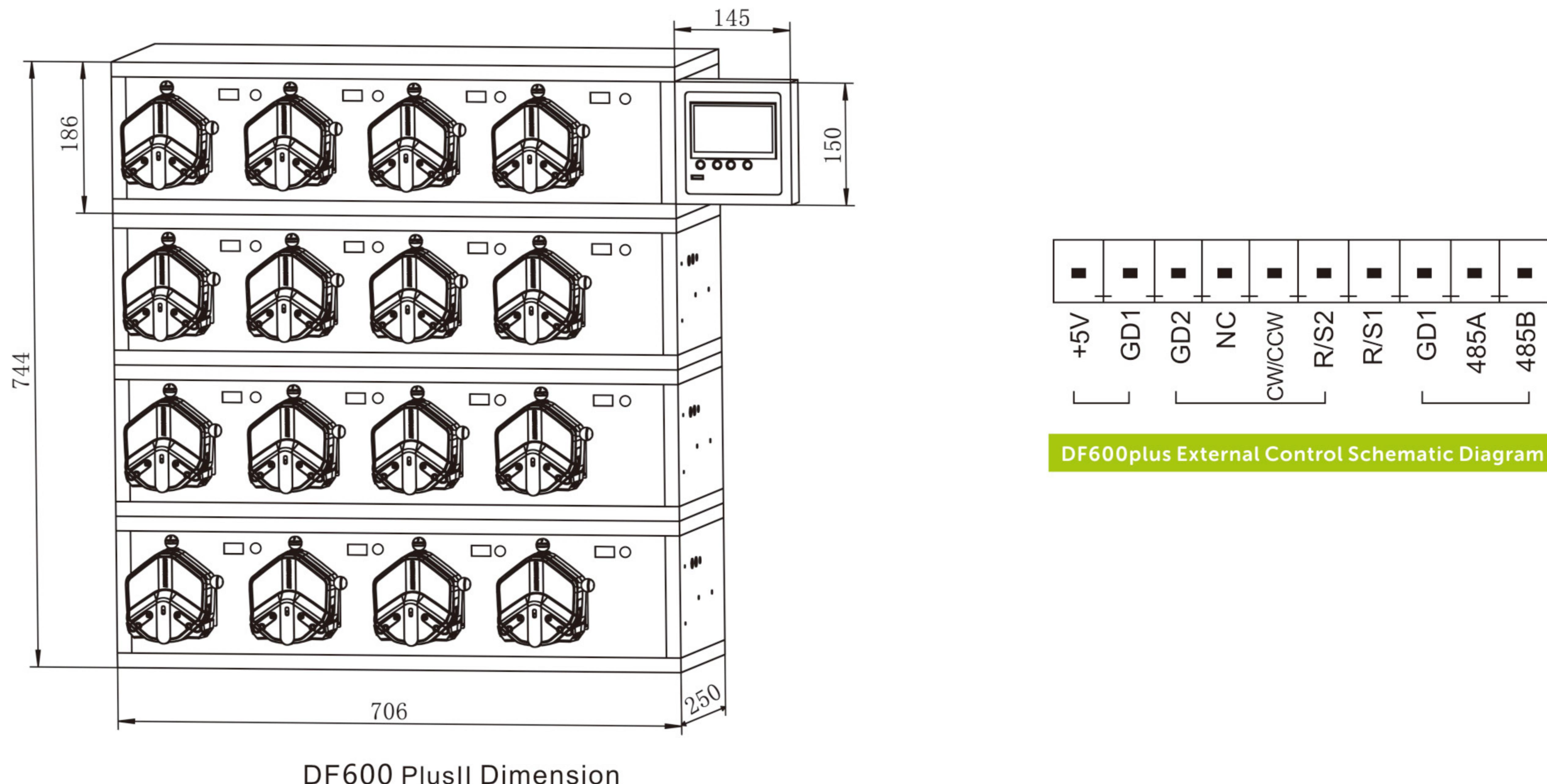

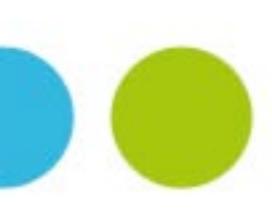

Product Composition and Flow Rate Range						
Drive	Pump Head	Pump Head Material	Tube Clamp Material	Tubing		(Motor Speed 0.1~600rpm) Flow Rate mL/min
				Tubing Size	ID×Wall Thickness	
DF600 Plus CF600 Plus	YZ1515x	PSF PPS		13#	0.8×1.6 (mm)	0.07~42
				14#	1.6×1.6 (mm)	0.27~162
				19#	2.4×1.6 (mm)	0.55~330
				16#	3.1×1.6 (mm)	0.82~492
				25#	4.8×1.6 (mm)	1.7~1020
				17#	6.4×1.6 (mm)	2.9~1740
				18#	7.9×1.6 (mm)	3.8~2280
	YZ2515x	PSF PPS		15#	4.8×2.4 (mm)	1.7~1020
				24#	6.4×2.4 (mm)	2.9~1740
DF600 PlusII CF600 PlusII	EasyPump I/III			13#	0.8×1.6 (mm)	0.0053~32
				14#	1.6×1.6 (mm)	0.027~162
				19#	2.4×1.6 (mm)	0.055~330
				16#	3.1×1.6 (mm)	0.093~560
				25#	4.8×1.6 (mm)	0.197~1180
				17#	6.4×1.6 (mm)	0.333~2000
				18#	7.9×1.6 (mm)	0.430~2580
	EasyPump II/IV			15#	4.8×2.4 (mm)	0.180~1080
				24#	6.4×2.4 (mm)	0.273~1640
				35#	7.9×2.4 (mm)	0.383~2300
				36#	6.4×2.4 (mm)	0.517~3100
				13#	0.8×1.6 (mm)	0.0053~32
				14#	1.6×1.6 (mm)	0.027~162
				19#	2.4×1.6 (mm)	0.055~330
DF600 PlusIII CF600 PlusIII	EasyPump V/VI			16#	3.1×1.6 (mm)	0.093~560
				25#	4.8×1.6 (mm)	0.197~1180
				15#	4.8×2.4 (mm)	2.11~1264
				24#	6.4×2.4 (mm)	3.85~2310
				35#	7.9×2.4 (mm)	5.08~3050
	DZ25-3L	Aluminum alloy/PPS	PP	36#	4.8×2.4 (mm)	6~3600
				15#	4.8×2.4 (mm)	3~1800
				24#	6.4×2.4 (mm)	5.5~3300
				35#	7.9×2.4 (mm)	8~4800
				36#	4.8×2.4 (mm)	10~6000
DF600 PlusIV CF600 PlusIV	YZ35-PPS	Aluminum alloy/PPS	PP	26#	6.4×3.3 (mm)	6.9~4200
				73#	9.6×3.3 (mm)	12.3~7400
				82#	12.7×3.3 (mm)	20~12000
Drive	Pump Head	Pump Head Material	Tube Clamp Material		Tubing	(Motor Speed 0.1~350rpm) Flow Rate mL/min
				Tubing Size	ID×Wall Thickness	
DF350 Plus CF350 Plus	Aluminum alloy	DY15	PP	13#	0.8×1.6 (mm)	0.1~48
				14#	1.6×1.6 (mm)	0.6~223
				19#	2.4×1.6 (mm)	1.3~448
				16#	3.1×1.6 (mm)	2~723
				25#	4.8×1.6 (mm)	4.7~1626
	Aluminum alloy	DY25	PP	17#	6.4×1.6 (mm)	6.4~2230
				18#	7.9×1.6 (mm)	9.5~3337
				15#	4.8×2.4 (mm)	4.2~1480
				24#	6.4×2.4 (mm)	7.6~2670
				35#	7.9×2.4 (mm)	10~3600
				36#	9.6×2.4 (mm)	12.4~4340



Filling Volume Reference (Media is water)							
Drive	Pump Head	Tubing	Filling Volume (mL)	Filling Time (s)	Accuracy (\pm %)	Output (pcs/min)	Motor Speed (rpm)
DF600 PlusII CF600 PlusII	EasyPump	13#	0.4	1.2	0.8	27	377.36
		13#	1	2.5	0.5	17	452.83
		14#	2	1	0.5	30	444.44
		19#	5	1.2	0.5	27	454.55
		16#	7	1	0.5	30	450.16
		25#	10	0.8	0.8	33	381.29
		25#	15	1	0.5	30	457.55
		25#	20	1.5	0.5	24	406.71
		17#	30	1.2	0.5	27	450.05
		18#	50	1.5	0.5	24	465.12
		15#	15	1.2	0.5	27	416.67
		24#	20	1.2	0.5	27	365.90
		35#	30	1.2	0.5	27	391.34
		36#	50	1.5	0.5	24	387.07
2*EasyPump	2*EasyPump	13#	0.5	0.8	0.8	33	353.77
		13#	1	1.2	0.5	27	471.70
		14#	2	0.6	1	38	370.37
		14#	5	1.5	0.5	24	370.37
		19#	10	1.2	0.5	27	454.55
		16#	15	1.2	0.5	27	401.93
		16#	20	1.5	0.5	24	428.72
		25#	30	1	0.5	30	457.55
		17#	50	1	0.5	30	450.05
		18#	70	1	0.5	30	488.37
		15#	30	1.2	0.5	27	416.67
		24#	50	1.2	0.5	27	457.37
		35#	70	1.2	0.5	27	456.56
		36#	100	1.5	0.5	24	387.07
DF600 PlusIII CF600 PlusIII	DZ25-3L	15#	16	1	0.5	30	454.98
		24#	30	1.2	0.5	27	389.61
		35#	150	4	0.5	12	442.91
		36#	200	4	0.5	12	500.00
DF600 PlusIII CF600 PlusIII	DZ25-6L	15#	30	1.2	0.5	27	500.00
		24#	50	1.2	0.5	27	454.55
		35#	70	1.2	0.5	27	437.50
		36#	100	1.5	0.5	24	400.00
DF350 Plus CF350 Plus	DY15	14#	0.5	0.5	0.8	40	94.19
		14#	1	0.6	0.8	38	156.99
		14#	3	1	0.5	30	282.57
		19#	5	1	0.8	30	234.38
		16#	15	1.5	0.5	24	290.42
		25#	20	1	0.5	30	258.29
		17#	30	1	0.5	30	282.53
		18#	100	2.5	0.5	17	251.73
DF600 PlusIV CF600 PlusIV	DY25	15#	10	0.8	0.8	33	177.35
		24#	30	1	0.8	30	235.94
		35#	70	1.5	0.5	24	272.21
		36#	100	2	0.5	20	241.94
DF600 PlusIV CF600 PlusIV	YZ35	26#	50	1	0.5	30	434.78
		26#	100	2	0.5	20	434.78
		73#	100	1.2	0.5	27	406.50
		73#	150	2	0.5	20	365.85
		82#	200	1.5	0.5	24	400.00
2*YZ35	2*YZ35	73#	100	1.2	0.5	27	203.25
		73#	200	1.5	0.5	24	325.20
		82#	500	2	0.5	20	375.00



Basic Peristaltic Pump


LabM Series



3 years warranty



Suitable Pump Head

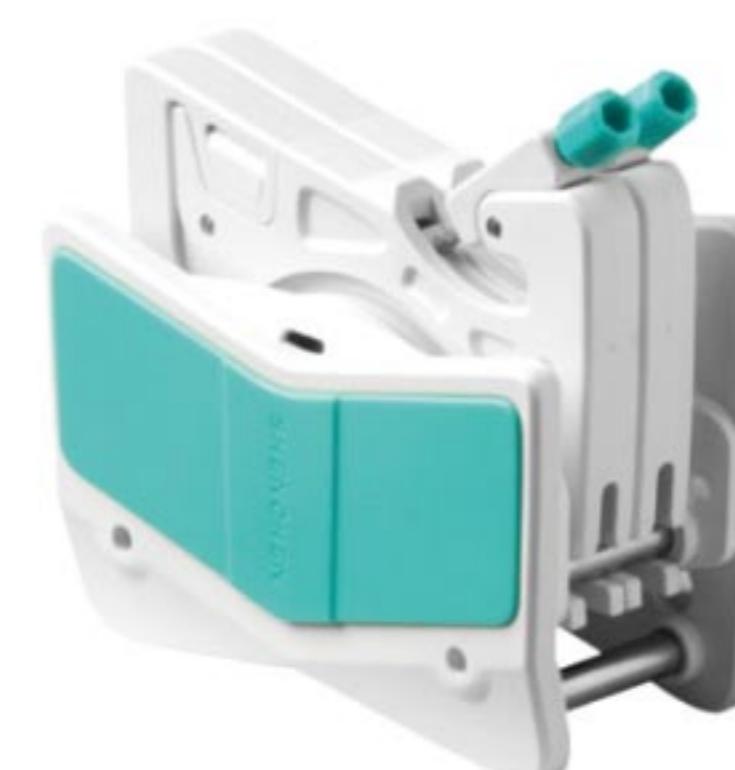EasyPump Series
(Pressure Adjustable)EasyPump Series
(Fixed Pressure)

AMC Series



YZ Series

Model Number | LabM1-III, LabM3-III, LabM6-III**Typical Application**

University laboratory.

Supporting ion chromatography and titrator.

Features

- 3 digital LED display motor speed, mechanical keypad control.
- Timing function, the time range of 0.5 seconds -999 seconds can be used as a simple dispensing.
- Support RS232 and RS485 communication, Modbus protocol.

Technical Specifications

Flow rate range	LabM1-III: 0.0053~775 mL/min	External speed control signal	0-5V, 4-20mA, 0-10V
	LabM3-III: 0.0053~1808 mL/min	Start/stop, direction signal	Passive switch signal, such as foot pedal switch; Active switch signal: 5-24V universal
	LabM6-III: 0.0053~3100 mL/min	Communication interface	RS232, RS458 communication
	LabM1: 0.0002~570 mL/min		Modbus protocol(RTU mode)
	LabM3: 0.0002~1330 mL/min	Power supply	AC 220V±10%, 50Hz/60Hz (standard)
Speed resolution	LabM6: 0.0002~2280 mL/min		AC 110V±10%, 50Hz/60Hz (optional)
	0.1rpm when the speed is 0-100rpm, 1rpm when the speed is 100-600rpm.	Drive dimension	261.4×157.3×236.9mm (L×W×H)
Testing time range	0.5s-999s	Drive weight	4.40 kg
Motor type	LabM-III series: Closed-loop stepper motor	Power consumption	<80W
	LabM series: Stepper motor	Condition temperature	0-40°C
Display	3 digital LED	Relative humidity	<80%
Control method	Mechanical keypad		
Keypad lifetime	300,000 times		

Product Composition and Flow Rate Range

Flow Rates Peristaltic Pump		Pump Head & Flow Rate (mL/min)		
		New Generation Quick Load Type Pump Head		
Drive&speed	Tubing	EasyPumpI/III 13#, 14#, 19#, 16#, 25#, 17#, 18# 0.0053~645	EasyPumpII/IV 15#, 24#, 35#, 36# 0.18~775	EasyPumpV/VI(dual channel) 13#, 14#, 19#, 16#, 25# 0.0053~295
	LabM1-III	0.1-150 rpm	0.0053~1505	0.0053~688
	LabM3-III	0.1-350 rpm	0.18~3100	0.0053~1180
Drive&speed	Tubing	YZ1515x 13#, 14#, 19#, 16# 25#, 17#, 18#	YZ2515x 15#, 24#	AMC1~AMC12(10) Inner diameter: 0.13-3mm Wall thickness: 0.8-1mm
	LabM1	0.1-150 rpm	0.17~435	0.0002-48(working speed≤150rpm)
	LabM3	0.1-350 rpm	0.17~1015	0.0002-65(working speed≤150rpm)
	LabM6	0.1-600 rpm	0.17~1740	



Basic Peristaltic Pump



Features

- Closed-loop stepper motor drive, powerful and maintenance-free.
- 3 digital LED display motor speed, rotary encoded switch control.
- Suitable for industrial sites, transfer liquid with large flow and high precision.

Model Number

- M6-3L/EasyPump
- M6-3L/DZ25-3L
- M6-6L/DZ25-6L
- M6-12L/YZ35

Technical Specifications

Flow rate range	M6-3L: 0.0053~3600 mL/min M6-6L: 0.3~6000 mL/min M6-12L: 0.69~12000 mL/min	Start/stop, direction signal	Passive switch signal, such as foot pedal Active switch signal: 24V default
Speed resolution	0.1-600 rpm	Communication interface	RS232, RS485 Modbus protocol(RTU mode)
Speed range	0-100rpm, 0.1rpm; 100-600rpm, 1rpm.	Power supply	AC 220V±10% 50Hz/60Hz (standard) AC 110V±10% 50Hz/60Hz (optional)
Flow rate accuracy	<±0.5%	Drive dimension (L×W×H)	M6-3L: 223×152×199mm M6-6L: 283×192×274mm M6-12L: 302×222×331mm
Motor type	Closed-loop stepper motor	Drive weight	M6-3L: 5.02kg; M6-6L: 7.85kg M6-12L: 13.14kg
Display	3 digital LED	Power consumption	M6-3L: <80W; M6-6L:<180W M6-12L: <300W
Control method	Digital knob control and mechanical keypad	Condition temperature	0-40°C
Keypad lifetime	300,000 times	Relative humidity	<80%
External speed control signal	0-5V, 0-10V, 4-20mA		

Product Composition and Flow Rate Range

Drive	Motor Type	Pump Head	Tubing Size	Speed Range(rpm)	Flow Rate(mL/min)
M6-3L	Closed-loop stepper motor	EasyPump	13#, 14#, 19#, 16#, 25#, 17# 18#/ 15#, 24#, 35#, 36#	0.1-600	0.0053~3100
		DZ25-3L	15#, 24#, 35#, 36#		0.211~3600
M6-6L		DZ25-6L	15#, 24#, 35#, 36#		0.3~6000
M6-12L		YZ35	26#, 73#, 82#		0.69~12000



Basic Peristaltic Pump

3 years warranty

BT100N, BT300N, BT600N




Suitable Pump Head



AMC Series



YZ Series

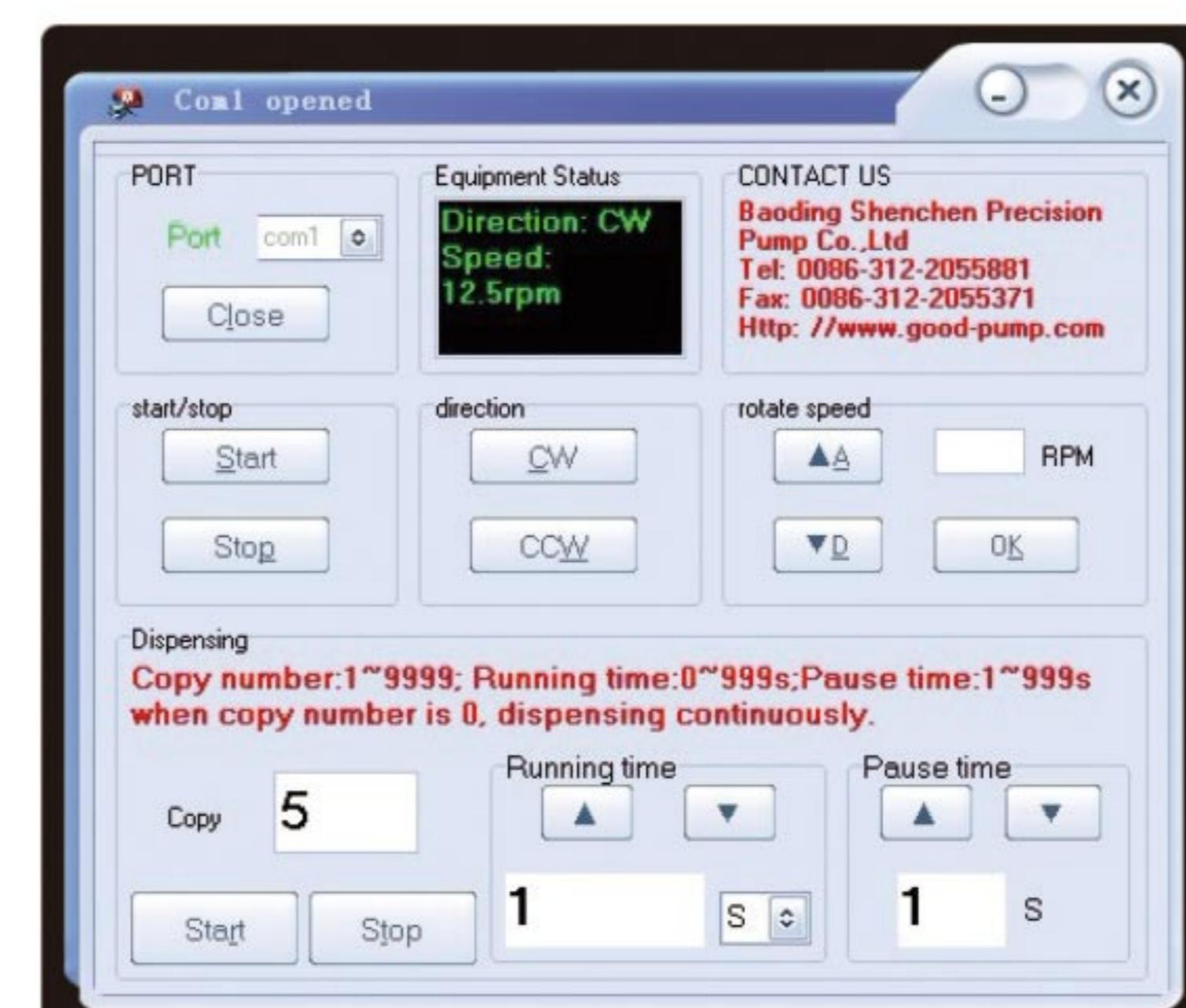

Features

- | Plastic coated metal housing, compact structure.
- | Timing function, time range 0.5s–999s, can be used for simple dispense.
- | RS232, RS485 Communication Interface.
- | Support Shenzhen communication protocol or standard Modbus communication protocol(RTU mode).

Technical Specifications

Flow rate range	BT100N: 0.0002~570 mL/min BT300N: 0.0002~1330 mL/min BT600N: 0.0002~2280 mL/min	Start/stop, reversing signal	Passive switch signal, such as foot pedal Active switch signal: 5V,12V,24V for option
Speed range	0.1-600rpm	Communication interface	RS232, RS458 communication Modbus protocol(RTU mode)
Speed resolution	0-100rpm, 0.1rpm; 100-600rpm, 1rpm	Power supply	AC 220V±10% 50Hz/60Hz (standard) AC 110V±10% 50Hz/60Hz (optional)
Testing time range	0.5 s-999 s	Drive dimension	183×131×194mm
Display	LED Display	Drive weight	4.20 kg
Control method	Mechanical keypad	Power consumption	<50W
Keypad lifetime	300,000 times	Condition temperature	0-40°C
External speed control signal	0-5V, 4-20mA, 0-10V for option	Relative humidity	<80%

Product Composition and Flow Rate Range

Flow Rates Peristaltic Pump		Pump Head & Flow Rate (mL/min)			
		YZ1515x	YZ2515x	AMC1~AMC12(10)	AMC1~AMC12(6)
Drive&speed	Tubing	13#, 14#, 19#, 16# 25#, 17#, 18#	15#, 24#	Inner diameter: 0.13-3mm Wall thickness: 0.8-1mm	
BT100N	0.1-150 rpm	0.007~570	0.17~435	0.0002-48(working speed≤150rpm)	0.0002-65(working speed≤150rpm)
BT300N	0.1-350 rpm	0.007~1330	0.17~1015		
BT600N	0.1-600 rpm	0.007~2280	0.17~1740		



Planetary Gear Type Industrial Peristaltic Pump

J025



DZ45-I



DZ45-II


Product Introduction

J025 peristaltic pump head use aluminum alloy shell, 304 stainless steel rollers, long lifetime, corrosion resistance; big flow rate, high pressure, suitable for transfer high viscosity and high lift liquid; Driven by AC motor, can use frequency adapter to adjust speed and flow rate, can also connect with PLC, IPC and computers.


Typical Application

- | Bio-medical
- | Chemical industry
- | Environmental protection

Roller Components

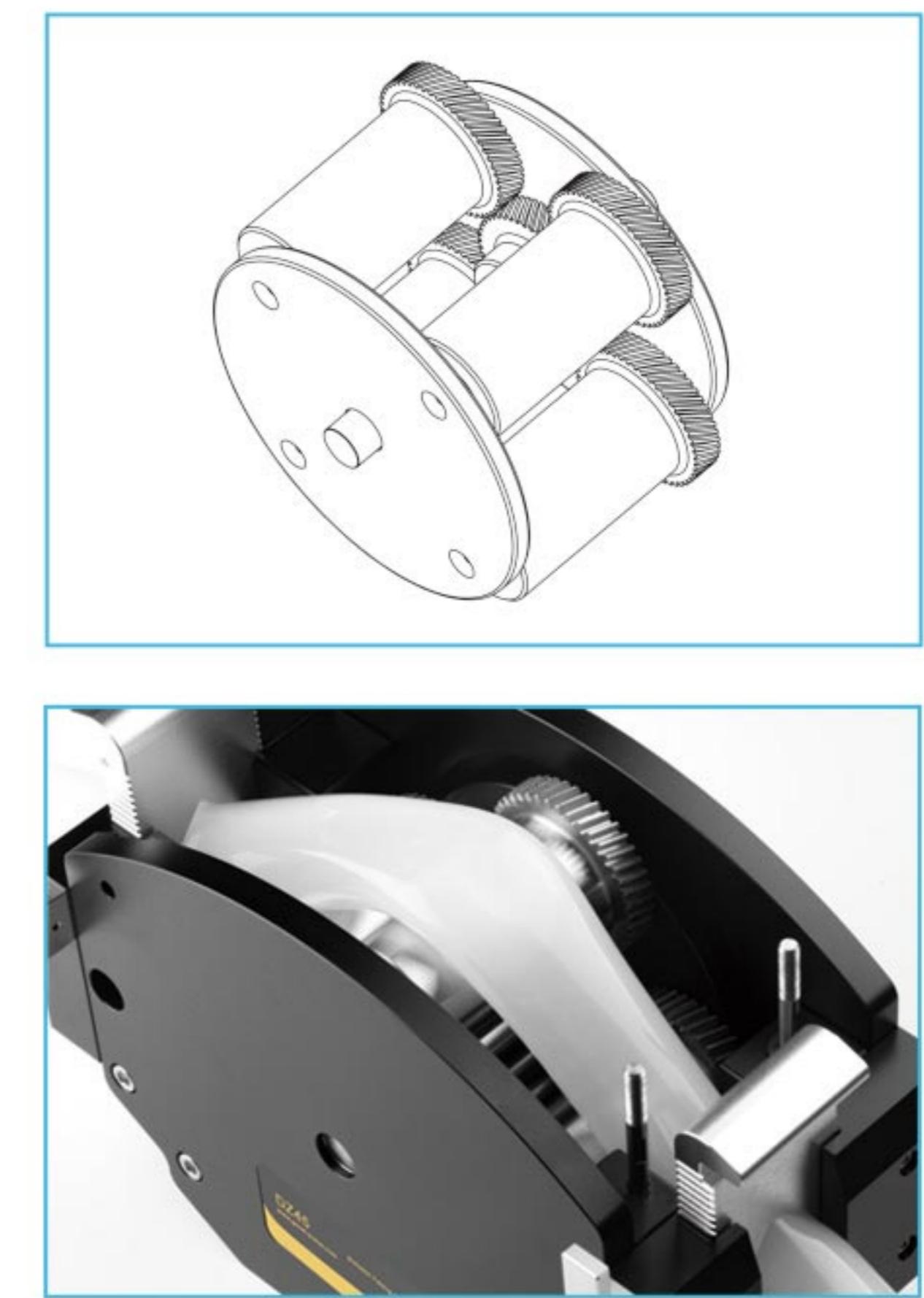

Dimension Drawing (Unit: mm)

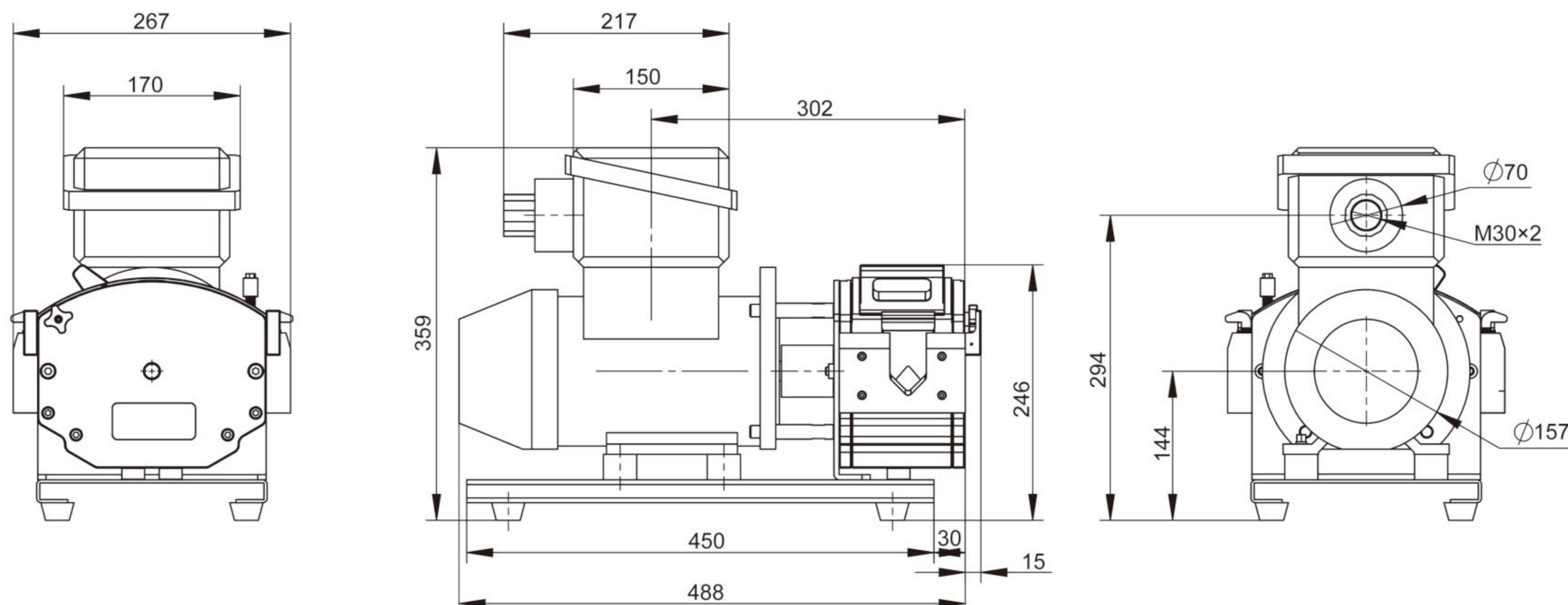

Product Composition and Flow Rate Range

Drive	Pump Head	Power Supply	Power	Speed (rpm)	Tubing		Flow Rate (L/min)	Pressure (Mpa)		Weight
					Size	ID*Wall thickness (mm)		Continuous	Intermittent	
J025	DZ45	AC380V/AC220V	370W	37.5-350	88#	12.7*4.8	1.0-12.5	0.25	0.3	30kg
				37.5-270	92#	25.4*4.8	3.68-28.15			



Explosion Proof Peristaltic Pump


EXP600/EasyPump



Features

- Driven by explosion-proof motor, the power is strong.
- Stepless speed regulation, reversible inversion.
- IP rate: IP55
- Explosion-proof grade: ExdIIBT4 Gb
- Suitable for inflammable and explosive industrial sites.

Dimension Drawing (Unit: mm)

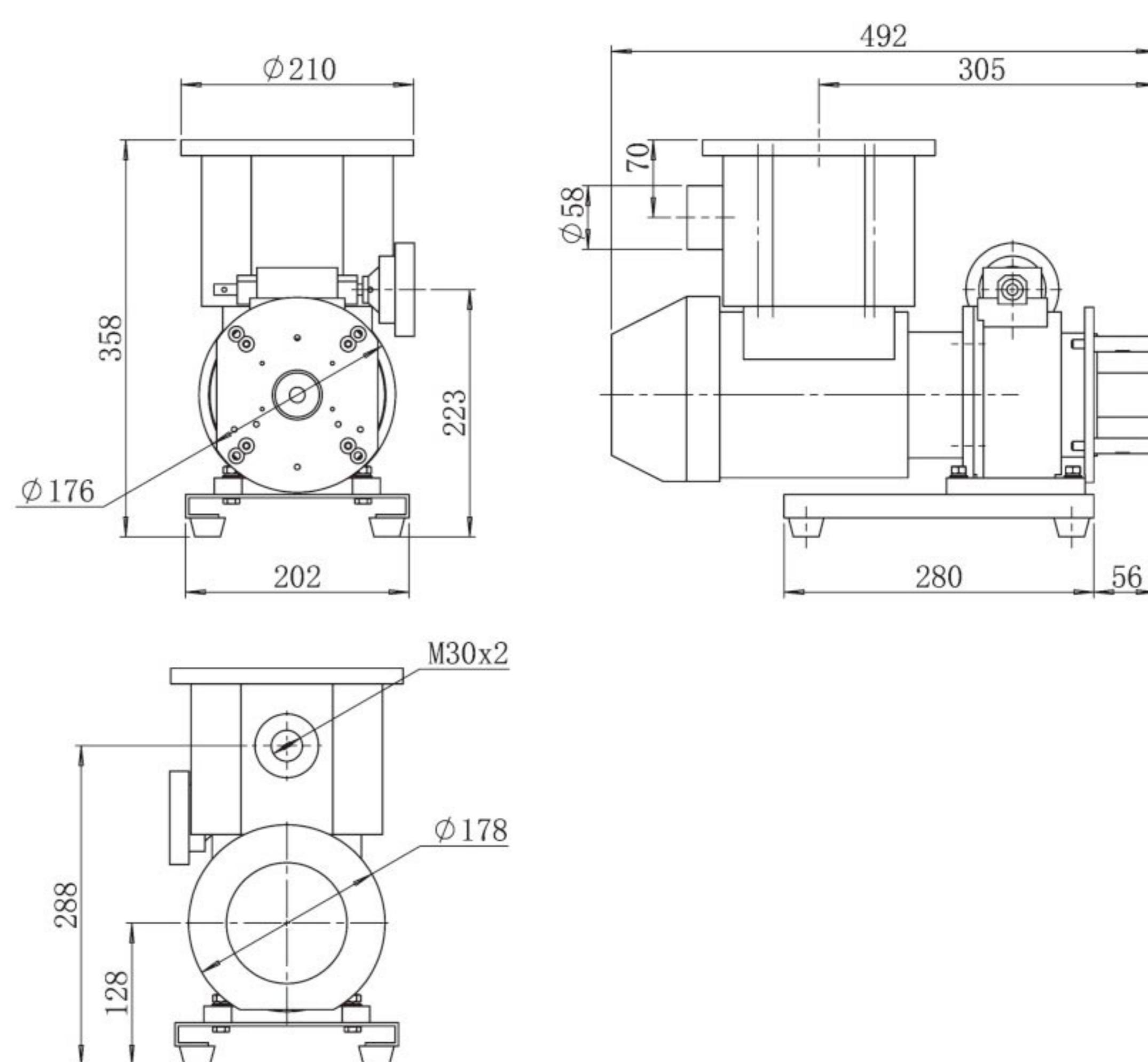

EXP600 with Frequency Converter



Technical Specifications

Speed range	Manual control: 180-600 rpm Frequency converter control: 60-600rpm	Condition temperature	0~40°C
		Relative humidity	<80%
Speed control	Manual stepless speed regulation/ Inverter speed regulation	Explosion-proof grade	ExdIIBT4 Gb, ExdIICt4
		IP rate	IP55
Voltage	3 phase 380V(standard)/ 3 phase 220V(optional)	Drive weight	30 kg
		Drive dimension(L×W×H)	424×230×330(mm)


Product Composition and Flow Rate Range

Drive	Speed(rpm)	Pump Head	Tubing Size	Flow Rate (mL/min)
EXP600	60-600	AMC	1×1, 2×1, 2.4×0.8, 3×1, 0.13×0.86, 0.19×0.86, 0.25×0.86, 0.51×0.86, 0.89×0.86, 1.14×0.86, 1.42×0.86, 2.06×0.86, 2.79×0.86	0.12~65.17
		EasyPump	13#, 14#, 19#, 16#, 25#, 17#, 18# 15#, 24#, 35#, 36#	3.18~2580 108~3100
		YZ1515x	13#, 14#, 19#, 16#, 25#, 17#, 18#	4.2~2280
		YZ2515x	15#, 24#	102~1740
		YZ35	26#, 73#, 82#	414~12000
		DZ25-3L	15#, 24#, 35#, 36#	126.6~3600
		DZ25-6L	15#, 24#, 35#, 36#	180~6000
		SN15	14#, 16#	14.4~528
		SN25	24#	150~1500
EXP300	37.5-350	DZ45	88#, 92#	1000~28150



Explosion Proof Peristaltic Pump

Suitable Pump Head

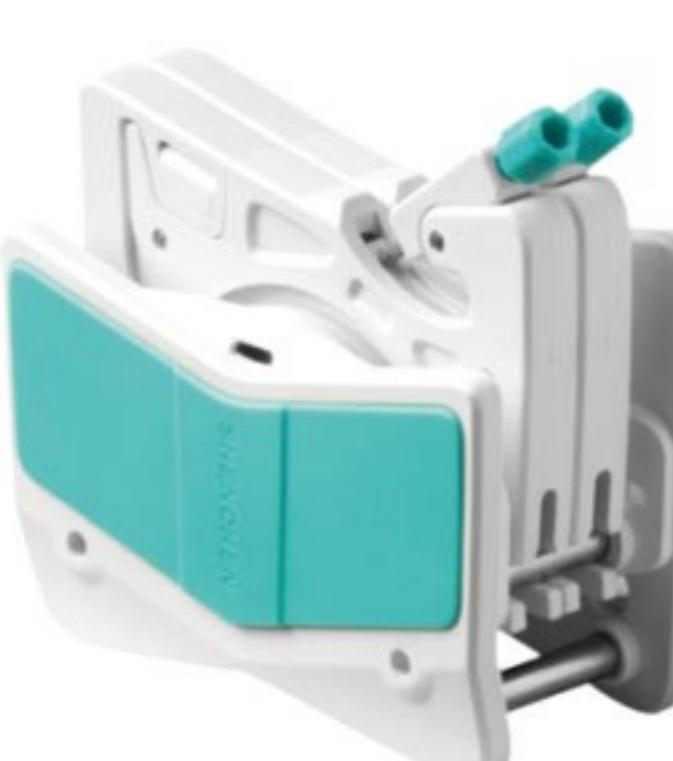

Features

Pneumatic motor driver, gas driving, explosion-proof and safe.


Stepless speed regulation, manual control flow valve. Start/stop, direction, overload protection functions.

Barometer, convenient for users to check the pressure value.


Suitable for inflammable, explosive industrial sites.

Dimension Drawing(Unit: mm)

Working status dimension drawing

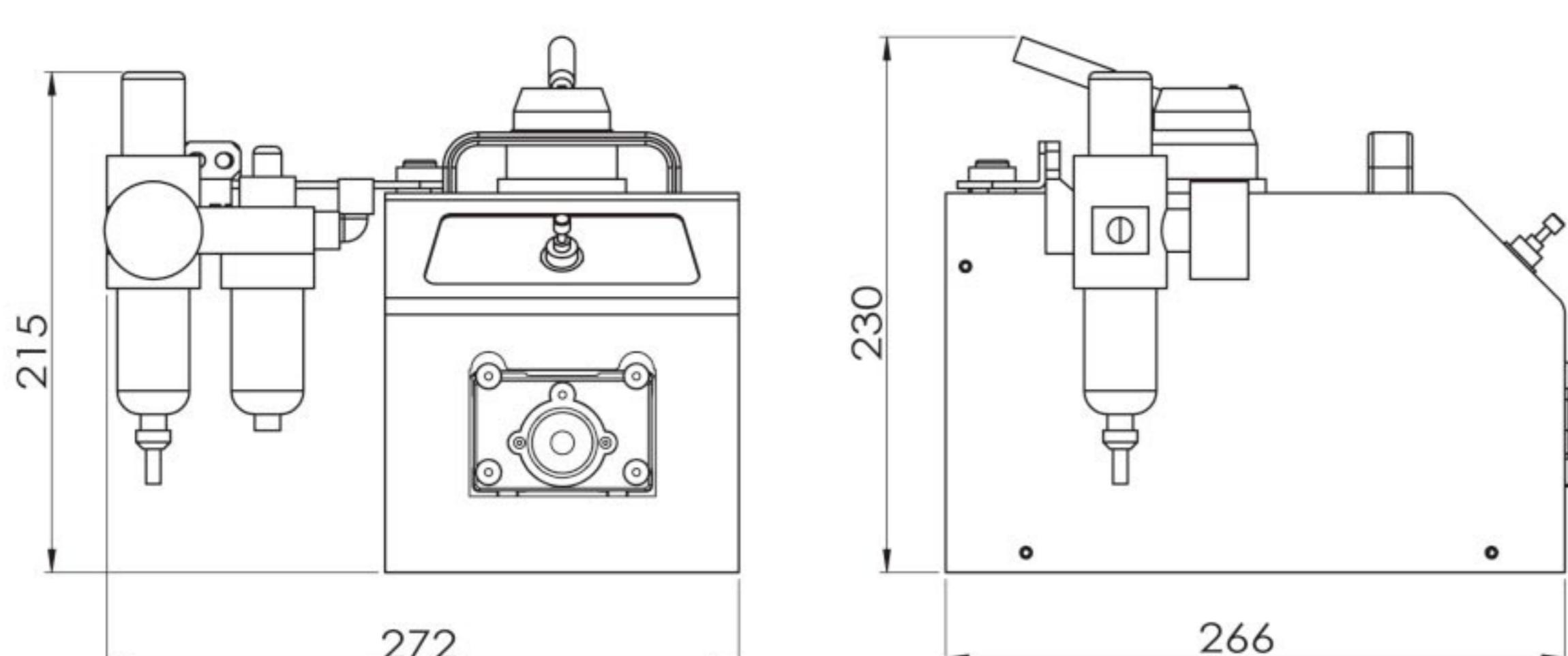

Non-working status dimension drawing

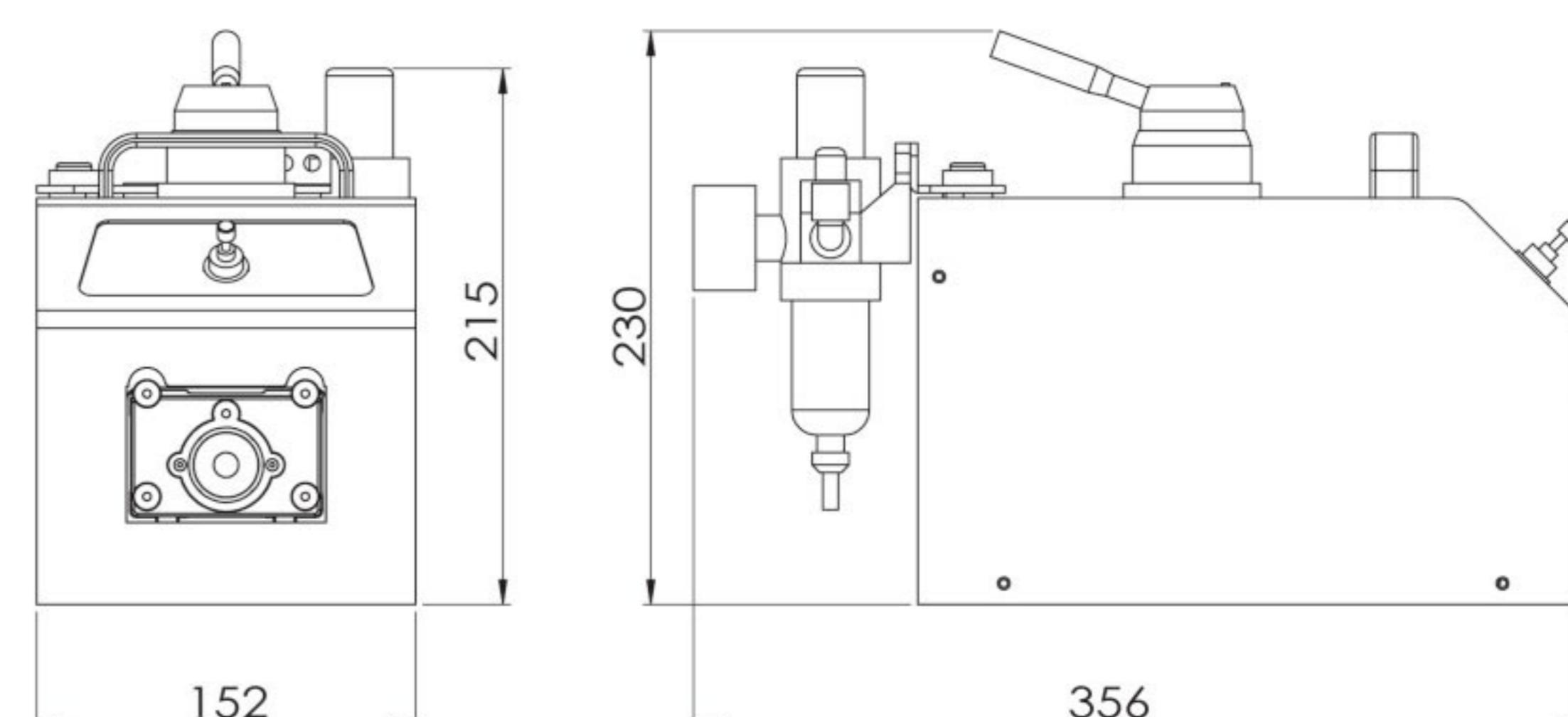

Technical Specifications

Speed range	60~600 rpm	Gas consumption	0.4-8L/Sec
Display	Display air pressure	Drive dimension (L x W x H)	Working:266x 272 x 230mm Non-working:356x 152 x 230mm
Speed control	stepless speed regulation (By adjusting the flow valve)	Drive weight	4.10 kg
Control function	Start/stop, reversing(gas reversing valve control)	Condition temperature	0-40°C
Working air pressure	0.4-0.7Mpa	Relative humidity	< 80%

Product Composition and Flow Rate Range

Peristaltic Pump	Tubing	Pump Head & Flow Rate (mL/min)			
		EasyPump		AMC1~AMC12	
Drive & Speed	QD600	EasyPumpI/III	EasyPumpII/IV	10 roller	6 roller
		13#, 14#, 19#, 16#, 25#, 17#, 18#	15#, 24#, 35#, 36#	Inner diameter: 0.13-3mm Wall thickness: 0.8-1mm	
		3.18~2580	108~3100	0.12~48.38	0.12~65.17
		YZ1515x	YZ2515x	DZ25-3L	
		13#, 14#, 19#, 16#, 25#, 17#, 18#	15#, 24#	15#, 24#, 35#, 36#	
		4.2~2280	102~1740	126.6~3600	



Pump Head Flow Meter

AMC Series Technical Specifications

Tubing	Speed	Flow rate of pump head with 10 rollers(mL/min)	Flow rate of pump head with 6 rollers(mL/min)	Tubing maximum pressure (Mpa)	
				Continuous	Intermittent
1×1	0.1~150rpm	0.0050~7.55	0.0062~9.36	0.1	0.1
2×1		0.0183~27.52	0.0220~33.06		
2.4×0.8		0.0254~38.13	0.0319~47.81		
3×1		0.0323~48.38	0.0434~65.17		
0.13×0.86		0.0002~0.29	0.0002~0.31		
0.19×0.86		0.0003~0.44	0.0003~0.46		
0.25×0.86		0.0005~0.76	0.0005~0.80		
0.51×0.86		0.0013~2.00	0.0014~2.05		
0.89×0.86		0.0030~4.47	0.0031~4.65		
1.14×0.86		0.0061~9.16	0.0065~9.74		
1.42×0.86		0.0125~18.75	0.0142~21.28		
2.06×0.86		0.0197~29.60	0.0234~35.17		
2.79×0.86		0.0286~42.86	0.0372~55.77		

EasyPump Series Technical Specifications

Pump Head	Tubing	ID×Wall thickness(mm)	mL / r	Flow Rate(mL/min) (0.1~600rpm)	Tubing Max.Pressure(Mpa)		Weight(kg)
					Intermittent	Continuous	
EasyPumpI/III	13#	0.8×1.6	0.053	0.0053~32	0.27	0.17	0.6
	14#	1.6×1.6	0.27	0.027~162			
	19#	2.4×1.6	0.55	0.055~330			
	16#	3.1×1.6	0.933	0.093~560			
	25#	4.8×1.6	1.967	0.197~1180			
	17#	6.4×1.6	3.333	0.333~2000			
	18#	7.9×1.6	4.3	0.430~2580			
EasyPumpII/IV	15#	4.8×2.4	1.8	0.180~1080	0.27	0.17	0.6
	24#	6.4×2.4	2.733	0.273~1640			
	35#	7.9×2.4	3.833	0.383~2300			
	36#	9.6×2.4	5.167	0.517~3100			
EasyPumpV/VI	13#	0.8×1.6	0.053	0.0053~32	0.27	0.17	0.6
	14#	1.6×1.6	0.27	0.027~162			
	19#	2.4×1.6	0.55	0.055~330			
	16#	3.1×1.6	0.933	0.093~560			
	25#	4.8×1.6	1.967	0.197~1180			

YZ Series Technical Specifications

Pump Head	Tubing	ID×Wall thickness(mm)	mL / r	Flow Rate(mL/min) (0.1~600rpm)	Tubing Max.Pressure(Mpa)		Material/Weight(kg)
					Intermittent	Continuous	
YZ1515x	13#	0.8×1.6	0.07	0.007~42	0.27	0.17	0.40 (3 rollers) 0.46 (3 rollers)
	14#	1.6×1.6	0.27	0.027~162			
	19#	2.4×1.6	0.55	0.055~330			
	16#	3.1×1.6	0.82	0.082~492			
	25#	4.8×1.6	1.7	0.17~1020			
	17#	6.4×1.6	2.9	0.29~1740			
	18#	7.9×1.6	3.8	0.38~2280			
YZ2515x	15#	4.8×2.4	1.7	0.17~1020	0.27	0.17	0.44 (6 rollers) 0.50 (6 rollers)
	24#	6.4×2.4	2.9	0.29~1740			

DZ25 Series Technical Specifications

Model No.	Housing Material	Tubing Clamp Material	Tubing		Flow Rate (mL/min) (0.1~600rpm)	Weight (kg)
			Tubing Sizes	ID × Wall Thickness		
DZ25-3L	PPS	PP	15#	4.8×2.4(mm)	0.211~1264	0.5
			24#	6.4×2.4(mm)	0.385~2310	
			35#	7.9×2.4(mm)	0.508~3050	
			36#	9.6×2.4(mm)	0.6~3600	
DZ25-6L	Aluminum alloy/PPS	PP	15#	4.8×2.4(mm)	0.3~1800	1.86/0.86
			24#	6.4×2.4(mm)	0.55~3300	
			35#	7.9×2.4(mm)	0.8~4800	
			36#	9.6×2.4(mm)	1~6000	



YZ35 Technical Specifications

Pump Head	Tubing	ID×Wall Thickness(mm)	mL / r	Flow Rate(mL/min) (0.1~600rpm)	Tubing Max.Pressure(Mpa)		Material/Weight(kg)	
					Intermittent	Continuous	Aluminum Alloy	PPS
YZ35	26#	6.4×3.3	6.9	0.69~4200	0.27	0.2	4.36	1.50
	73#	9.6×3.3	12.3	1.23~7400				
	82#	12.7×3.3	20	2~12000	0.14	0.1		

MiniPump Technical Specifications

Pump Head	Tube Size	ID × Wall Thickness(mm)	Speed (rpm)	Flow Rate (mL/min)
MiniPump01	13#	0.8×1.6	0.1~300	0.0024-8.28
	14#	1.6×1.6		0.0112-33.88
	19#	2.4×1.6		0.0252-77.23
	16#	3.1×1.6		0.0394-114.31
	25#	4.8×1.6		0.0652-190.00
MiniPump01/ MiniPump02	1×1	1×1	0.1~300	0.005-15.01
	2×1	2×1		0.018-54.63
	2.5×1	2.5×1		0.0256-76.84
	3×1	3×1		0.0356-108.39

HandyPump Technical Specifications

Pump Head	Channel number	Tubing	ID×Wall thickness(mm)	mL/r (3 rollers)	mL/r (6 rollers)	Speed (rpm)	Flow Rate(mL/min) (3 rollers)	Flow Rate(mL/min) (6 rollers)	Weight (kg)
HandyPump01	Single channel	13#	0.8×1.6	0.033	0.031	0.1~300	0.0033~10.03	0.0031~29.36	0.224 (3 rollers)
		14#	1.6×1.6	0.187	0.126		0.0187~56.09	0.0208~37.68	
		19#	2.4×1.6	0.371	0.215		0.0371~111.17	0.0036~64.51	0.302 (6 rollers)
		16#	3.1×1.6	0.636	0.345		0.0636~190.76	0.0059~103.51	
		25#	4.8×1.6	1.219	0.636		0.1219~365.69	0.01038~190.81	
HandyPump02	Dual channel	13#	0.8×1.6	0.033	—	0.1~300	0.0033~10.03	—	0.224 (3 rollers)
		14#	1.6×1.6	0.187	—		0.0187~56.09	—	
		19#	2.4×1.6	0.371	—		0.0371~111.17	—	
		16#	3.1×1.6	0.636	—		0.0636~190.76	—	

KT15 Technical Specifications

Pump Head	Channel number	Tubing	ID×Wall thickness(mm)	mL / r	Speed(rpm)	Flow Rate(mL/min)	Weight(kg)
KT15	Single channel	13#	0.8×1.6	0.033	0.1~600	0.0033~19.97	0.096
		14#	1.6×1.6	0.156		0.0156~93.42	
		19#	2.4×1.6	0.286		0.0286~171.6	
		16#	3.1×1.6	0.477		0.0477~286.32	
		25#	4.8×1.6	0.933		0.0933~560.04	

UD15 Technical Specifications

Model No.	Housing Material		Tubing		Speed Range(rpm)	Flow Rate (mL/min)	Weight (kg)
	Body	Protective cover	Tubing Sizes	ID×Wall Thickness			
UD15	PSF	Transparent PC	16#	3.1×1.6(mm)	0.1~350	0.08~280	0.12
			25#	4.8×1.6(mm)		0.16~580	
			17#	6.4×1.6(mm)		0.26~930	

UC25 Technical Specifications

Model No.	Housing Material		Tubing		Speed Range(rpm)	Flow Rate(mL/min)	Weight (kg)
	Base	Protective cover	Tubing Sizes	ID×Wall(mm)			
UC25	PSF	PC	15#	4.8×2.4	0.1~600	0.3423~2054	1.39
			24#	6.4×2.4		0.5033~3020	
			35#	7.9×2.4		0.8588~5153	
			36#	9.6×2.4		1.1105~6663	



Peristaltic Pump Accessories

A Filling Nozzle

Name	Material	Picture
Reducer anti-splash filling nozzle	SS316	
Flat filling nozzle	SS304/316	

B One Way Checkvalve



Avoid liquid drop off after filling and transferring.


C Filling Countersunk



Used for the output of tube, preventing the tube floating or absorbing the wall of container.

Name	Material	Tube
Counter sunk	304/316 stainless steel	13#, 14#, 19#, 16#, 25#, 17#, 18#, 15#, 24#, 35#, 36#, 26#, 73#, 82#

D Fluid Pulse Damper

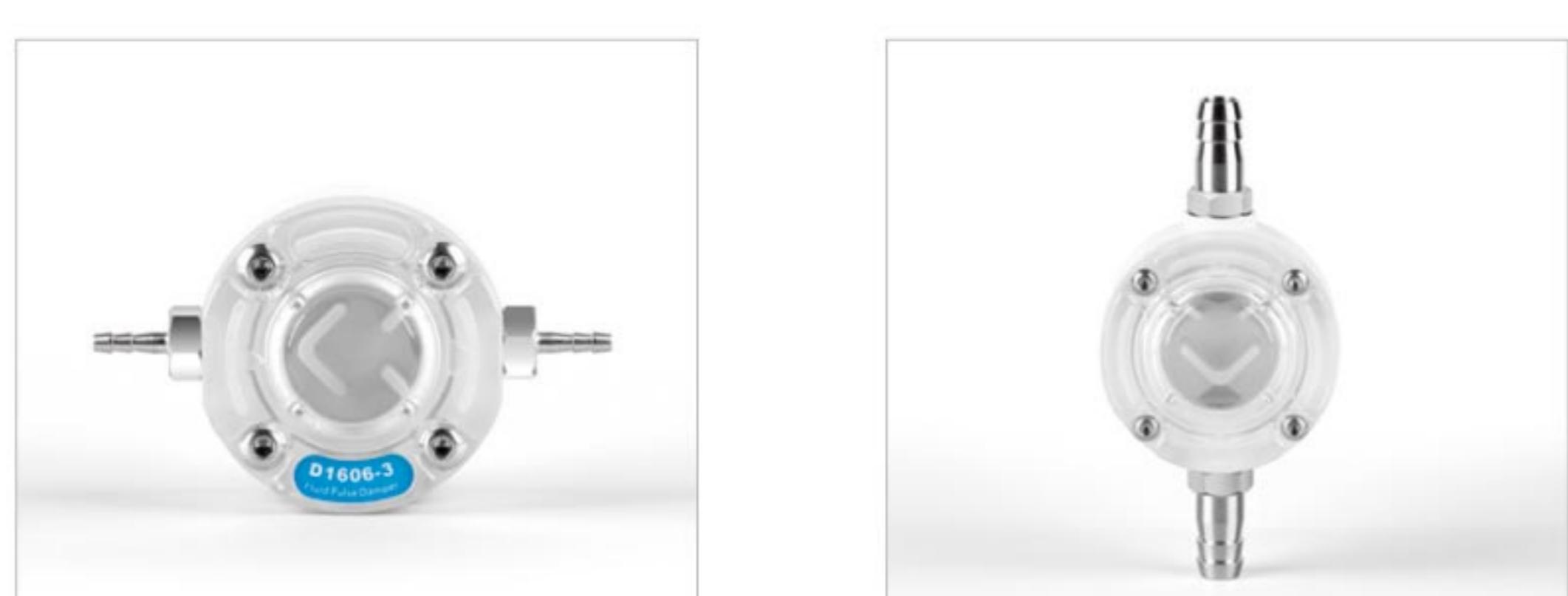

Special design for peristaltic pump, effectively suppress the peristaltic pump pulsation and improve the flow rate accuracy. The pulsation suppression rate can reach more than 95%.

E Handling Dispenser



Filling nozzle and tubing cap


Filling nozzle size	13#	14#	19#
Inner diameter	3mm	3.5mm	4.5mm
Picture			
Filling nozzle size	16#	15#/25#	17#/24#
Inner diameter	5mm	7mm	9mm
Picture			
Tubing size	17#	18#	Plum blossom cap
Inner diameter	9.6mm	11.1mm	Plum blossom cap
Picture			

Based on ergonomics design, elegant appearance, grip feeling comfortable, easy operation. Connect to peristaltic pump external control interface, with start/stop and full speed control, can realize transferring and dispensing function. Power supply and working indicator, show the dispenser working status. With hanging hole, can be hang up when do not use.

F Foot Pedal Switch



Control the pump start/stop with foot pedal switch.

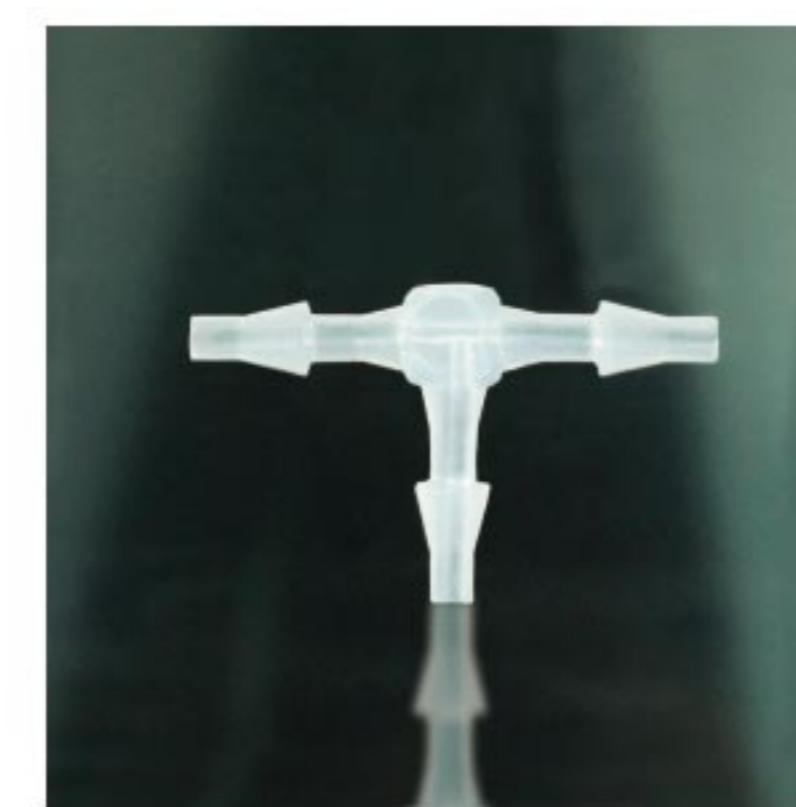G Tube Connector**a** Straight tube connector**b** "Y" tube connector**c** Tee tube connector**d** Reducer tube connector**e** "L" tube connector**f** Flexible/Hard tube connector**H PH Controller**

Work with peristaltic pump, can control the liquid PH value, add acid or alkali automatically.

Function:

1. Liquid: Acid–Base Solutions
2. PH value : 0–14PH
3. Set up target PH value
4. Add acid or alkali liquid automaticall
5. Control: RS485 , 4–20mA
6. Power supply: DC24V (AC220V for option)
7. Suitable temperature: 0–60°C

J Benchtop Tubing Cutter

Stainless steel blade, makes right-angle cuts in several sizes of plastic tubing.

K Support Stand**I 5V Sensor**

When applied in the dispensing line, it can detect weather there is filling bottle in the production line.

When the bottle approach the sensor side, the switch action will be made without any mechanical contact or pressure, thereby providing filling control order to the pump. In the same way, when no filling bottle is detected, the stop filling control order is provided to the pump.



The multiple filling stand is suitable for more than 2 channels filling. It can hold 2–8 filling nozzles. We can customize the suitable one according to your request.



Peristaltic Pump Tubing

Silicone Tubing


- | Platinum-cured silicone tubing
- | Slightly clarity, smooth surface, low protein binding levels, fewer potential leachable .
- | Ideal for pharmaceutical and biotechnology use, suitable temperature range -51~238 °C.

Micro Flow Rate Tubing										
Tubing Size	0.13×0.86	0.5×0.86	0.86×0.86	1.52×0.86	2.06×0.86	2.79×0.86	1×1	2×1	3×1	2.4×0.8
Tubing cross sections (1:1)	●	●	●	○	○	○	●	○	○	○
Wall thickness (mm)				0.86				1.0		0.8
Inside diameter (mm)	0.13	0.5	0.86	1.52	2.06	2.79	1.0	2.0	3.0	2.4
Maximum pressure (Mpa)	Continuous				0.1					
	Intermittent				0.1					

Basic Flow Rate Tubing											
Tubing Size	13#	14#	19#	16#	25#	17#	18#	15#	24#	35#	36#
Tubing cross sections (1:1)	●	●	○	○	○	○	○	○	○	○	○
Wall thickness	mm			1.6						2.4	
	inch			1/16						3/32	
Inside diameter	mm	0.8	1.6	2.4	3.2	4.8	6.4	7.9	4.8	6.4	9.6
	inch	1/32	1/16	3/32	1/8	3/16	1/4	5/16	3/16	1/4	5/16
Maximum pressure (Mpa)	Continuous		0.17		0.14	0.1	0.07		0.17		0.14
	Intermittent		0.27		0.24	0.14	0.1		0.27		0.24

Industrial Tubing							
Tubing Size	26#	73#	82#	86#	90#	88#	92#
Tubing cross sections (1:1)							
Wall thickness	mm		3.3		6.4		4.8
	inch		1/8		1/4		3/16
Inside diameter	mm	6.4	9.6	12.7	9.5	19	12.7
	inch	1/4	3/8	1/2	3/8	3/4	1/2
Maximum pressure (Mpa)	Continuous		0.2			0.25	
	Intermittent		0.27			0.3	



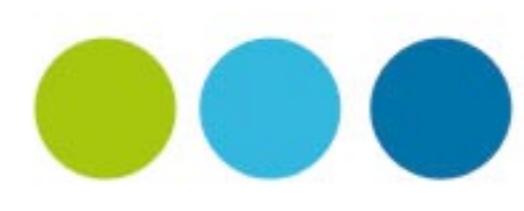


Peristaltic Pump Tubing

SAINT-GOBAIN Tubing: Tygon, PharMed BPT, Norprene etc

	A Tygon3350	B Tygon E-3603	C Norprene Chemical	D PharMed	E Norprene A-60-F
Formulation	Tygon3350	Tygon R-3603	Norprene Chemical	PharMed	Norprene A-60-F
Application	Pharmaceutical, cosmetic, medical and auto-analysis application.	General laboratory, food & beverage, biopharmaceutical, analytical instruments.	Excellent for chemical processing and general industrial applications. Food and beverage applications where extractables are a concern.	Cell and tissue culture work and pharmaceutical uses. Also good for light-sensitive samples.	Ideal for the food, dairy and beverage.
Advantages	Ultra-smooth; minimizes bacterial growth. Good for mild to medium concentration bases, salts and alcohols; odorless, tasteless, and nontoxic. Transparent.	Inexpensive tubing for general lab application. Nonaging, nonoxidizing. Clear for easy flow monitoring. Handles virtually all inorganic chemicals. Low gas permeability. Smooth bore; good for viscous fluids. High dielectric constant.	Norprene thermoplastic elastomer outer jacket with chemically inert Tygon® 2075 inner bore for excellent chemical resistance. Plasticizer-free to guard against extractables. Long flex life. Opaque beige.	Great for tissue and cell work-nontoxic and nonhemolytic; long service life minimizes risk of fluid exposure; reduces tubing costs and pump downtime. Opaque to UV and visible light to protect light-sensitive fluids. Heat sealable, bondable, and formable. Extremely low gas permeability.	Heat, ozone, and UV light resistant. Nonaging; nonoxidizing; superior acid and alkali resistance. Opaque beige.
Application Suitability	ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE GOOD VACUUM GOOD VISCOUS FLUIDS EXCELLENT STERILE FLUIDS GOOD	ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE GOOD VACUUM EXCELLENT VISCOUS FLUIDS GOOD STERILE FLUIDS EXCELLENT	ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE GOOD VACUUM EXCELLENT VISCOUS FLUIDS GOOD STERILE FLUIDS EXCELLENT		
Physical characteristics	Thermoplastic. PVC-based material with plasticizer. Firm (stiff) material. Transparent, clear.			Thermoplastic elastomer. Polypropylene-based material with USP mineral oil. Excellent tensile strength. Firm (stiff) material. Opaque, beige.	
Temp. range	-75 to 450° F (-60~232° C)	-58 to 165° F (-50~74° C)	-76 to 165° F (-60~74° C)	-60 to 270° F (-59~135° C)	-60 to 275° F (-51~135° C)
Meets classifications	FDA 21 CFR 177.2600 USP Class VI EP 3.1.9. Exceeds 3A standards Manufactured according to GMP.	FDA 21 CFR 175.300	None.	None.	FDA 21 CFR 177.2600 NSF listed (Standard 51) Manufactured according to GMP.
Cleaning/ Sterilization	Ethylene oxide gamma irradiation, or autoclave for 30 min, 15psi (1 bar).	Unaffected by commercial sanitizers (with recommended procedures) Sterilize with ethylene oxide (ETO) or autoclave. To autoclave: Coil loosely in nonlinting cloth or paper, autoclave at 121°C (250°F). 1KG/cm³ (15psi) for 30 minutes (tubing will appear milky); air dry at max 66°C (150°F) for 2 to 2 ½ hours until clear.	Sterilize with ethylene-oxide(ETO), autoclave or gamma irradiation up to 2.5Mrad. Repeated autoclaving will not affect overall life.	Autoclave, ethylene oxide, or gamma irradiation.	Autoclave.



Peristaltic Pump Tubing

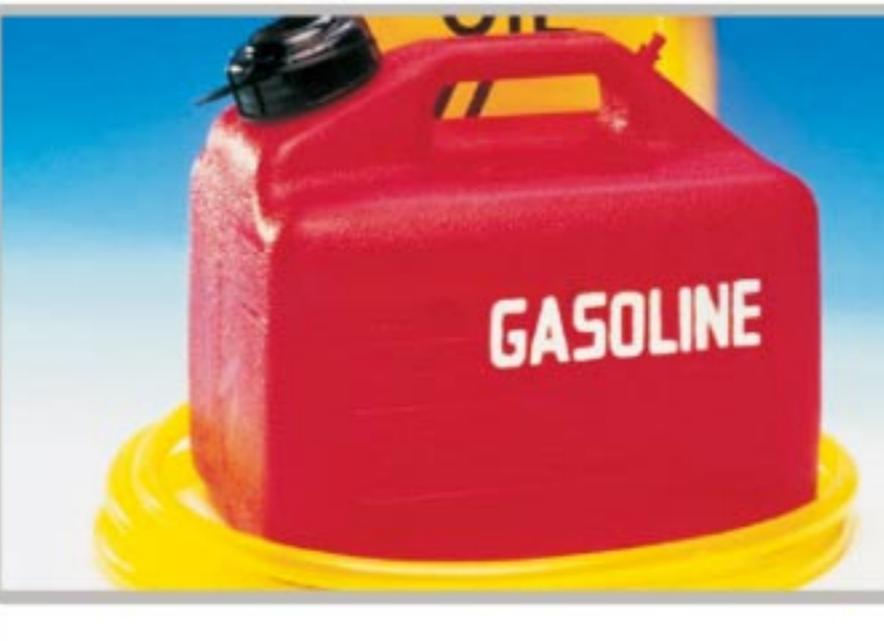	 F Norprene A-60-G	 G Tygon F-4040-A	 H Tygon LFL	 I T Y G O N 2475	 K Viton
Formulation	Norprene A-60-G	Tygon F-4040-A	Tygon LFL	TYGON 2475	Viton
Application	For applications requiring excellent chemical, heat, ozone, and ultra-violet (UV) light resistance.	Fuels and industrial lubricants-gasoline, kerosene, heating oils, cutting compounds, and glycol-based coolants. Resists most hydrocarbons.	General laboratory use, provides longer life with peristaltic tubing pumps.	Sensitive fluid transfer applications requiring high purity.	Acid and solvent transfer, high-temperature.
Advantages	Best choice for vacuum/pressure applications. Offers longest life with good flow consistency. Heat and ambient ozone resistant. Good resistance to acids/alkalies. Black color hides dirt and dust. Heat sealable, nonaging, and nonoxidizing. High dielectric constant.	Resists embrittlement and swelling, ozone-and UV-resistant, with low-extractability. Translucent yellow.	Longest life of all Tygon® peristaltic tubing (1000hrs). Nonaging, nonoxidizing. Clear for easy flow monitoring. Broad chemical resistance; low gas permeability. Smooth bore. Good for viscous fluids. High dielectric constant.	Plasticizer free, smooth inner surface (inhibits particulate buildup and bacterial growth), safely disposed of through incineration and nontoxic. Transparent.	The most chemical resistant tubing. Registand to corrosives, solvents, and oils at elevated temperatures. Low gas permeability.
Application Suitability	ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE EXCELLENT VACUUM EXCELLENT VISCOUS FLUIDS EXCELLENT STERILE FLUIDS NO		ACIDS GOOD ALKALIES GOOD ORGANIC SOLVENTS NO PRESSURE GOOD VACUUM GOOD VISCOUS FLUIDS EXCELLENT STERILE FLUIDS POOR		ACIDS EXCELLENT ALKALIES EXCELLENT ORGANIC SOLVENTS EXCELLENT PRESSURE GOOD VACUUM GOOD VISCOUS FLUIDS GOOD STERILE FLUIDS FAIR
Physical characteristics	Thermoplastic elastomer. Polypropylene-based material with USP mineral oil. Excellent tensile strength. Firm (stiff) material. Opaque, black. Manufactured according to GMP.		Thermoplastic. PVC-based material with plasticizer. Firm (stiff) material. Transparent, clear.		Thermal set rubber. Viton B (67% fluorine) Firm (stiff) material Opaque, black. Manufactured according to GMP.
Temp. range	-60 to 270° F (-59~135° C)	-35 to 165° F (-37~74° C)	-58 to 165° F (-50~74° C)	-94 to 125° F (-70~52° C)	-25 to 400° F (-32~205° C)
Meets classifications	None.	Meets NSF-51 and 3A sanitary standards.	USP Class VI, FDA 21 CFR 175.300	FDA 21 CFR 177.1520, USP 23 Class VI, Manufactured according to GMP.	None.
Cleaning/ Sterilization	Sterilize by autoclave only.	Not recommended.	Sterilize by ETO/autoclave. Coil loosely in nonlinting cloth or paper; autoclave at 250°F(121°C), 15 psi (1kg/cm²), 30 minutes (tubing will appear milky); air dry at max 150°F (66°C) for 2 to 2 ½ hrs until clear.	Ethylene oxide or gamma irradiation.	Sterilization is not recommended.





BLANC - LABO S.A.

Chemin des Mouettes 1
Centre du Bief,
CH - 1027 Lonay

T. +41 (0)21 804 18 50
F. +41 (0)21 804 18 59
E. sales@blanc-labo.com
I. www.blanc-labo.com