




COLD THAT CARES

evermed

WEB APP

DASHBOARD

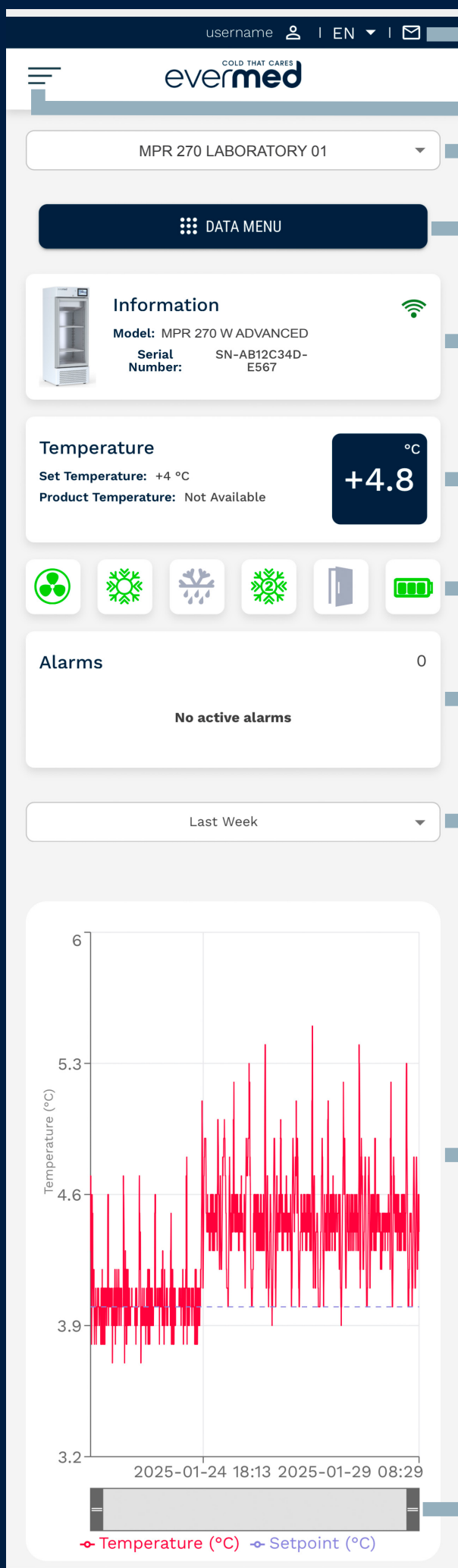

User menu (fig. 2)

Menu (fig. 3)

Name of the appliance (customizable)

Data menu (fig. 4)

Appliance informations:
Model
Serial Number
Wi-Fi connection status

Temperature informations:
Set temperature
Product temperature
Real time temperature

Appliance status icons.
Green icons: in function
Grey icons: not in function


- Fan
- Compressor 1
- Defrosting
- Compressor 2
- Door ajar
- Battery status

Real time active alarms

Select the time interval
(day - last week - last month)


Temperature graph
Touching a point on the graph it will be displayed the temperature and time

By moving the sliders on the bar you can zoom the graph

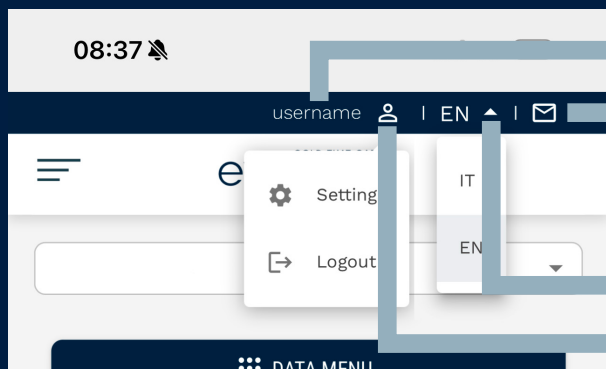


(fig. 2)

- Name of the user registered
- Request assistance via email
- Select language
- Logout and user settings:
Name
Email
Phone number
Date of birth
City
Change password

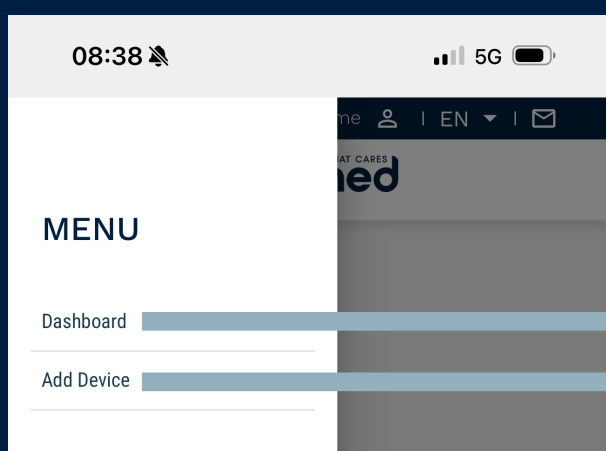

(fig. 3)

- Return to Dashboard
- Add a new Device

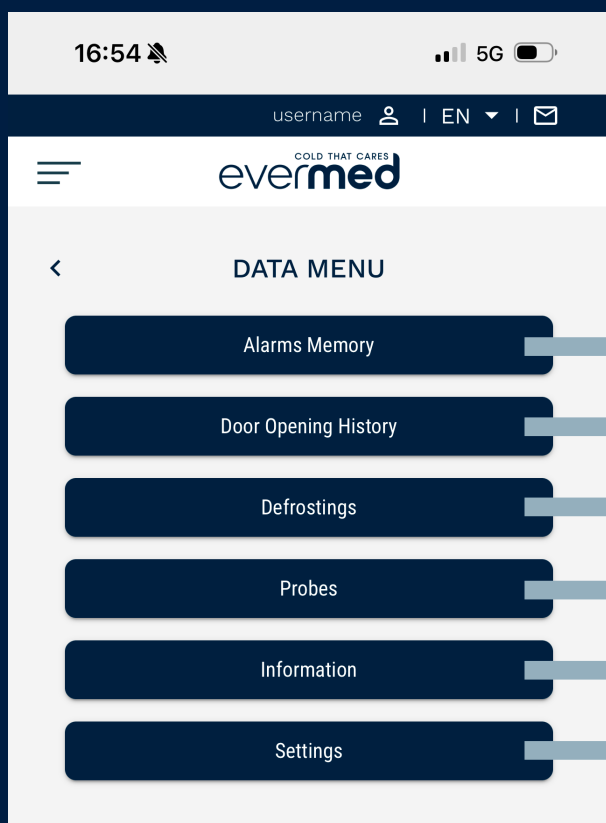

(fig. 4)


- Alarms Memory (fig. 5)
- Door Opening History (fig. 6)
- Defrosting (fig. 7)
- Probes (fig. 8)
- Information (fig. 9)
- Settings (fig. 10)

USER MENU

MENU

DATA MENU


ALARMS MEMORY

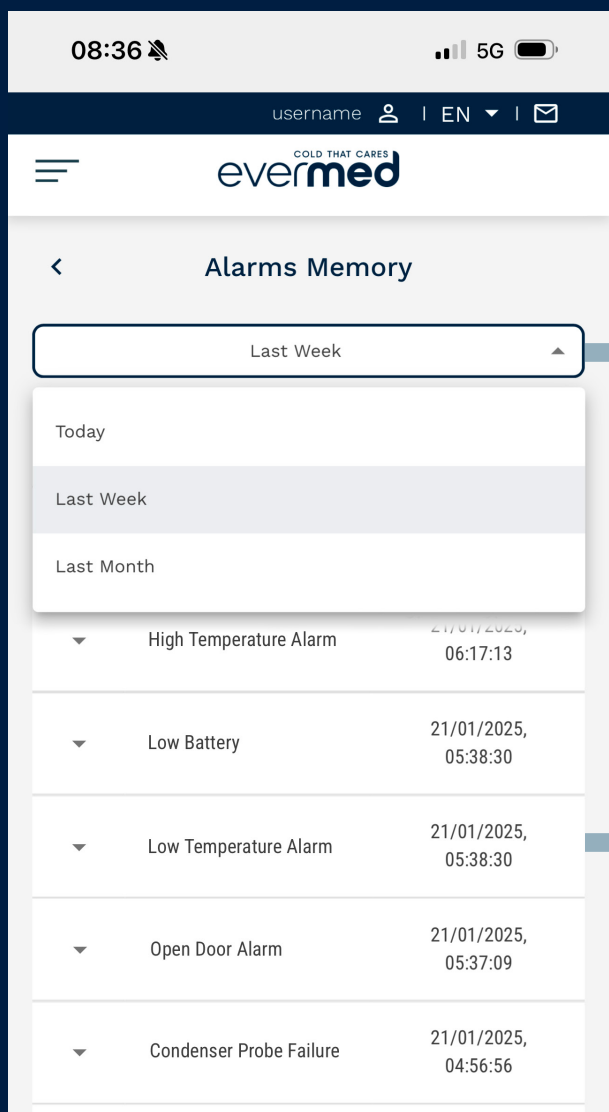

(fig. 5)

Select the time interval
(day - last week - last month)

List of stored alarms

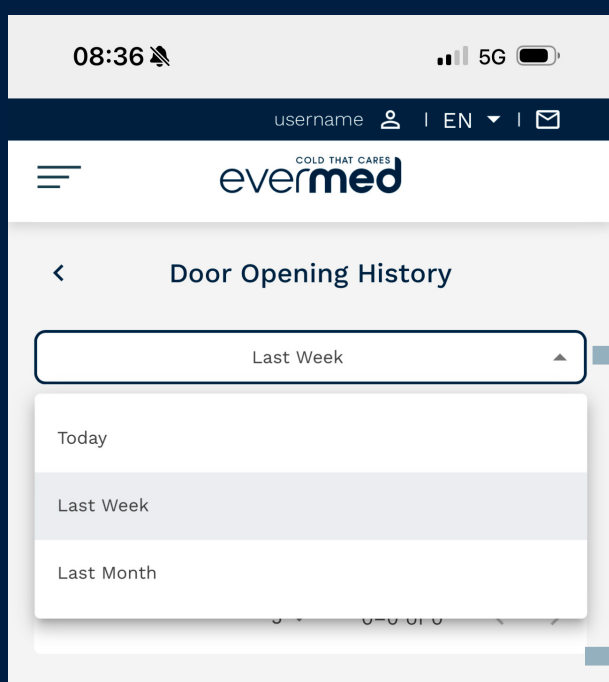


(fig. 6)

Select the time interval
(day - last week - last month)

List of door opening

DOOR OPENING HISTORY

DEFROSTING

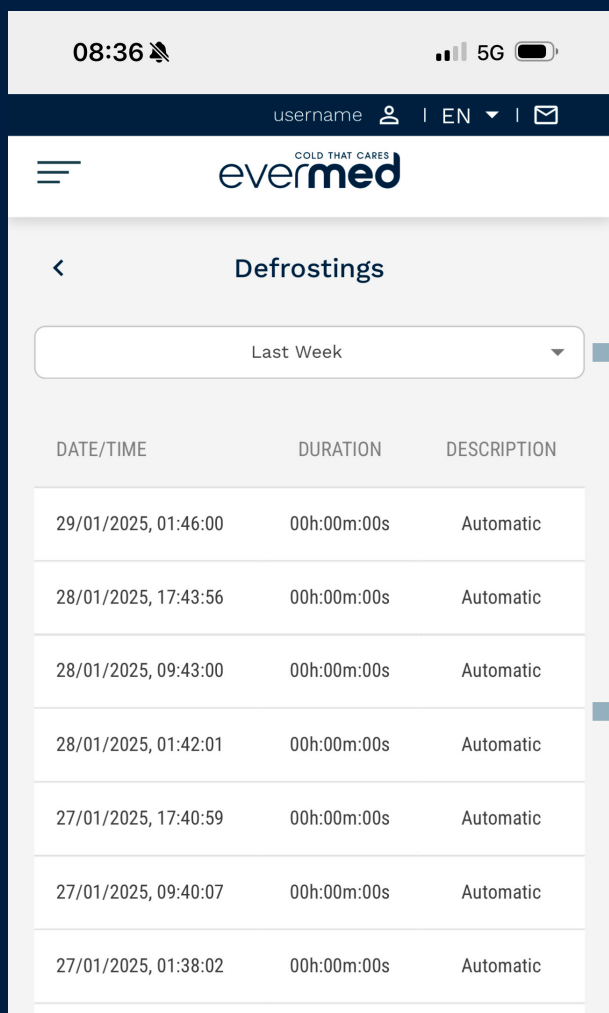


(fig. 7)

Select the time interval
(day - last week - last month)

List of stored defrostings

PROBES

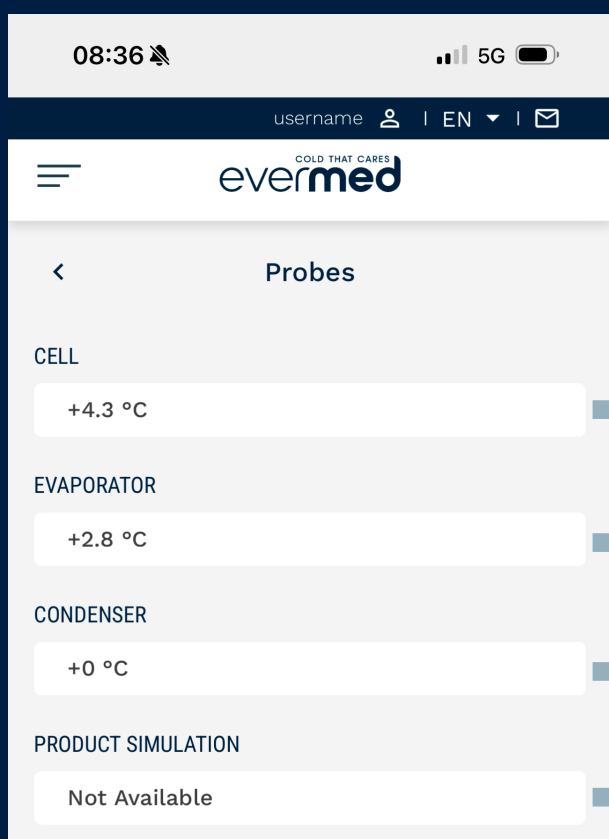


(fig. 8)

Temperature of the CELL PROBE

Temperature of the EVAPORATOR PROBE


Temperature of the CONDENSER PROBE

Temperature of the PRODUCT SIMULATION
PROBE

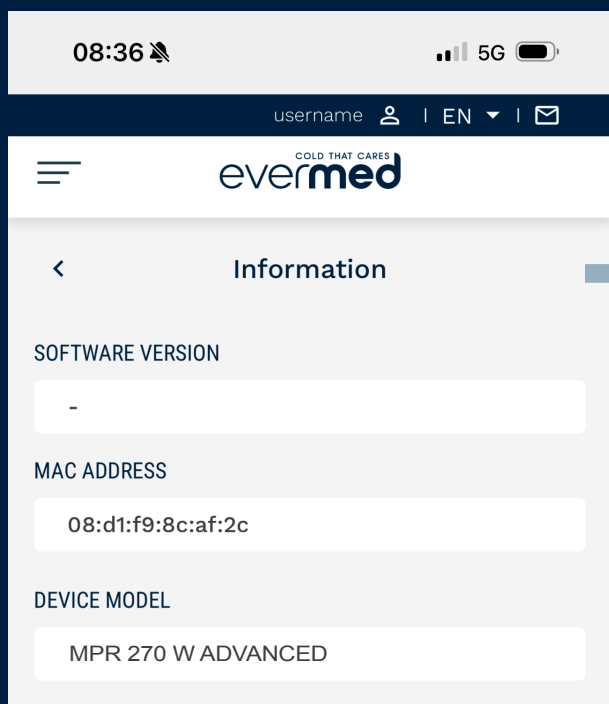

(fig. 9)

Information about the software and the device model

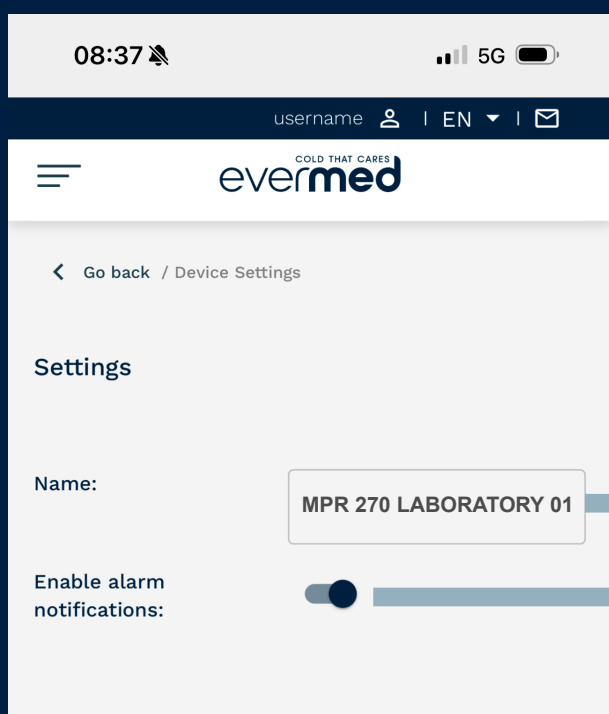

(fig. 10)

Name of the appliance (customizable)

Enable alarm notifications by e-mail

INFORMATION

SETTINGS